

GMP PE/Dazzle™ 594 anti-human CD64 Antibody

Catalog# / Size	260248 / 100 tests
Clone	10.1
Workshop	VI MA36
Other Names	FcγRI, FcR I
Isotype	Mouse IgG1, κ
Description	CD64 is a 72 kD single chain type I glycoprotein also known as FcγRI and FcR I. CD64 is a member of the immunoglobulin superfamily and is expressed on monocytes/macrophages, dendritic cells, and activated granulocytes. The expression can be upregulated by IFN-γ stimulation. CD64 binds IgG immune complex. It plays a role in antigen capture, phagocytosis of IgG/antigen complexes, and antibody-dependent cellular cytotoxicity (ADCC).

Product Details

Reactivity	Human
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	Human rheumatoid synovial fluid cells and fibronectin-purified monocytes.
Formulation	Phosphate-buffered solution, pH 7.2, containing True-Stain Monocyte Blocker™, 0.09% sodium azide and 0.2% (w/v) BSA (origin USA) and a stabilizer.
Preparation	The antibody was purified by affinity chromatography and conjugated with PE/Dazzle™ 594 under optimal conditions.
Concentration	100 µg/mL
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application. * PE/Dazzle™ 594 has a maximum excitation of 566 nm and a maximum emission of 610 nm.
Excitation Laser	Blue Laser (488 nm) Green Laser (532 nm)/Yellow-Green Laser (561 nm)
Application Notes	Clone 10.1 recognizes the EC3 epitope of CD64. While both contain the EC3 domain, in-house testing suggests that clone 10.1 preferentially binds to CD64A (FcγRIA), but not CD64B (FcγRIB). Additional reported applications (for the relevant formats) include: blocking of human IgG3 and murine IgG2a binding to FcγRI ^{2,5,6,11} and immunohistochemical staining of acetone-fixed frozen tissue sections ¹² .
Application References	<ol style="list-style-type: none"> McMichael A, <i>et al.</i> Eds. 1987. Leucocyte Typing III. Oxford University Press. New York. Schlossman S, <i>et al.</i> Eds. 1995. Leucocyte Typing V. Oxford University Press. New York. p. 874. Kishimoto T, <i>et al.</i> Eds. 1997. Leucocyte Typing VI. Garland Publishing Inc. London. Holl V, <i>et al.</i> 2004. <i>J. Immunol.</i> 173:6274. Hober D, <i>et al.</i> 2002. <i>J. Gen. Virol.</i> 83:2169. Cho HJ, <i>et al.</i> 2007. <i>Physiol Genomics</i> 149:60. van Tits L, <i>et al.</i> 2005. <i>Arterioscler Thromb Vasc Biol.</i> 25:717. PubMed Bruhns P, <i>et al.</i> 2008. <i>Blood</i> 113:3716. PubMed Yoshino N, <i>et al.</i> 2000. <i>Exp. Anim. (Tokyo)</i> 49:97. (FC) Carter DL, <i>et al.</i> 1999. <i>Cytometry</i> 37:41. (FC) Dougherty GJ, <i>et al.</i> 1987. <i>Eur. J. Immunol.</i> 17:1453.
(PubMed link indicates BioLegend citation)	

Disclaimer

GMP RUO Flow Cytometry Antibodies. BioLegend GMP RUO fluorophore conjugated antibodies are manufactured in a dedicated GMP facility and compliant with ISO 13485:2016. For research use only. Not for use in diagnostic or therapeutic procedures. Our processes include:

- Batch-to-batch consistency
- Material traceability
- Documented procedures
- Documented employee training
- Equipment maintenance and monitoring records
- Lot-specific certificates of analysis
- Quality audits per ISO 13485:2016
- QA review of released products

Antigen Details

Structure	Ig superfamily, type I glycoprotein, 72 kD
Distribution	Monocytes, macrophages, dendritic cells, activated granulocytes
Function	Phagocytosis, ADCC
Ligand/Receptor	IgG receptor
Cell Type	Dendritic cells, Granulocytes, Macrophages, Monocytes
Biology Area	Immunology, Innate Immunity
Molecular Family	CD Molecules, Fc Receptors
Antigen References	1. Hulett M, et al. 1994. <i>Adv. Immunol.</i> 57:1. 2. van de Winkel J, et al. 1993. <i>Immunol. Today</i> 14:215.
Gene ID	2209


Related Protocols

[Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

Biotin anti-human CD64, FITC anti-human CD64, PE anti-human CD64, Purified anti-human CD64, Alexa Fluor® 488 anti-human CD64, Alexa Fluor® 647 anti-human CD64, APC anti-human CD64, Pacific Blue™ anti-human CD64, Brilliant Violet 421™ anti-human CD64, PE/Cyanine7 anti-human CD64, PerCP/Cyanine5.5 anti-human CD64, APC/Cyanine7 anti-human CD64, Brilliant Violet 510™ anti-human CD64, Purified anti-human CD64 (Maxpar® Ready), PE/Dazzle™ 594 anti-human CD64, Brilliant Violet 605™ anti-human CD64, APC/Fire™ 750 anti-human CD64, TotalSeq™-A0162 anti-human CD64, Brilliant Violet 711™ anti-human CD64, Alexa Fluor® 700 anti-human CD64, Brilliant Violet 785™ anti-human CD64, TotalSeq™-C0162 anti-human CD64, Ultra-LEAF™ Purified anti-human CD64, TotalSeq™-B0162 anti-human CD64, TotalSeq™-D0162 anti-human CD64, GMP PE anti-human CD64, GMP FITC anti-human CD64, Brilliant Violet 650™ anti-human CD64

Product Data

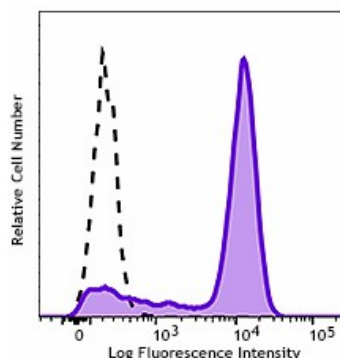

Typical results from human peripheral blood monocytes stained either with 10.1 PE/Dazzle™ 594 used at 5 µL/test (filled histogram) or with an isotype control (open histogram).

For research use only. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587