

## Brilliant Violet 421™ anti-mouse CD62L Antibody

<b>Catalog# / Size</b>	104435 / 125 µL 104436 / 50 µg
<b>Clone</b>	MEL-14
<b>Regulatory Status</b>	RUO
<b>Other Names</b>	L-selectin, LECAM-1, Ly-22, LAM-1, MEL-14
<b>Isotype</b>	Rat IgG2a, κ
<b>Description</b>	CD62L is a 74-95 kD glycoprotein also known as L-selectin, LECAM-1, Ly-22, LAM-1, and MEL-14. It is a member of the selectin family and is expressed on the majority of B and naïve T cells, a subset of memory T cells, monocytes, granulocytes, most thymocytes, and a subset of NK cells. CD62L is important in lymphocyte homing to high endothelial venules (HEV) in peripheral lymph nodes and leukocyte "rolling" on activated endothelium. CD62L also contributes to neutrophil emigration at inflammatory sites. CD62L is rapidly shed from lymphocytes and neutrophils upon cellular activation and the expression levels of CD62L (in conjunction with other markers) have been used to distinguish naïve, effector, and memory T cells. CD62L has been reported to interact with CD34, GlyCAM-1, and MAdCAM-1.

### Product Details

---

<b>Verified Reactivity</b>	Mouse
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Rat
<b>Immunogen</b>	C3H/eb mouse B lymphoma 38C-13
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA).
<b>Preparation</b>	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 421™ under optimal conditions.
<b>Concentration</b>	µg sizes: 0.2 mg/mL µL sizes: lot-specific (to obtain lot-specific concentration, please enter the lot number in our <a href="#">Concentration and Expiration Lookup</a> or <a href="#">Certificate of Analysis</a> online tools.)
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	<p>Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a>. For immunofluorescent staining using the µg size, the suggested use of this reagent is ≤0.125 µg per million cells in 100 µl volume. For immunofluorescent staining using the µl size, the suggested use of this reagent is 5 µl per million cells in 100 µl staining volume or 5 µl per 100 µl of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.</p> <p>Brilliant Violet 421™ excites at 405 nm and emits at 421 nm. The standard bandpass filter 450/50 nm is recommended for detection. Brilliant Violet 421™ is a trademark of Sirigen Group Ltd.</p> <p><a href="#">Learn more about Brilliant Violet™.</a></p> <p>This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.</p>
<b>Excitation Laser</b>	Violet Laser (405 nm)

**Application Notes** Additional reported applications (for the relevant formats) include: immunoprecipitation<sup>1-3</sup>, complement-dependent cytotoxicity<sup>4</sup>, *in vivo* and *in vitro* blocking of adhesion<sup>1-3,5</sup>, and immunohistochemical staining of acetone-fixed frozen sections and zinc-fixed paraffin-embedded sections<sup>6</sup>. The Ultra-LEAF™ purified antibody (Endotoxin < 0.01 EU/μg, Azide-Free, 0.2 μm filtered) is recommended for functional assays (Cat. Nos. 104457-104462).

#### Application References

(PubMed link indicates BioLegend citation)

1. Gallatin WM, *et al.* 1983. *Nature* 304:30. (IP, Block)
2. Siegelman MH, *et al.* 1990. *Cell* 61:611. (IP, Block)
3. Lewinsohn DM, *et al.* 1987. *J. Immunol.* 138:4313. (IP, Block)
4. Iwabuchi K, *et al.* 1991. *Immunobiology* 182:161. (CMCD)
5. Pizcueta P, *et al.* 1994. *Am. J. Pathol.* 145:461.
6. Reichert RA, *et al.* 1986. *J. Immunol.* 136:3535. (IHC, FC)
7. Olver S, *et al.* 2006. *Cancer Res.* 66:571.
8. Fukushima A, *et al.* 2006. *Invest. Ophthalmol. Vis. Sci.* 47:657. [PubMed](#)
9. Benson MJ, *et al.* 2007. *J. Exp. Med.* doi:10.1084/jem.20070719. (FC) [PubMed](#)
10. Chappaz S, *et al.* 2007. *Blood* doi:10.1182/blood-2007-02-074245. (FC) [PubMed](#)
11. Lee JW, *et al.* 2006. *Nature Immunol.* 8:181.
12. Shigeta A, *et al.* 2008. *Blood* 112:4915 (FC) [PubMed](#)
13. de Vries VC, *et al.* 2009. *Am. J. Transplant.* 9:2270 [PubMed](#)

#### Product Citations

1. Bergot AS, *et al.* 2020. *J Immunol.* 204:1787. [PubMed](#)
2. Valanparambil RM, *et al.* 2017. *PLoS Pathog.* 13:e1006647. [PubMed](#)
3. Skinner DD, *et al.* 2022. *Front Immunol.* 13:931388. [PubMed](#)
4. Aldon Y, *et al.* 2020. *J Immunol.* 204:903. [PubMed](#)
5. He Y, *et al.* 2020. *Sci Adv.* 6:eaba7589. [PubMed](#)
6. Limon JJ *et al.* 2019. *Cell host & microbe.* 25(3):377-388. [PubMed](#)
7. Wiede F, *et al.* 2021. *Cancer Discov.* Online ahead of print. [PubMed](#)
8. Li Z, *et al.* 2022. *Elife.* 11:. [PubMed](#)
9. Isakova-Sivak I, *et al.* 2020. *Vaccines (Basel).* 8:00. [PubMed](#)
10. Levine LS, *et al.* 2021. *Immunity.* 54(4):829-844.e5. [PubMed](#)
11. Perry JSA, *et al.* 2018. *Immunity.* 48:923. [PubMed](#)
12. Pfenninger P, *et al.* 2022. *Front Immunol.* 13:777113. [PubMed](#)
13. Mamedov MR, *et al.* 2018. *Immunity.* 48:350. [PubMed](#)
14. Lee CG, *et al.* 2019. *JCI Insight.* 4:5. [PubMed](#)
15. Yue X, *et al.* 2019. *Nat Commun.* 10:2011. [PubMed](#)
16. Kemp V, *et al.* 2018. *Cancer Gene Ther.* 26:268. [PubMed](#)
17. Goh PK, *et al.* 2022. *Sci Adv.* 8:eabk3338. [PubMed](#)
18. Gomez S, *et al.* 2022. *J Immunother Cancer.* 10:. [PubMed](#)
19. Jandke A, *et al.* 2020. *Nat Commun.* 3.075694444. [PubMed](#)
20. Yue X, *et al.* 2021. *EMBO Rep.* 22:e52716. [PubMed](#)
21. Zhang J, *et al.* 2022. *Biomater Res.* 26:44. [PubMed](#)
22. Lu X, *et al.* 2019. *Circ Res.* 125:1055. [PubMed](#)
23. Chen Y, *et al.* 2022. *Nat Commun.* 13:4468. [PubMed](#)
24. Damgaard RB *et al.* 2016. *Cell.* 166(5):1215-1230. [PubMed](#)
25. Muri J, *et al.* 2020. *Cell Reports.* 29(9):2731-2744.e4.. [PubMed](#)
26. Zhong C, *et al.* 2021. *J Virol.* 95:e0092521. [PubMed](#)
27. Boulch M, *et al.* 2021. *Sci Immunol.* 6:. [PubMed](#)
28. Culina S, *et al.* 2015. *Diabetes.* 64: 3532 - 3542. [PubMed](#)
29. Kong IY, *et al.* 2020. *Cell Rep.* 33:108290. [PubMed](#)

**RRID** AB\_10900082 (BioLegend Cat. No. 104435)  
AB\_2562560 (BioLegend Cat. No. 104436)

## Antigen Details

---

<b>Structure</b>	Selectin, 95 kD (neutrophils) or 74 kD (lymphocytes)
<b>Distribution</b>	Subsets of B and T cells, monocytes, granulocytes, subset of NK cells
<b>Function</b>	Lymphocyte homing to HEV, rolling on activated endothelium
<b>Ligand/Receptor</b>	CD34, GlyCAM-1, MAdCAM-1
<b>Cell Type</b>	B cells, Granulocytes, Monocytes, Neutrophils, NK cells, T cells, Tregs
<b>Biology Area</b>	Cell Adhesion, Cell Biology, Costimulatory Molecules, Immunology, Innate Immunity
<b>Molecular Family</b>	Adhesion Molecules, CD Molecules
<b>Antigen References</b>	1. Barclay AN, <i>et al.</i> 1997. <i>The Leukocyte Antigen FactsBook</i> Academic Press. 2. Kishimoto TK, <i>et al.</i> 1990. <i>P. Natl. Acad. Sci. USA</i> 87:2244. 3. Tedder TF, <i>et al.</i> 1995. <i>J. Exp. Med.</i> 181:2259.
<b>Gene ID</b>	<a href="#">20343</a>


## Related Protocols

[Cell Surface Flow Cytometry Staining Protocol](#)

## Other Formats

APC anti-mouse CD62L, Biotin anti-mouse CD62L, FITC anti-mouse CD62L, PE anti-mouse CD62L, PE/Cyanine5 anti-mouse CD62L, Purified anti-mouse CD62L, PE/Cyanine7 anti-mouse CD62L, Alexa Fluor® 488 anti-mouse CD62L, Alexa Fluor® 647 anti-mouse CD62L, Pacific Blue™ anti-mouse CD62L, Alexa Fluor® 700 anti-mouse CD62L, APC/Cyanine7 anti-mouse CD62L, PerCP/Cyanine5.5 anti-mouse CD62L, PerCP anti-mouse CD62L, Brilliant Violet 421™ anti-mouse CD62L, Brilliant Violet 570™ anti-mouse CD62L, Brilliant Violet 605™ anti-mouse CD62L, Brilliant Violet 510™ anti-mouse CD62L, Purified anti-mouse CD62L (Maxpar® Ready), Brilliant Violet 711™ anti-mouse CD62L, Brilliant Violet 785™ anti-mouse CD62L, PE/Dazzle™ 594 anti-mouse CD62L, APC/Fire™ 750 anti-mouse CD62L, TotalSeq™-A0112 anti-mouse CD62L, Brilliant Violet 650™ anti-mouse CD62L, TotalSeq™-C0112 anti-mouse CD62L, Ultra-LEAF™ Purified anti-mouse CD62L, KIRAVIA Blue 520™ anti-mouse CD62L, TotalSeq™-B0112 anti-mouse CD62L

## Product Data

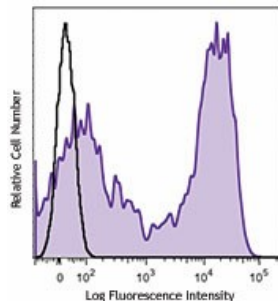

C57BL/6 mouse bone marrow cells were stained with CD62L (clone MEL-14) Brilliant Violet 421™ (filled histogram) or rat IgG2a, κ Brilliant Violet 421™ isotype control (open histogram). Data shown was gated on total cell population.

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

\*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, [www.biolegend.com/ordering#license](http://www.biolegend.com/ordering#license)). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587