

Brilliant Violet 510™ anti-human CD123 Antibody

| | |
|--------------------------|--|
| Catalog# / Size | 306021 / 25 tests 306022 / 100 tests |
| Clone | 6H6 |
| Regulatory Status | RUO |
| Other Names | IL-3R α , IL-3 Receptor alpha |
| Isotype | Mouse IgG1, κ |
| Description | CD123 is the 70 kD transmembrane α chain of the IL-3 receptor. Alone, CD123 binds IL-3 with low affinity; when CD123 associates with CDw131 (common β chain), it binds IL-3 with high affinity. CD123 does not transduce intracellular signals upon binding IL-3 and requires the β chain for this function. CD123 is expressed by myeloid precursors, macrophages, dendritic cells, mast cells, basophils, megakaryocytes, and some B cells. |

Product Details

| | |
|-------------------------------|--|
| Verified Reactivity | Human |
| Reported Reactivity | Rhesus |
| Antibody Type | Monoclonal |
| Host Species | Mouse |
| Immunogen | Human IL-3R α transfected COS cells. |
| Formulation | Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA). |
| Preparation | The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 510™ under optimal conditions. |
| Concentration | Lot-specific (to obtain lot-specific concentration, please enter the lot number in our Concentration and Expiration Lookup or Certificate of Analysis online tools.) |
| Storage & Handling | The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze. |
| Application | FC - Quality tested |
| Recommended Usage | <p>Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 μl per million cells in 100 μl staining volume or 5 μl per 100 μl of whole blood.</p> <p>Brilliant Violet 510™ excites at 405 nm and emits at 510 nm. The bandpass filter 510/50 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 510™ is a trademark of Sirigen Group Ltd.</p> <p>Learn more about Brilliant Violet™.</p> <p>This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.</p> |
| Excitation Laser | Violet Laser (405 nm) |
| Application Notes | Clone 6H6 does not inhibit IL-3 binding to low- or high-affinity IL-3Rs. Additional reported applications (for the relevant formats) include: Western blotting ¹ , immunoprecipitation ¹ , and immunohistochemical staining of acetone-fixed frozen sections ² and also paraformaldehyde fixed |

paraffin embedded tissue⁷.

Application References

(PubMed link indicates BioLegend citation)

1. Sun Q, *et al.* 1996. *Blood* 87:83. (IP, WB)
2. Herling M, *et al.* 2003. *Blood* 101:5007. (IHC)
3. Charles N, *et al.* 2010. *Nat. Med.* 16:701. (FC) [PubMed](#)
4. Martin-Gayo E, *et al.* 2010. *Blood* 115:5366. [PubMed](#)
5. Chen SC, *et al.* 2010. *Arch Dermatol Res.* 302:113. [PubMed](#)
6. Liu Y, *et al.* 2012. *Food Chem Toxicol.* 50:1920. [PubMed](#)
7. Peduzzi E, *et al.* 2007. *J. Invest. Dermatol.* 127:638. (IHC)

Product Citations

1. Carre C, *et al.* 2021. *iScience.* 24:102970. [PubMed](#)
2. Ols S, *et al.* 2020. *Cell Reports.* 30(12):3964-3971. [PubMed](#)
3. Romero-Olmedo AJ, *et al.* 2022. *Nat Microbiol.* 7:195. [PubMed](#)
4. Briceño O, *et al.* 2016. *PLoS One.* 11:e0166496. [PubMed](#)
5. Silver J, *et al.* 2020. *Methods Mol Biol.* 2121:115. [PubMed](#)
6. Grootens J *et al.* 2019. *EBioMedicine.* 43:150-158. [PubMed](#)
7. Lee J, *et al.* 2015. *J Exp Med.* 212:385. [PubMed](#)
8. Breton G, *et al.* 2015. *J Exp Med.* 212:401. [PubMed](#)
9. Abd Hamid M *et al.* 2019. *Cancer Immunol Res.* 7(8):1293-1306. [PubMed](#)

RRID

AB_2562067 (BioLegend Cat. No. 306021)
AB_2562068 (BioLegend Cat. No. 306022)

Antigen Details

| | |
|---------------------------|--|
| Structure | Ig superfamily, type I transmembrane glycoprotein, associates with CDw131, 70 kD |
| Distribution | Myeloid precursors, basophils, mast cells, macrophages, dendritic cells, megakaryocytes, subset of lymphocytes |
| Function | Hematopoietic cell proliferation, differentiation |
| Ligand/Receptor | IL-3 |
| Cell Type | Basophils, Dendritic cells, Hematopoietic stem and progenitors, Lymphocytes, Macrophages, Mast cells, Megakaryocytes |
| Biology Area | Immunology |
| Molecular Family | CD Molecules, Cytokine/Chemokine Receptors |
| Antigen References | 1. Miyajima A, <i>et al.</i> 1993. <i>Blood</i> 82:1960. |
| Gene ID | 3563 |



Related Protocols

[Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

Biotin anti-human CD123, PE anti-human CD123, Purified anti-human CD123, PE/Cyanine5 anti-human CD123, PE/Cyanine7 anti-human CD123, APC anti-human CD123, FITC anti-human CD123, PerCP/Cyanine5.5 anti-human CD123, Brilliant Violet 421™ anti-human CD123, Brilliant Violet 650™ anti-human CD123, Brilliant Violet 510™ anti-human CD123, Alexa Fluor® 647 anti-human CD123, Brilliant Violet 605™ anti-human CD123, Purified anti-human CD123 (Maxpar® Ready), Brilliant Violet 711™ anti-human CD123, Brilliant Violet 785™ anti-human CD123, PE/Dazzle™ 594 anti-human CD123, Alexa Fluor® 488 anti-human CD123, TotalSeq™-A0064 anti-human CD123, Alexa Fluor® 700 anti-human CD123, APC/Fire™ 750 anti-human CD123, Pacific Blue™ anti-human CD123, TotalSeq™-C0064 anti-human CD123, TotalSeq™-B0064 anti-human CD123, TotalSeq™-D0064 anti-human CD123, PerCP anti-human CD123, GMP PE/Cyanine7 anti-human CD123, APC/Fire™ 810 anti-human CD123 Antibody

Product Data

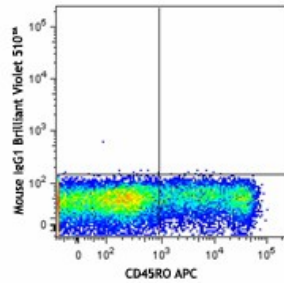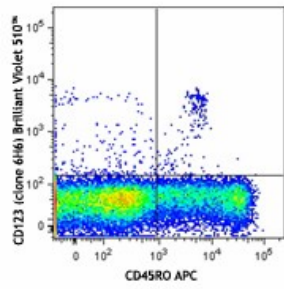

Human peripheral blood lymphocytes were stained with CD45RO APC and CD123 (clone 6H6) Brilliant Violet 510™ (top) or mouse IgG1 Brilliant Violet 510™ isotype control (bottom).

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
 Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587