

Brilliant Violet 510™ anti-human IL-17A Antibody

Catalog# / Size	512329 / 25 tests 512330 / 100 tests
Clone	BL168
Regulatory Status	RUO
Other Names	Interleukin IL-17A, IL-17, Cytotoxic T lymphocyte-associated antigen 8 (CTLA-8)
Isotype	Mouse IgG1, κ
Description	IL-17A is the founding member of the IL-17 family, a group of six structurally related pro-inflammatory cytokines. IL-17A, secreted by activated CD4 ⁺ Th17 cell subpopulation, elicits multiple biological activities on a variety of cells including: the induction of IL-6, IL-8, G-CSF, and PGE2 production in epithelial, endothelial or fibroblasts; the enhancement of surface expression of ICAM-1 in fibroblasts; activation of NF- κ B and costimulation of T cell proliferation. Recent studies demonstrated that, in mice, activated IL-17-secreting CD4 ⁺ helper T cells (Th17 cells) mediate an autoimmune arthritis that clinically and immunologically resembles rheumatoid arthritis (RA). Human IL-17A shows 63%, 63%, and 72% amino acid sequence identity to rat IL-17A, mouse IL-17A, and a protein encoded by the ORF13 gene of herpesvirus Saimiri (HVS), respectively.

Product Details

Verified Reactivity	Human
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	Recombinant full length human IL-17A
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA).
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 510™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration, please enter the lot number in our Concentration and Expiration Lookup or Certificate of Analysis online tools.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	ICFC - Quality tested
Recommended Usage	<p>Each lot of this antibody is quality control tested by intracellular immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 μl per million cells in 100 μl staining volume or 5 μl per 100 μl of whole blood.</p> <p>Brilliant Violet 510™ excites at 405 nm and emits at 510 nm. The bandpass filter 510/50 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 510™ is a trademark of Sirigen Group Ltd.</p> <p>Learn more about Brilliant Violet™.</p> <p>This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.</p>
Excitation Laser	Violet Laser (405 nm)

Product Citations 1. Wang J, *et al.* 2020. Cell. 183(7):1867-1883.e26. [PubMed](#)

RRID AB_2562744 (BioLegend Cat. No. 512329)
AB_2562745 (BioLegend Cat. No. 512330)

Antigen Details

Structure	Cytokine; dimer; ~15 kD (Mammalian)
Bioactivity	Secretion of IL-6, IL-8, G-CSF, prostaglandin E2 by epithelial, endothelial or fibroblastic cells; stimulates cell migration, cord formation, and IL-6 secretion by stromal cells
Cell Sources	Activated T cell subpopulation
Cell Targets	Fibroblasts, epithelial and endothelial cells, stromal cells
Receptors	IL-17R
Cell Type	Tregs
Biology Area	Cell Biology, Immunology, Neuroinflammation, Neuroscience
Molecular Family	Cytokines/Chemokines
Antigen References	1. Hirota K, <i>et al.</i> 2007. <i>J. Exp. Med.</i> 204:41. 2. Furuzawa-Carballeda J, <i>et al.</i> 2007. <i>Autoimmun. Rev.</i> 6:169. 3. Witowski J, <i>et al.</i> 2007. <i>Kidney Int.</i> 71:514. 4. Gaffen SL, <i>et al.</i> 2006. <i>Vitam. Horm.</i> 74:255. 5. Hymowitz S, <i>et al.</i> 2001. <i>EMBO J.</i> 20:5332.
Gene ID	3605

Related Protocols


[Surface and Intracellular Cytokine Staining for Flow Cytometry - Video](#)

[Intracellular Flow Cytometry Staining Protocol](#)


Other Formats

Brilliant Violet 421™ anti-human IL-17A, Purified anti-human IL-17A, FITC anti-human IL-17A, PE anti-human IL-17A, Alexa Fluor® 488 anti-human IL-17A, Alexa Fluor® 647 anti-human IL-17A, Pacific Blue™ anti-human IL-17A, PerCP/Cyanine5.5 anti-human IL-17A, PE/Cyanine7 anti-human IL-17A, Alexa Fluor® 700 anti-human IL-17A, APC/Cyanine7 anti-human IL-17A, Brilliant Violet 570™ anti-human IL-17A, Brilliant Violet 605™ anti-human IL-17A, Brilliant Violet 711™ anti-human IL-17A, Brilliant Violet 510™ anti-human IL-17A, APC anti-human IL-17A, Purified anti-human IL-17A (Maxpar® Ready), PE/Dazzle™ 594 anti-human IL-17A, Brilliant Violet 785™ anti-human IL-17A

Product Data

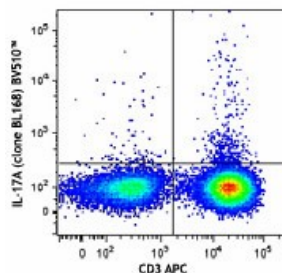

PMA+ ionomycin-stimulated human peripheral blood lymphocytes (in the presence of monensin) were stained with CD3 APC, fixed, permeabilized, and then stained with IL-17A (clone BL168) Brilliant Violet 510™ (top) or mouse IgG1, κ Brilliant Violet 510™ isotype control (bottom).

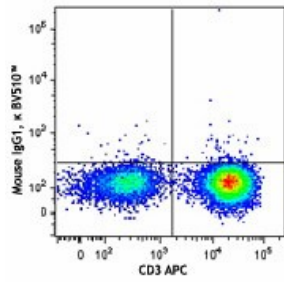

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587