

Brilliant Violet 650™ anti-human CD27 Antibody

Catalog# / Size	302827 / 25 tests 302828 / 100 tests
Clone	O323
Regulatory Status	RUO
Workshop	IV T-186
Other Names	S152, T14, TNFRSF7
Isotype	Mouse IgG1, κ
Description	CD27 is a 50-55 kD type I membrane protein also known as S152 and T14. It is a lymphocyte-specific member of the TNF-receptor superfamily. CD27 is expressed on medullary thymocytes, virtually all mature T cells, some B cells, and NK cells. CD27 binds to CD70 and plays an important role in costimulation of T cell activation, and regulation of B cell differentiation and proliferation. The cytoplasmic domains of CD27 have also been shown to interact with TRAF2 and TRAF5 to elicit NF- κ B and SAPK/JNK activation.

Product Details

Verified Reactivity	Human, Cynomolgus, Rhesus
Reported Reactivity	African Green, Baboon, Squirrel Monkey
Antibody Type	Monoclonal
Host Species	Mouse
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA).
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 650™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration, please enter the lot number in our Concentration and Expiration Lookup or Certificate of Analysis online tools.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	<p>Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 μl per million cells in 100 μl staining volume or 5 μl per 100 μl of whole blood.</p> <p>Brilliant Violet 650™ excites at 405 nm and emits at 645 nm. The bandpass filter 660/20 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 650™ is a trademark of Sirigen Group Ltd.</p> <p>Learn more about Brilliant Violet™.</p> <p>This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.</p>
Excitation Laser	Violet Laser (405 nm)
Application References	1. Knapp W, <i>et al.</i> Eds. 1989. Leucocyte Typing IV. Oxford University Press. New York.

(PubMed link indicates BioLegend citation)

Product Citations

2. Correia DV, *et al.* 2011. *Blood* 118:992. (FC) [PubMed](#)
1. Ye Y, *et al.* 2022. *Nat Commun.* 13:6458. [PubMed](#)
2. Dudley D, *et al.* 2016. *Nat Commun.* 7:12204. [PubMed](#)
3. Woodham A, *et al.* 2016. *Papillomavirus Res.* 2:21-30. [PubMed](#)
4. Gagne M, *et al.* 2022. *Cell.* 185:113. [PubMed](#)
5. Fajgenbaum DC, *et al.* 2019. *J Clin Invest.* 130:4451. [PubMed](#)
6. Kuri-Cervantes L, *et al.* 2020. *STAR Protoc.* 1:100154. [PubMed](#)
7. Szabo PA, *et al.* 2021. *Immunity.* 54(4):797-814.e6. [PubMed](#)
8. Aliota M, *et al.* 2016. *PLoS Negl Trop Dis.* 10:e0005168. [PubMed](#)
9. Jegaskanda S, *et al.* 2019. *JCI Insight.* 4:e129035. [PubMed](#)
10. Vanshylla K, *et al.* 2022. *Cell Host Microbe.* 30:69. [PubMed](#)
11. Ravimohan S, *et al.* 2016. *Clin Infect Dis.* 62: 795 - 803. [PubMed](#)
12. Fromentin R, *et al.* 2019. *Nat Commun.* 10:814. [PubMed](#)
13. Phad GE, *et al.* 2022. *Nat Immunol.* 23:1. [PubMed](#)
14. Japp AS, *et al.* 2021. *Cell.* 184(3):827-839.e14. [PubMed](#)
15. Kiniry B, *et al.* 2016. *Mucosal Immunol.* 10.1038/mi.2016.100. [PubMed](#)
16. McCracken M, *et al.* 2017. *PLoS Pathog.* 10.1371/journal.ppat.1006487. [PubMed](#)
17. Buggert M, *et al.* 2020. *Cell.* 183(7):1946-1961.e15. [PubMed](#)
18. Kiniry BE, *et al.* 2018. *J Immunol.* 200:1876. [PubMed](#)
19. Pinto D, *et al.* 2020. *Cell Host & Microbe.* 26(5):623-637.e8. [PubMed](#)
20. McBrien JB, *et al.* 2020. *Nature.* 578:154. [PubMed](#)
21. Fromentin R, *et al.* 2016. *PLoS Pathog.* 12: 1005761. [PubMed](#)
22. Pombo C, *et al.* 2015. *J Infect Dis.* 12: 1376-1386. [PubMed](#)
23. Burel JG, *et al.* 2019. *Elife.* 8. [PubMed](#)
24. Pardons M, *et al.* 2020. *Cell Reports.* 29(9):2783-2795.e5. [PubMed](#)

RRID

AB_11124941 (BioLegend Cat. No. 302827)
AB_2562096 (BioLegend Cat. No. 302828)

Antigen Details

Structure	TNF-R superfamily, type I transmembrane glycoprotein, 50-55 kD
Distribution	Medullary thymocytes, T and B cell subsets, NK cells
Function	T cell costimulation
Ligand/Receptor	CD70
Cell Type	B cells, NK cells, T cells, Thymocytes, Tregs
Biology Area	Costimulatory Molecules, Immunology
Molecular Family	CD Molecules
Antigen References	1. Hintzen R, <i>et al.</i> 1994. <i>Immunol. Today</i> 15:307. 2. Agematsu K, <i>et al.</i> 1995. <i>J. Immunol.</i> 154:3627.
Gene ID	939


Related Protocols

[Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

APC anti-human CD27, Biotin anti-human CD27, FITC anti-human CD27, PE anti-human CD27, Purified anti-human CD27, Alexa Fluor® 647 anti-human CD27, Alexa Fluor® 700 anti-human CD27, APC/Cyanine7 anti-human CD27, PerCP anti-human CD27, PerCP/Cyanine5.5 anti-human CD27, Pacific Blue™ anti-human CD27, Brilliant Violet 421™ anti-human CD27, Brilliant Violet 570™ anti-human CD27, Brilliant Violet 650™ anti-human CD27, Brilliant Violet 605™ anti-human CD27, Brilliant Violet 711™ anti-human CD27, Brilliant Violet 785™ anti-human CD27, Brilliant Violet 510™ anti-human CD27, PE/Cyanine7 anti-human CD27, Purified anti-human CD27 (Maxpar® Ready), PE/Dazzle™ 594 anti-human CD27, APC/Fire™ 750 anti-human CD27, TotalSeq™-A0154 anti-human CD27, Brilliant Violet 750™ anti-human CD27, TotalSeq™-B0154 anti-human CD27, TotalSeq™-C0154 anti-human CD27, Spark NIR™ 685 anti-human CD27, PE/Fire™ 810 anti-human CD27, TotalSeq™-D0154 anti-human CD27, APC/Fire™ 810 anti-human CD27, PE/Cyanine5 anti-human CD27 Antibody, Spark UV™ 387 anti-human CD27

Product Data

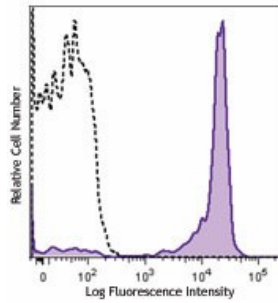

Human peripheral blood lymphocytes were stained with CD27 (clone O323) Brilliant Violet 650™.

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587