

Brilliant Violet 750™ anti-human CD56 (NCAM) Antibody

Catalog# / Size	362555 / 25 tests 362556 / 100 tests
Clone	5.1H11
Regulatory Status	RUO
Other Names	Leu-19, NKH1, NCAM-1
Isotype	Mouse IgG1, κ
Description	CD56 is a single transmembrane glycoprotein also known as NCAM (neural cell adhesion molecule), Leu-19, or NKH1. It is a member of the Ig superfamily. The 140 kD isoform is expressed on NK and NKT cells. CD56 is also expressed in the brain (cerebellum and cortex) and at neuromuscular junctions. Certain large granular lymphocyte (LGL) leukemias, small-cell lung carcinomas, neuronal derived tumors, myelomas, and myeloid leukemias also express CD56. CD56 plays a role in homophilic and heterophilic adhesion via binding to itself or heparan sulfate.

Product Details

Verified Reactivity	Human
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	Human myotube cells
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA).
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 750™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration, please enter the lot number in our Concentration and Expiration Lookup or Certificate of Analysis online tools.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 µl per million cells in 100 µl staining volume or 5 µl per 100 µl of whole blood.

Brilliant Violet 750™ excites at 405 nm and emits at 750 nm. The bandpass filter 780/60 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 750™ is a trademark of Sirigen Group Ltd.

[Learn more about Brilliant Violet™.](#)

This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.

Application References

(PubMed link indicates BioLegend citation)

1. Walsh FS, *et al.* 1981. *Nature* 289:60. (FC)
2. Pavlath GK, *et al.* 1986. *J. Cell Biol.* 102:124. (FC)
3. Pavlath GK, *et al.* 1989. *Nature* 337:570. (FC)
4. Pulido R, *et al.* 1988. *J. Immunol.* 140:3851. (FC)

Product Citations

1. Cohen KW, *et al.* 2021. Cell Reports Medicine. 2(7):100354. [PubMed](#)
2. Obradovic A, *et al.* 2021. Cell. 184(11):2988-3005.e16. [PubMed](#)
3. Qing G, *et al.* 2021. Front Cell Dev Biol. 9:761300. [PubMed](#)

RRID

AB_2734396 (BioLegend Cat. No. 362555)
AB_2801001 (BioLegend Cat. No. 362556)


Antigen Details

Structure	Ig superfamily, single transmembrane or GPI-anchored glycoprotein
Distribution	NK cells, T subset, neural tissue, some LGL and myeloid leukemias
Function	Adhesion
Ligand/Receptor	Heparan sulfate
Antigen References	1. Lanier L, <i>et al.</i> 1991. <i>J. Immunol.</i> 146:4421 2. Hemperly J, <i>et al.</i> 1990. <i>J. Mol. Neurosci.</i> 2:71 3. Cremer H, <i>et al.</i> 1994. <i>Nature</i> 367:455.
Gene ID	4684

Other Formats

Purified anti-human CD56 (NCAM), APC anti-human CD56 (NCAM), PerCP/Cyanine5.5 anti-human CD56 (NCAM), PE anti-human CD56 (NCAM), PE/Cyanine7 anti-human CD56 (NCAM), APC/Cyanine7 anti-human CD56 (NCAM), PE/Cyanine5 anti-human CD56 (NCAM), Alexa Fluor® 647 anti-human CD56 (NCAM), Alexa Fluor® 488 anti-human CD56 (NCAM), Pacific Blue™ anti-human CD56 (NCAM), Alexa Fluor® 700 anti-human CD56 (NCAM), PerCP anti-human CD56 (NCAM), Brilliant Violet 650™ anti-human CD56 (NCAM), Brilliant Violet 510™ anti-human CD56 (NCAM), Brilliant Violet 421™ anti-human CD56 (NCAM), Biotin anti-human CD56 (NCAM), Brilliant Violet 605™ anti-human CD56 (NCAM), Brilliant Violet 570™ anti-human CD56 (NCAM), Brilliant Violet 711™ anti-human CD56 (NCAM), PE/Dazzle™ 594 anti-human CD56 (NCAM), FITC anti-human CD56 (NCAM), Brilliant Violet 785™ anti-human CD56 (NCAM), Ultra-LEAF™ Purified anti-human CD56 (NCAM), APC/Fire™ 750 anti-human CD56 (NCAM), Brilliant Violet 750™ anti-human CD56 (NCAM), TotalSeq™-A0047 anti-human CD56 (NCAM), TotalSeq™-C0047 anti-human CD56 (NCAM), TotalSeq™-B0047 anti-human CD56 (NCAM), Spark NIR™ 685 anti-human CD56 (NCAM), KIRAVIA Blue 520™ anti-human CD56 (NCAM), GMP PE anti-human CD56 (NCAM), TotalSeq™-D0047 anti-human CD56 (NCAM), GMP APC anti-human CD56 (NCAM)

Product Data

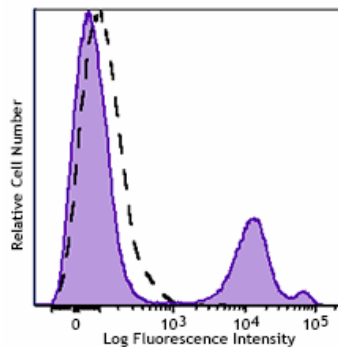

Human peripheral blood lymphocytes stained with CD56 (clone 5.1H11) Brilliant Violet 750™ (filled histogram) or mouse IgG1, κ Brilliant Violet 750™ isotype control (open histogram).

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587

