

## Purified anti-mouse CD69 Antibody

<b>Catalog# / Size</b>	104502 / 500 µg
<b>Clone</b>	H1.2F3
<b>Regulatory Status</b>	RUO
<b>Other Names</b>	Very Early Activation Antigen (VEA), AIM, EA1, MLR3, gp34/28
<b>Isotype</b>	Armenian Hamster IgG
<b>Description</b>	CD69 is a 60 kD type II membrane protein composed of a 27/33 kD disulfide-linked homodimer, also known as Very Early Activation Antigen (VEA), AIM, EA1, MLR3, and gp34/28. It is expressed on a subset of thymocytes and platelets. CD69 is rapidly induced on activated T and B cells, neutrophils, and NK cells. It is a C-type lectin, closely related to the NKR-P1 and Ly-49 NK cell activation molecules. CD69 is involved in the early events of cell activation and thymocyte positive selection.

### Product Details

---

<b>Verified Reactivity</b>	Mouse
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Armenian Hamster
<b>Immunogen</b>	Mouse dendritic epidermal T cell line Y245
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide.
<b>Preparation</b>	The antibody was purified by affinity chromatography.
<b>Concentration</b>	0.5 mg/ml
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C.
<b>Application</b>	<a href="#">FC - Quality tested</a> <a href="#">IHC-F, IP - Reported in the literature, not verified in house</a>
<b>Recommended Usage</b>	Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a> . For flow cytometric staining, the suggested use of this reagent is ≤ 1.0 µg per million cells in 100 µl volume. It is recommended that the reagent be titrated for optimal performance for each application.
<b>Application Notes</b>	The H1.2F3 antibody has been reported to augment T cell activation. Additional reported applications (for the relevant formats) include: <i>in vitro</i> T cell and NK cell activation <sup>1-3</sup> , immunohistochemistry <sup>4,5</sup> , and immunoprecipitation <sup>1</sup> .  This antibody has been characterized in the literature as containing a lambda (?) light chain.

### Application References

(PubMed link indicates BioLegend citation)

1. Yokoyama WM, *et al.* 1988. *J. Immunol.* 141:369. (IP)
2. Sobel ES, *et al.* 1993. *J. Immunol.* 150:673.
3. Karlhofer FM, *et al.* 1991. *J. Immunol.* 146:3662.
4. Zhou X, *et al.* 2005. *J. Biol. Chem.* 280:31240. (IHC)
5. Podd BS, *et al.* 2006. *J. Immunol.* 176:6532. (IHC)
6. Lawson BR, *et al.* 2007. *J. Immunol.* 178:5366.
7. Lee JW, *et al.* 2006. *Nature Immunol.* 8:181.
8. Epardaud M, *et al.* 2008. *Cancer Res.* 15:2972. [PubMed](#)
9. Jordan JM, *et al.* 2008. 76:3717. [PubMed](#)
10. Kenna TJ, *et al.* 2008. *Blood* 111:2091. [PubMed](#)
11. Ishikawa C, *et al.* 2013. *Biochim Biophys Acta.* 167:99. [PubMed](#)

### Product Citations

1. Arnold IC, *et al.* 2019. *PLoS Pathog.* 15:e1007866. [PubMed](#)
2. Guo R, *et al.* 2020. *Cell Res.* 30:21. [PubMed](#)
3. Khan KA, *et al.* 2020. *NPJ Breast Cancer.* 6:29. [PubMed](#)

4. Krovi SH, *et al.* 2019. Proc Natl Acad Sci U S A. 116:22252. [PubMed](#)
5. Burrello C, *et al.* 2018. Front Med (Lausanne). 5:21. [PubMed](#)
6. Chartrand K, *et al.* 2018. Front Immunol. 1.642361111. [PubMed](#)
7. Xu L, *et al.* 2017. Int J Biol Macromol. . 10.1016/j.ijbiomac.2017.07.090. [PubMed](#)
8. Rieck M, *et al.* 2017. Eur J Immunol. 47:677. [PubMed](#)
9. Wang W, *et al.* 2018. Cancer Cell. 34:757. [PubMed](#)
10. Oda SK, *et al.* 2020. J Exp Med. 217: . [PubMed](#)
11. Carpenter SM, *et al.* 2017. PLoS Pathog. 13:e1006704. [PubMed](#)
12. Niss Arfelt K, *et al.* 2017. Blood. 129:866. [PubMed](#)
13. Kroczek AL, *et al.* 2019. Front Immunol. 9:2806. [PubMed](#)
14. Pawaria S, *et al.* 2020. Arthritis Rheumatol. 72:359. [PubMed](#)
15. Morabito KM, *et al.* 2017. Mucosal Immunol. 0.795138889. [PubMed](#)
16. Chen L, *et al.* 2020. Front Immunol. 11:584458. [PubMed](#)
17. Kiripolsky J, *et al.* 2021. J Autoimmun. 118:102608. [PubMed](#)
18. Khan KA, *et al.* 2020. NPJ Breast Cancer. 6:29. [PubMed](#)
19. Seelige R, *et al.* 2018. Sci Rep. 8:13670. [PubMed](#)
20. Snell LM, *et al.* 2018. Immunity. 49:678. [PubMed](#)
21. Xie D, *et al.* 2020. Eur J Immunol. 50:1729. [PubMed](#)
22. Liu H, *et al.* 2020. J Immunol. 205:1207. [PubMed](#)
23. Abboud G, *et al.* 2018. Front Immunol. 9:1973. [PubMed](#)
24. Hsiao CC, *et al.* 2021. Cells. 10: . [PubMed](#)
25. Tsubaki T, *et al.* 2018. Oncotarget. 9:11209. [PubMed](#)
26. Hoves S, *et al.* 2018. J Exp Med. 215:859. [PubMed](#)
27. Sun Y, *et al.* 2020. J Immunol. 205:3358. [PubMed](#)
28. Schuler CF, *et al.* 2020. Allergy. 75:2279. [PubMed](#)
29. Ishikawa C, *et al.* 2013. Biochim Biophys Acta. 167:99. [PubMed](#)
30. Park GY, *et al.* 2020. Immune Netw. e43:20. [PubMed](#)
31. Lee JY, *et al.* 2018. Front Immunol. 0.678472222. [PubMed](#)
32. Kenna T, *et al.* 2008. Blood. 111:2091. [PubMed](#)
33. Clark K, *et al.* 2019. Front Immunol. 10:2355. [PubMed](#)
34. Karmaus P, *et al.* 2011. J Leukoc Biol. 90:983. [PubMed](#)
35. Brownlie D, *et al.* 2021. J Immunother Cancer. 9: . [PubMed](#)
36. Tsai MS, *et al.* 2021. Int J Mol Sci. 22: . [PubMed](#)
37. Roy S, *et al.* 2018. Cancer Res. 78:5600. [PubMed](#)
38. Wagner AK, *et al.* 2017. Nat Commun. 8:15627. [PubMed](#)
39. Cui X, *et al.* 2017. J Immunol. 199:4066. [PubMed](#)
40. Scutts SR, *et al.* 2018. Cell Rep. 25:1953. [PubMed](#)
41. Cané S, *et al.* 2021. J Immunother Cancer. 9: . [PubMed](#)
42. Xu W, *et al.* 2021. Immunity. 54(3):526-541.e7. [PubMed](#)
43. Sivapatham S, *et al.* 2019. Front Immunol. 1.865277778. [PubMed](#)
44. Richards KA, *et al.* 2018. Front Immunol. 0.829861111. [PubMed](#)
45. Englezou PC, *et al.* 2018. Mol Ther Nucleic Acids. 12:118. [PubMed](#)
46. Diao J, *et al.* 2018. J Immunol. 201:1306. [PubMed](#)
47. Rios-Doria J, *et al.* 2017. Cancer Res. 77:2686. [PubMed](#)
48. Jordan J, *et al.* 2008. Infect Immun. 76:3717. [PubMed](#)
49. Kästele V, *et al.* 2021. Mucosal Immunol. 14:717. [PubMed](#)
50. Sutiwisesak R, *et al.* 2020. PLoS Pathog. 16:e1009000. [PubMed](#)
51. Epardaud M, *et al.* 2008. Cancer Res. 68:2972. [PubMed](#)

**RRID** AB\_313105 (BioLegend Cat. No. 104502)

## Antigen Details

---

<b>Structure</b>	C-type lectin, 27/33 kD
<b>Distribution</b>	Activated T cells and B cells, NK cells, granulocytes, thymocytes, platelets
<b>Function</b>	Lymphocyte activation
<b>Cell Type</b>	B cells, Granulocytes, NK cells, Platelets, T cells, Thymocytes, Tregs
<b>Biology Area</b>	Costimulatory Molecules, Immunology, Innate Immunity
<b>Molecular Family</b>	CD Molecules
<b>Antigen References</b>	<ol style="list-style-type: none"> <li>1. Barclay AN, <i>et al.</i> 1997. The Leukocyte Antigen FactsBook Academic Press.</li> <li>2. Testi R, <i>et al.</i> 1994. <i>Immunol. Today</i> 15:479.</li> <li>3. Moretta A, <i>et al.</i> 1991. <i>J. Exp. Med.</i> 174:1393.</li> <li>4. Yokoyama WM, <i>et al.</i> 1988. <i>J. Immunol.</i> 141:369.</li> </ol>
<b>Gene ID</b>	<a href="#">12515</a>


## Related Protocols

---

## Other Formats

Biotin anti-mouse CD69, FITC anti-mouse CD69, PE anti-mouse CD69, PE/Cyanine5 anti-mouse CD69, Purified anti-mouse CD69, PE/Cyanine7 anti-mouse CD69, APC anti-mouse CD69, Alexa Fluor® 488 anti-mouse CD69, Alexa Fluor® 647 anti-mouse CD69, PerCP anti-mouse CD69, PerCP/Cyanine5.5 anti-mouse CD69, Pacific Blue™ anti-mouse CD69, Brilliant Violet 421™ anti-mouse CD69, APC/Cyanine7 anti-mouse CD69, Brilliant Violet 605™ anti-mouse CD69, Brilliant Violet 510™ anti-mouse CD69, Purified anti-mouse CD69 (Maxpar® Ready), PE/Dazzle™ 594 anti-mouse CD69, Brilliant Violet 711™ anti-mouse CD69, Alexa Fluor® 700 anti-mouse CD69, Brilliant Violet 650™ anti-mouse CD69, Brilliant Violet 785™ anti-mouse CD69, TotalSeq™-A0197 anti-mouse CD69, APC/Fire™ 750 anti-mouse CD69, TotalSeq™-C0197 anti-mouse CD69, TotalSeq™-B0197 anti-mouse CD69, KIRAVIA Blue 520™ anti-mouse CD69, Spark NIR™ 685 anti-mouse CD69, Spark Red™ 718 anti-mouse CD69

## Product Data

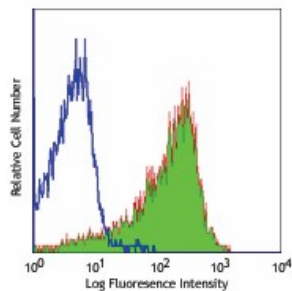

Con A-stimulated Balb/c mouse splenocytes (24 hours) stained with purified H1.2F3, followed by anti-Armenian hamster IgG FITC

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

\*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, [www.biolegend.com/ordering#license](http://www.biolegend.com/ordering#license)). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587