

## Purified anti-mouse F4/80 (Maxpar<sup>®</sup> Ready) Antibody

<b>Catalog# / Size</b>	123143 / 100 µg
<b>Clone</b>	BM8
<b>Regulatory Status</b>	RUO
<b>Other Names</b>	EMR1, Ly71
<b>Isotype</b>	Rat IgG2a, κ
<b>Description</b>	F4/80 is a 160 kD glycoprotein. It is characterized as a member of the epidermal growth factor (EGF)-transmembrane 7 (TM7) family. F4/80, also known as EMR1 or Ly71, has been widely used as a murine macrophage marker, which is expressed on the majority of tissue macrophages including peritoneal macrophages, macrophages in lung, gut, thymus and red pulp of spleen (but not on the macrophages located in T cell areas of the spleen, lymph node and Peyer's patch), Kuffer cells, Langerhans cells, and bone marrow stromal cells. F4/80 has also been shown on a subset of dendritic cells. The biological ligand of F4/80 has not been identified, but it has been reported that F4/80 is required for induction of CD8 <sup>+</sup> T cells-mediated peripheral tolerance.

### Product Details

<b>Verified Reactivity</b>	Mouse
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Rat
<b>Immunogen</b>	Murine macrophages
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and EDTA.
<b>Preparation</b>	The antibody was purified by affinity chromatography.
<b>Concentration</b>	1.0 mg/ml
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C.
<b>Application</b>	<a href="#">FC - Quality tested</a> <a href="#">CyTOF<sup>®</sup> - Verified</a>
<b>Recommended Usage</b>	This product is suitable for use with the <a href="#">Maxpar<sup>®</sup> Metal Labeling Kits</a> . For metal labeling using Maxpar <sup>®</sup> Ready antibodies, proceed directly to the step to Partially Reduce the Antibody by adding 100 µl of Maxpar <sup>®</sup> Ready antibody to 100 µl of 4 mM TCEP-R in a 50 kDa filter and continue with the protocol. Always refer to the latest version of Maxpar <sup>®</sup> User Guide when conjugating Maxpar <sup>®</sup> Ready antibodies.
<b>Application Notes</b>	Additional reported applications (for the relevant formats) include: immunohistochemical staining of acetone-fixed frozen sections <sup>1,2</sup> and formalin-fixed paraffin-embedded sections <sup>6,7</sup> , Western blotting, and spatial biology (IBEX) <sup>12,13</sup> .
<b>Additional Product Notes</b>	Maxpar <sup>®</sup> is a registered trademark of Standard BioTools Inc.
<b>Application References</b>	<ol style="list-style-type: none"> <li>Schaller E, <i>et al.</i> 2002. <i>Mol. Cell. Biol.</i> 22:8035. (IHC)</li> <li>Stevceva L, <i>et al.</i> 2001. <i>BMC Clin Pathol.</i> 1:3. (IHC)</li> <li>Kobayashi M, <i>et al.</i> 2008. <i>J. Leukoc. Biol.</i> 83:1354. <a href="#">PubMed</a></li> <li>Poeckel D, <i>et al.</i> 2009. <i>J. Biol Chem.</i> 284:21077. <a href="#">PubMed</a></li> <li>Glass AM, <i>et al.</i> 2013. <i>J. Immunol.</i> 190:4830. <a href="#">PubMed</a></li> <li>Koehm S, <i>et al.</i> 2007. <i>J. Allergy Clin. Immunol.</i> 120:570. (IHC)</li> <li>Rankin AL, <i>et al.</i> 2010. <i>J. Immunol.</i> 184:1526. (IHC)</li> <li>Sasi SP, <i>et al.</i> 2014. <i>J Biol Chem.</i> 289:14178. <a href="#">PubMed</a></li> <li>Thakus VS, <i>et al.</i> 2014. <i>Toxicol Lett.</i> 230:322. <a href="#">PubMed</a></li> <li>Watson NB, <i>et al.</i> 2015. <i>J Immunol.</i> 194:2796. <a href="#">PubMed</a></li> <li>Hirakawa H, <i>et al.</i> 2015. <i>PLoS One.</i> 10:119360. <a href="#">PubMed</a></li> <li>Radtke AJ, <i>et al.</i> 2020. <i>Proc Natl Acad Sci U S A.</i> 117:33455-65. (SB) <a href="#">PubMed</a></li> </ol>
<b>(PubMed link indicates BioLegend citation)</b>	

13. Radtke AJ, et al. 2022. *Nat Protoc.* 17:378-401. (SB) [PubMed](#)

## Product Citations

1. McDonald B, et al. 2020. *Cell Host Microbe.* 28(5):660-668.e4. [PubMed](#)
2. Jordan S, et al. 2020. *Cell.* 178(5):1102-1114.e17.. [PubMed](#)

## RRID

AB\_2563767 (BioLegend Cat. No. 123143)

## Antigen Details

<b>Structure</b>	EGF-TM7 family member, 160 kD glycoprotein
<b>Distribution</b>	Majority of tissue macrophages including peritoneal macrophages, macrophages in lung, gut, thymus and red pulp of spleen, Kuffer cells, Langerhans cells, bone marrow stromal cells, and a subset of dendritic cells
<b>Function</b>	Induction of immunological tolerance
<b>Cell Type</b>	Dendritic cells, Langerhans cells, Macrophages, Tregs
<b>Biology Area</b>	Cell Biology, Immunology, Innate Immunity, Neuroinflammation, Neuroscience
<b>Antigen References</b>	<ol style="list-style-type: none"><li>1. Austy JM and Gordon S. 1981. <i>Eur. J. Immunol.</i> 11:805.</li><li>2. Hume DA, et al. 1983. <i>J. Exp. Med.</i> 158:1522.</li><li>3. Ruedl C, et al. 1996. <i>Eur. J. Immunol.</i> 26:1801.</li><li>4. McKnight AJ, et al. 1996. <i>J. Biol. Chem.</i> 271:486.</li><li>5. Lin HH, et al. 2005. <i>J. Exp. Med.</i> 201:1615.</li></ol>
<b>Gene ID</b>	<a href="#">13733</a>


## Related Protocols

[Cell Surface Flow Cytometry Staining Protocol](#)

## Other Formats

Brilliant Violet 605™ anti-mouse F4/80, Purified anti-mouse F4/80, Biotin anti-mouse F4/80, FITC anti-mouse F4/80, PE anti-mouse F4/80, PE/Cyanine5 anti-mouse F4/80, PE/Cyanine7 anti-mouse F4/80, APC anti-mouse F4/80, APC/Cyanine7 anti-mouse F4/80, Alexa Fluor® 488 anti-mouse F4/80, Alexa Fluor® 647 anti-mouse F4/80, Pacific Blue™ anti-mouse F4/80, PerCP anti-mouse F4/80, PerCP/Cyanine5.5 anti-mouse F4/80, Alexa Fluor® 700 anti-mouse F4/80, Brilliant Violet 421™ anti-mouse F4/80, Brilliant Violet 510™ anti-mouse F4/80, Alexa Fluor® 594 anti-mouse F4/80, Brilliant Violet 785™ anti-mouse F4/80, Purified anti-mouse F4/80 (Maxpar® Ready), PE/Dazzle™ 594 anti-mouse F4/80, Brilliant Violet 650™ anti-mouse F4/80, Brilliant Violet 711™ anti-mouse F4/80, APC/Fire™ 750 anti-mouse F4/80, TotalSeq™-A0114 anti-mouse F4/80, TotalSeq™-B0114 anti-mouse F4/80, TotalSeq™-C0114 anti-mouse F4/80, Spark YG™ 570 anti-mouse F4/80, KIRAVIA Blue 520™ anti-mouse F4/80, Ultra-LEAF™ Purified anti-mouse F4/80, APC/Fire™ 810 anti-mouse F4/80, Spark NIR™ 685 anti-mouse F4/80

## Product Data

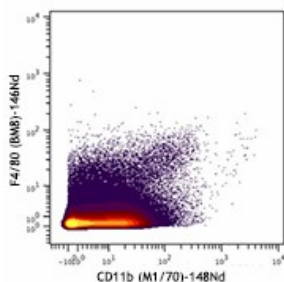

Mouse bone-marrow stained with 148Nd-anti-CD11b (M1/70) and 146Nd-anti-F4/80 (BM8). Data provided by DVS Sciences.

use of our products.

\*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, [www.biolegend.com/ordering#license](http://www.biolegend.com/ordering#license)). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587