

Alexa Fluor[®] 594 anti-mouse/human CD45R/B220 Antibody

Catalog# / Size	103254 / 100 µg
Clone	RA3-6B2
Regulatory Status	RUO
Other Names	B220
Isotype	Rat IgG2a, κ

Description CD45R, also known as B220, is an isoform of CD45. It is a member of the protein tyrosine phosphatase (PTP) family with a molecular weight of approximately 180-240 kD. CD45R is expressed on B cells (at all developmental stages from pro-B cells through mature B cells), activated B cells, and subsets of T and NK cells. CD45R (B220) is also expressed on a subset of abnormal T cells involved in the pathogenesis of systemic autoimmunity in MRL-*Fas*^{lpr} and MRL-*Fas*^{gld} mice. It plays a critical role in TCR and BCR signaling. The primary ligands for CD45 are galectin-1, CD2, CD3, and CD4. CD45R is commonly used as a pan-B cell marker; however, CD19 may be more appropriate for B cell specificity.

Product Details

Verified Reactivity	Mouse, Human
Reported Reactivity	Cat
Antibody Type	Monoclonal
Host Species	Rat
Immunogen	Abelson murine leukemia virus-induced pre-B tumor cells
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide.
Preparation	The antibody was purified by affinity chromatography and conjugated with Alexa Fluor [®] 594 under optimal conditions.
Concentration	0.5 mg/mL
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	IHC-F - Quality tested FC, 3D IHC - Verified
Recommended Usage	<p>Each lot of this antibody is quality control tested by immunofluorescence staining. For immunohistochemistry, a concentration range of 2.5-5 µg/mL is recommended. For immunofluorescence microscopy, a concentration range of 1.25 - 10 µg/mL is recommended. For flow cytometric staining, the suggested use of this reagent is ≤0.25 µg per million cells in 100 µL volume. For 3D immunohistochemistry on formalin-fixed tissues, a concentration of 5.0 µg/mL is suggested. It is recommended that the reagent be titrated for optimal performance for each application.</p> <p>* Alexa Fluor[®] 594 has an excitation maximum of 590 nm, and a maximum emission of 617 nm.</p> <p>Alexa Fluor[®] and Pacific Blue™ are trademarks of Life Technologies Corporation.</p> <p>View full statement regarding label licenses</p>
Application Notes	Clone RA3-6B2 has been described to react with an epitope on the extracellular domain of the transmembrane CD45 glycoprotein which is dependent upon the expression of exon A and specific carbohydrate residues. Additional reported applications (for the relevant formats) include: immunoprecipitation ¹ , <i>in vitro</i> and <i>in vivo</i> modulation of B cell responses ^{2,4} , immunohistochemistry of acetone-fixed frozen sections and formalin-fixed paraffin-embedded sections ^{5,6} , and spatial biology (IBEX) ^{14,15} .

Application References

**(PubMed link indicates
BioLegend citation)**

1. Coffman RL. 1982. *Immunol. Rev.* 69:5. (IP)
2. George A, *et al.* 1994. *J. Immunol.* 152:1014. (Activ)
3. Asensi V, *et al.* 1989. *Immunology* 68:204. (Activ)
4. Domiati-Saad R, *et al.* 1993. *J. Immunol.* 151:5936. (Activ)
5. Hata H, *et al.* 2004. *J. Clin. Invest.* 114:582. (IHC)
6. Monteith CE, *et al.* 1996. *Can. J. Vet. Res.* 60:193. (IHC)
7. Shih FF, *et al.* 2006. *J. Immunol.* 176:3438. (FC)
8. Chang C L-T, *et al.* 2007. *J. Immunol.* 178:6984.
9. Fazilleau N, *et al.* 2007. *Nature Immunol.* 8:753.
10. Lang GL, *et al.* 2008. *Blood* 111:2158. [PubMed](#)
11. Charles N, *et al.* 2010. *Nat. Med.* 16:701. (FC) [PubMed](#)
12. del Rio ML, *et al.* 2011. *Transpl. Int.* 24:501. (FC) [PubMed](#)
13. Murakami R, *et al.* 2013. *PLoS One.* 8:73270. [PubMed](#)
14. Radtke AJ, *et al.* 2020. *Proc Natl Acad Sci U S A.* 117:33455-65. (SB) [PubMed](#)
15. Radtke AJ, *et al.* 2022. *Nat Protoc.* 17:378-401. (SB) [PubMed](#)

Product Citations

1. Baptista AP *et al.* 2019. *Immunity.* 50(5):1188-1201 . [PubMed](#)
2. Li S, *et al.* 2020. *Cell Death Dis.* 0.9875. [PubMed](#)
3. Pham THM, *et al.* 2020. *Cell Host & Microbe.* 27(1):54-67.e5.. [PubMed](#)
4. Counoupas C, *et al.* 2020. *NPJ Vaccines.* 0.28125. [PubMed](#)
5. Matryba P, *et al.* 2020. *J Immunol.* 1395:204. [PubMed](#)
6. Wang X, *et al.* 2021. *EMBO J.* 40:e105926. [PubMed](#)
7. Kobia FM, *et al.* 2020. *PLoS Biol.* 18:e3000850. [PubMed](#)
8. Luo L, *et al.* 2021. *J Neuroinflammation.* 18:27. [PubMed](#)
9. Shi H *et al.* 2018. *Immunity.* 49(5):899-914 . [PubMed](#)
10. Song W, *et al.* 2022. *Immunity.* 55:290. [PubMed](#)
11. Goldberg MF *et al.* 2018. *Immunity.* 49(6):1090-1102 . [PubMed](#)
12. Kuhn NF *et al.* 2019. *Cancer cell.* 35(3):473-488 . [PubMed](#)
13. Bidet K, *et al.* 2019. *NPJ Vaccines.* 4:27. [PubMed](#)

RRID

AB_2563229 (BioLegend Cat. No. 103254)

Antigen Details

Structure	Protein tyrosine phosphatase (PTP) family, 180-240 kD
Distribution	B cells, T cell subset, NK cell subset
Function	Phosphatase, T and B cell activation
Ligand/Receptor	Galectin-1, CD2, CD3, CD4
Cell Type	B cells, NK cells, T cells
Biology Area	Cell Biology, Immunology, Inhibitory Molecules, Neuroscience, Neuroscience Cell Markers
Molecular Family	CD Molecules
Antigen References	<ol style="list-style-type: none">1. Barclay A, <i>et al.</i> 1997. <i>The Leukocyte Antigen FactsBook</i> Academic Press.2. Trowbridge IS, <i>et al.</i> 1993. <i>Annu. Rev. Immunol.</i> 12:85.3. Kishihara K, <i>et al.</i> 1993. <i>Cell</i> 74:143.4. Pulido R, <i>et al.</i> 1988. <i>J. Immunol.</i> 140:3851.
Gene ID	19264 5788

Related Protocols

[Immunohistochemistry Protocol for Frozen Sections](#)

[Cell Surface Flow Cytometry Staining Protocol](#)


[Immunocytochemistry Staining Protocol](#)

[Ce3D™ Tissue Clearing Kit](#)


Other Formats

Alexa Fluor® 594 anti-mouse/human CD45R/B220, APC anti-mouse/human CD45R/B220, Biotin anti-mouse/human CD45R/B220, FITC anti-mouse/human CD45R/B220, PE anti-mouse/human CD45R/B220, PE/Cyanine5 anti-mouse/human CD45R/B220, Purified anti-mouse/human CD45R/B220, PE/Cyanine7 anti-mouse/human CD45R/B220, APC/Cyanine7 anti-mouse/human CD45R/B220, Alexa Fluor® 488 anti-mouse/human CD45R/B220, Alexa Fluor® 647 anti-mouse/human CD45R/B220, Pacific Blue™ anti-mouse/human CD45R/B220, Alexa Fluor® 700 anti-mouse/human CD45R/B220, PerCP anti-mouse/human CD45R/B220, PerCP/Cyanine5.5 anti-mouse/human CD45R/B220, Brilliant Violet 421™ anti-mouse/human CD45R/B220, Brilliant Violet 570™ anti-mouse/human CD45R/B220, Brilliant Violet 650™ anti-mouse/human CD45R/B220, Brilliant Violet 605™ anti-mouse/human CD45R/B220, Brilliant Violet 785™ anti-mouse/human CD45R/B220, Brilliant Violet 510™ anti-mouse/human CD45R/B220, Purified anti-mouse/human CD45R/B220 (Maxpar® Ready), Brilliant Violet 711™ anti-mouse/human CD45R/B220, PE/Dazzle™ 594 anti-mouse/human CD45R/B220, APC/Fire™ 750 anti-mouse/human CD45R/B220, Brilliant Violet 750™ anti-mouse/human CD45R/B220, TotalSeq™-A0103 anti-mouse/human CD45R/B220, Spark Blue™ 550 anti-mouse/human CD45R/B220, Spark NIR™ 685 anti-mouse/human CD45R/B220, TotalSeq™-B0103 anti-mouse/human CD45R/B220, Ultra-LEAF™ Purified anti-mouse/human CD45R/B220, TotalSeq™-C0103 anti-mouse/human CD45R/B220, PE/Fire™ 640 anti-mouse/human CD45R/B220, APC/Fire™ 810 anti-mouse/human CD45R/B220, PE/Fire™ 700 anti-mouse/human CD45R/B220, Spark Violet™ 538 anti-mouse/human CD45R/B220, Spark YG™ 581 anti-mouse/human CD45R/B220, Spark YG™ 570 anti-mouse/human CD45R/B220, PE/Fire™ 810 anti-mouse/human CD45R/B220, Spark Blue™ 574 anti-mouse/human CD45R/B220 Antibody, Spark Violet™ 423 anti-mouse/human CD45R/B220 Antibody, Spark Red™ 718 anti-mouse/human CD45R/B220


Product Data

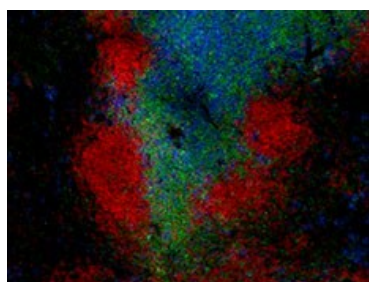


C57BL/6 mouse frozen lymph node section was fixed with 4% paraformaldehyde (PFA) for 10 minutes at room temperature and blocked with 5% FBS plus 5% rat serum for 1 hour at room temperature. Then the section was stained with 2.5 µg/mL of B220 (clone RA3-6B2) Alexa Fluor® 594 (red), 2.5 µg/mL of CD4 (clone RM4-5) Alexa Fluor® 488 (green), and 2.5 µg/mL of CD8 (clone 53-6.7) Alexa Fluor® 647 (blue) overnight at 4°C. The image was captured by 10X objective.

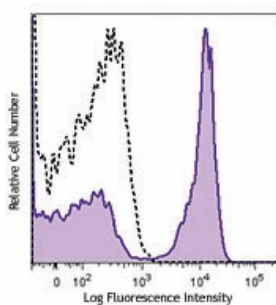

C57BL/6 mouse splenocytes were stained with CD45R/B220 (clone RA3-6B2) Alexa Fluor® 594 (filled histogram) or rat IgG2a, κ Alexa Fluor® 594 isotype control (open histogram). The data was acquired by BD LSRFortessa™ cell analyzer equipped with Yellow-Green Laser (561 nm).

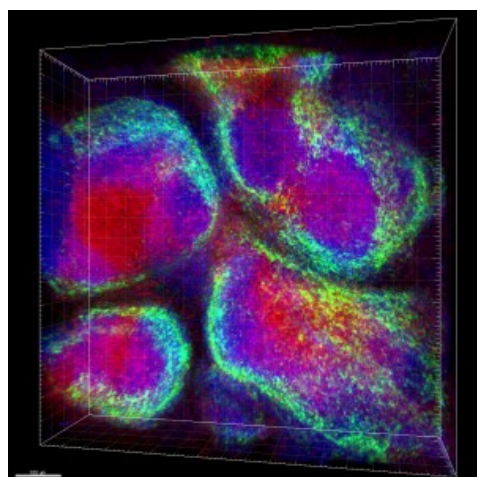

Formalin-fixed, 300 micron-thick mouse spleen section was blocked, permeabilized and stained overnight with CD4 (clone RM4-5) Alexa Fluor® 647 (red), CD169 (Siglec-1)(clone 3D6.112) Alexa Fluor® 488 (green), and CD45R/B220 (clone RA3-6B2) Alexa Fluor® 594 (blue) all at 5 µg/mL, optically cleared, then analyzed at 220 µm imaging depth on a confocal microscope. [Watch the video.](#)

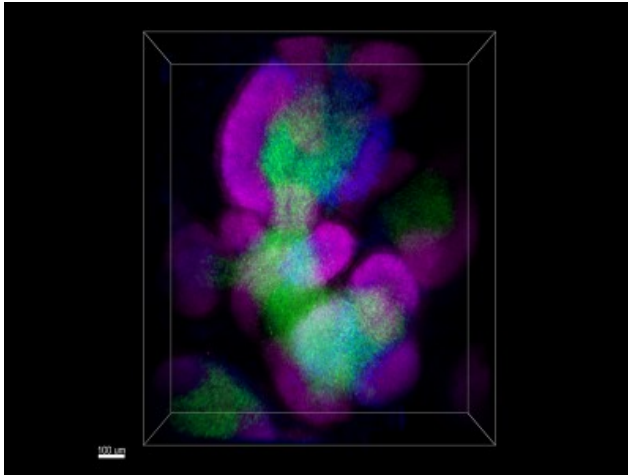

Paraformaldehyde-fixed (4%), 500 μm -thick mouse spleen section was processed according to the Ce3DTM Tissue Clearing Kit protocol (Cat. No. 427701). The section was costained with anti-mouse CD3 ϵ Antibody (clone 500A2) Alexa Fluor $^{\text{®}}$ 488 at 5 $\mu\text{g}/\text{mL}$ (green), and anti-mouse/human CD45R/B220 Antibody (clone RA3-6B2) Alexa Fluor $^{\text{®}}$ 594 at 5 $\mu\text{g}/\text{mL}$ (magenta) and counterstained with DAPI (blue). The section was then optically cleared and mounted in a sample chamber. The image was captured with a 10X objective using Zeiss 780 confocal microscope and processed by Imaris image analysis software.

[Watch the video.](#)

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587