

Brilliant Violet 421™ anti-mouse/human CD45R/B220 Antibody

Catalog# / Size	103239 / 125 µL 103251 / 50 µg 103240 / 500 µL
Clone	RA3-6B2
Regulatory Status	RUO
Other Names	B220
Isotype	Rat IgG2a, κ
Description	CD45R, also known as B220, is an isoform of CD45. It is a member of the protein tyrosine phosphatase (PTP) family with a molecular weight of approximately 180-240 kD. CD45R is expressed on B cells (at all developmental stages from pro-B cells through mature B cells), activated B cells, and subsets of T and NK cells. CD45R (B220) is also expressed on a subset of abnormal T cells involved in the pathogenesis of systemic autoimmunity in MRL- <i>Fas^{lpr}</i> and MRL- <i>Fas^{gld}</i> mice. It plays a critical role in TCR and BCR signaling. The primary ligands for CD45 are galectin-1, CD2, CD3, and CD4. CD45R is commonly used as a pan-B cell marker; however, CD19 may be more appropriate for B cell specificity.

Product Details

Verified Reactivity	Mouse, Human
Reported Reactivity	Cat
Antibody Type	Monoclonal
Host Species	Rat
Immunogen	Abelson murine leukemia virus-induced pre-B tumor cells
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA).
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 421™ under optimal conditions.
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested IHC-F - Verified
Recommended Usage	<p>Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 µl per million cells in 100 µl staining volume or 5 µl per 100 µl of whole blood. For immunohistochemistry on frozen tissue sections, a concentration range of 2.5 - 10 µg per ml. It is recommended that the reagent be titrated for optimal performance for each application.</p> <p>Brilliant Violet 421™ excites at 405 nm and emits at 421 nm. The standard bandpass filter 450/50 nm is recommended for detection. Brilliant Violet 421™ is a trademark of Sirigen Group Ltd.</p> <p>Learn more about Brilliant Violet™.</p> <p>This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.</p>
Excitation Laser	Violet Laser (405 nm)
Application Notes	Clone RA3-6B2 has been described to react with an epitope on the extracellular domain of the transmembrane CD45 glycoprotein which is dependent upon the expression of exon A and specific

carbohydrate residues. Additional reported applications (for the relevant formats) include: immunoprecipitation¹, *in vitro* and *in vivo* modulation of B cell responses^{2,4}, immunohistochemistry of acetone-fixed frozen sections and formalin-fixed paraffin-embedded sections^{5,6}, and spatial biology (IBEX)^{14,15}.

Application References

(PubMed link indicates BioLegend citation)

1. Coffman RL. 1982. *Immunol. Rev.* 69:5. (IP)
2. George A, *et al.* 1994. *J. Immunol.* 152:1014. (Activ)
3. Asensi V, *et al.* 1989. *Immunology* 68:204. (Activ)
4. Domiati-Saad R, *et al.* 1993. *J. Immunol.* 151:5936. (Activ)
5. Hata H, *et al.* 2004. *J. Clin. Invest.* 114:582. (IHC)
6. Monteith CE, *et al.* 1996. *Can. J. Vet. Res.* 60:193. (IHC)
7. Shih FF, *et al.* 2006. *J. Immunol.* 176:3438. (FC)
8. Chang C L-T, *et al.* 2007. *J. Immunol.* 178:6984.
9. Fazilleau N, *et al.* 2007. *Nature Immunol.* 8:753.
10. Lang GL, *et al.* 2008. *Blood* 111:2158. [PubMed](#)
11. Charles N, *et al.* 2010. *Nat. Med.* 16:701. (FC) [PubMed](#)
12. del Rio ML, *et al.* 2011. *Transpl. Int.* 24:501. (FC) [PubMed](#)
13. Murakami R, *et al.* 2013. *PLoS One.* 8:73270. [PubMed](#)
14. Radtke AJ, *et al.* 2020. *Proc Natl Acad Sci U S A.* 117:33455-65. (SB) [PubMed](#)
15. Radtke AJ, *et al.* 2022. *Nat Protoc.* 17:378-401. (SB) [PubMed](#)

Product Citations

1. Argüello RJ, *et al.* 2020. *Cell Metab.* 32:1063. [PubMed](#)
2. Wang X, *et al.* 2019. *Cell Res.* 29:787. [PubMed](#)
3. Baptista AP *et al.* 2019. *Immunity.* 50(5):1188-1201. [PubMed](#)
4. Mogilenko DA, *et al.* 2020. *Immunity.* 54(1):99-115.e12. [PubMed](#)
5. Sano T, *et al.* 2021. *Cell Reports.* 36(8):109608. [PubMed](#)
6. Chen J, *et al.* 2022. *J Nanobiotechnology.* 20:283. [PubMed](#)
7. Sugimoto C, *et al.* 2022. *Elife.* 11:. [PubMed](#)
8. Powell WE, *et al.* 2018. *Diabetologia.* 61:1794. [PubMed](#)
9. Hiramoto T, *et al.* 2018. *Mol Ther.* 26:1255. [PubMed](#)
10. Utley A, *et al.* 2020. *Cell Rep.* 31:107815. [PubMed](#)
11. Huang YL, *et al.* 2021. *Commun Biol.* 4:229. [PubMed](#)
12. Liu X, *et al.* 2022. *Nat Cancer.* .: [PubMed](#)
13. Liao T, *et al.* 2017. *Front Immunol.* 8:1334. [PubMed](#)
14. Marchalot A, *et al.* 2020. *Front Immunol.* 11:780. [PubMed](#)
15. Li J, *et al.* 2020. *Elife.* 9:00. [PubMed](#)
16. Venturutti L, *et al.* 2020. *Cell.* 182(2):297-316.e27. [PubMed](#)
17. Wittner J, *et al.* 2022. *Mucosal Immunol.* 15:668. [PubMed](#)
18. Bates PD, *et al.* 2021. *Front Immunol.* 12:668307. [PubMed](#)
19. Lindner J, *et al.* 2014. *J Immunol.* 192:4483. [PubMed](#)
20. Argüello RJ, *et al.* 2018. *J Cell Sci.* 131:jcs214346. [PubMed](#)
21. Siegemund S, *et al.* 2015. *PLoS One.* 10:124661. [PubMed](#)
22. Robertson TF, *et al.* 2021. *J Cell Biol.* 220:. [PubMed](#)
23. Bi E, *et al.* 2018. *J Clin Invest.* 128:4821. [PubMed](#)
24. Zwick M, *et al.* 2019. *Front Immunol.* 10:222. [PubMed](#)
25. Baomei Wang *et al.* 2019. *Cell reports.* 26(6):1614-1626. [PubMed](#)
26. Verheijen M, *et al.* 2020. *Cell Rep.* 108376:33. [PubMed](#)
27. Xu G, *et al.* 2020. *Cell Rep.* 32:108170. [PubMed](#)
28. Alexander Mildner *et al.* 2017. *Immunity.* 46(5):849-862. [PubMed](#)
29. Lazarian G, *et al.* 2021. *Cancer Cell.* 39(3):380-393.e8. [PubMed](#)
30. Girón-Pérez DA, *et al.* 2020. *J Cell Sci.* 133:00:00. [PubMed](#)
31. Heil J, *et al.* 2021. *Nat Commun.* 12:6963. [PubMed](#)
32. Koohy H, *et al.* 2018. *Genome Biol.* 0.879166667. [PubMed](#)
33. Yoshida H, *et al.* 2019. *Cell.* 176:897. [PubMed](#)
34. Han P, *et al.* 2020. *Sci Adv.* 6:eaaz1580. [PubMed](#)
35. Nagatake T, *et al.* 2019. *Mucosal Immunol.* 12:p1082. [PubMed](#)
36. Lightman SM, *et al.* 2020. *Cell Reports.* 1.086805556. [PubMed](#)
37. Tan DQ, *et al.* 2019. *Cell Rep.* 26:2316. [PubMed](#)
38. Yáñez A *et al.* 2017. *Immunity.* 47(5):890-902. [PubMed](#)

RRID

AB_10933424 (BioLegend Cat. No. 103239)
AB_2562905 (BioLegend Cat. No. 103251)
AB_11203896 (BioLegend Cat. No. 103240)

Antigen Details

Structure	Protein tyrosine phosphatase (PTP) family, 180-240 kD
Distribution	B cells, T cell subset, NK cell subset
Function	Phosphatase, T and B cell activation
Ligand/Receptor	Galectin-1, CD2, CD3, CD4

Cell Type	B cells, NK cells, T cells
Biology Area	Cell Biology, Immunology, Inhibitory Molecules, Neuroscience, Neuroscience Cell Markers
Molecular Family	CD Molecules
Antigen References	<ol style="list-style-type: none"> 1. Barclay A, <i>et al.</i> 1997. The Leukocyte Antigen FactsBook Academic Press. 2. Trowbridge IS, <i>et al.</i> 1993. <i>Annu. Rev. Immunol.</i> 12:85. 3. Kishihara K, <i>et al.</i> 1993. <i>Cell</i> 74:143. 4. Pulido R, <i>et al.</i> 1988. <i>J. Immunol.</i> 140:3851.
Gene ID	19264 5788

Related Protocols


[Immunohistochemistry Protocol for Frozen Sections](#)

[Cell Surface Flow Cytometry Staining Protocol](#)


Other Formats

Alexa Fluor® 594 anti-mouse/human CD45R/B220, APC anti-mouse/human CD45R/B220, Biotin anti-mouse/human CD45R/B220, FITC anti-mouse/human CD45R/B220, PE anti-mouse/human CD45R/B220, PE/Cyanine5 anti-mouse/human CD45R/B220, Purified anti-mouse/human CD45R/B220, PE/Cyanine7 anti-mouse/human CD45R/B220, APC/Cyanine7 anti-mouse/human CD45R/B220, Alexa Fluor® 488 anti-mouse/human CD45R/B220, Alexa Fluor® 647 anti-mouse/human CD45R/B220, Pacific Blue™ anti-mouse/human CD45R/B220, Alexa Fluor® 700 anti-mouse/human CD45R/B220, PerCP anti-mouse/human CD45R/B220, PerCP/Cyanine5.5 anti-mouse/human CD45R/B220, Brilliant Violet 421™ anti-mouse/human CD45R/B220, Brilliant Violet 570™ anti-mouse/human CD45R/B220, Brilliant Violet 650™ anti-mouse/human CD45R/B220, Brilliant Violet 605™ anti-mouse/human CD45R/B220, Brilliant Violet 785™ anti-mouse/human CD45R/B220, Brilliant Violet 510™ anti-mouse/human CD45R/B220, Purified anti-mouse/human CD45R/B220 (Maxpar® Ready), Brilliant Violet 711™ anti-mouse/human CD45R/B220, PE/Dazzle™ 594 anti-mouse/human CD45R/B220, APC/Fire™ 750 anti-mouse/human CD45R/B220, Brilliant Violet 750™ anti-mouse/human CD45R/B220, TotalSeq™-A0103 anti-mouse/human CD45R/B220, Spark Blue™ 550 anti-mouse/human CD45R/B220, Spark NIR™ 685 anti-mouse/human CD45R/B220, TotalSeq™-B0103 anti-mouse/human CD45R/B220, Ultra-LEAF™ Purified anti-mouse/human CD45R/B220, TotalSeq™-C0103 anti-mouse/human CD45R/B220, PE/Fire™ 640 anti-mouse/human CD45R/B220, APC/Fire™ 810 anti-mouse/human CD45R/B220, PE/Fire™ 700 anti-mouse/human CD45R/B220, Spark Violet™ 538 anti-mouse/human CD45R/B220, Spark YG™ 581 anti-mouse/human CD45R/B220, Spark YG™ 570 anti-mouse/human CD45R/B220, PE/Fire™ 810 anti-mouse/human CD45R/B220, Spark Blue™ 574 anti-mouse/human CD45R/B220 Antibody, Spark Violet™ 423 anti-mouse/human CD45R/B220 Antibody, Spark Red™ 718 anti-mouse/human CD45R/B220

Product Data

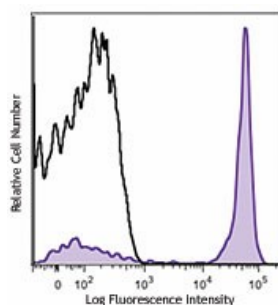

C57BL/6 mouse splenocytes were stained with CD45R/B220 (clone RA3-6B2) Brilliant Violet 421™ (filled histogram) or rat IgG2a, κ Brilliant Violet 421™ isotype control (open histogram).

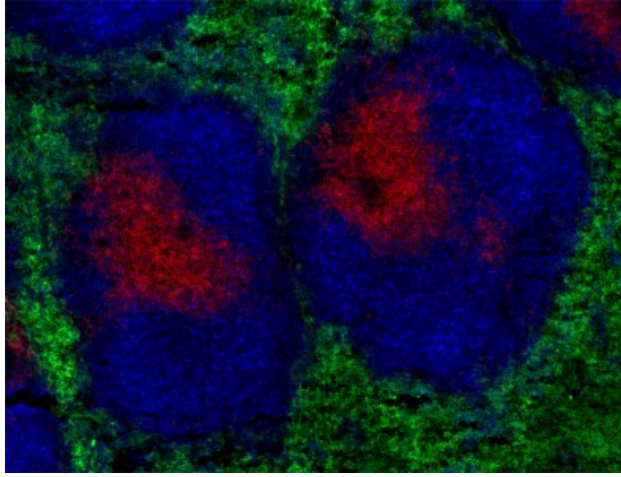

C57BL/6 mouse frozen spleen sections were fixed with 4% paraformaldehyde (PFA) for ten minutes at room temperature and blocked with 5% FBS plus 5% rat/mouse serum for 30 minutes at room temperature. Then the section was stained with 1.25 µg/ml of anti-mouse/human CD45R/B220 (clone RA3-6B2) Brilliant Violet 421™ (blue), anti-mouse CD5 (clone 53-7.3) Alexa Fluor® 594 (red), and anti-mouse F4/80 (clone BM8) Alexa Fluor® 647 (green) overnight at 4°C. The image was captured with a 10X objective.

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587