

Brilliant Violet 605™ anti-human CD95 (Fas) Antibody

Catalog# / Size	305627 / 25 tests 305628 / 100 tests
Clone	DX2
Regulatory Status	RUO
Workshop	VI C-64
Other Names	Fas, APO-1, TNFRSF6
Isotype	Mouse IgG1, κ
Description	CD95 is a 45 kD single chain type I glycoprotein also known as Fas, APO-1, and TNFRSF6. It is a member of the TNF receptor superfamily. CD95 is expressed on T and B lymphocytes, monocytes, neutrophils, and fibroblasts. CD95 expression is upregulated by activation. The extracellular region of CD95 binds to CD178 (Fas ligand). CD178 binding to CD95 induces apoptosis and has been shown to play a role in the maintenance of peripheral tolerance.

Product Details

Verified Reactivity	Human, Cynomolgus, Rhesus
Reported Reactivity	African Green, Baboon, Capuchin Monkey, Chimpanzee, Common Marmoset, Cotton-topped Tamarin, Pigtailed Macaque, Sooty Mangabey
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	CD95 transfected L cells
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA).
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 605™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration, please enter the lot number in our Concentration and Expiration Lookup or Certificate of Analysis online tools.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 µl per million cells in 100 µl staining volume or 5 µl per 100 µl of whole blood. Brilliant Violet 605™ excites at 405 nm and emits at 603 nm. The bandpass filter 610/20 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 605™ is a trademark of Sirigen Group Ltd. Learn more about Brilliant Violet™. This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.
Excitation Laser	Violet Laser (405 nm)

Application Notes

The DX2 antibody is useful for inducing apoptosis of Fas-positive cells. Additional reported applications (for the relevant formats) include: *in vitro* induction of apoptosis³ (DX2 antibody is required to be cross-linked for effective induction of apoptosis) and immunohistochemical staining^{4,5} of acetone-fixed frozen tissue sections and formalin-fixed paraffin-embedded tissue sections. The Ultra-LEAF™ purified antibody (Endotoxin < 0.01 EU/μg, Azide-Free, 0.2 μm filtered) is recommended for functional assays (Cat. No. 305655 and 305656).

Note: EOS9.1 antibody can induce apoptosis without cross-linking.

Application References

(PubMed link indicates BioLegend citation)

1. Schlossman S, *et al.* Eds. 1995. Leucocyte Typing V. Oxford University Press. New York.
2. Kishimoto T, *et al.* Eds. 1997. Leucocyte Typing VI. Garland Publishing Inc. New York.
3. Cifone M, *et al.* 1994. *J. Exp. Med.* 180:1547. (Apop)
4. Zietz C, *et al.* 2001. *Am. J. Pathol.* 159:963. (IHC)
5. Sergi C, *et al.* 2000. *Am. J. Pathol.* 156:1589. (IHC)
6. Xie S, *et al.* 2010. *J. Immunol.* 184:2289. (FC) [PubMed](#)
7. Yoshino N, *et al.* 2000. *Exp. Anim. (Tokyo)* 49:97. (FC)
8. Sestak K, *et al.* 2007. *Vet. Immunol. Immunopathol.* 119:21.
9. Rout N, *et al.* 2010. *PLoS One* 5:e9787. (FC)
10. Dixit N, *et al.* 2012. *J. Immunol.* 189:5954. [PubMed](#)

Product Citations

1. Hoang TN, *et al.* 2021. *Cell.* 184:460. [PubMed](#)
2. Meehan S, *et al.* 2019. *Commun Biol.* 2:229. [PubMed](#)
3. Colleen S McGary *et al.* 2017. *Immunity.* 47(4):776-788 . [PubMed](#)
4. Hoang TN, *et al.* 2020. *bioRxiv.* . [PubMed](#)
5. Barcelo H, *et al.* 2018. *Curr Protoc Cytom.* 84:e35. [PubMed](#)
6. Harper JL, *et al.* 2020. *Nat Med.* 519:26. [PubMed](#)
7. Violari A, *et al.* 2019. *Nat Commun.* 10:412. [PubMed](#)
8. Zhang P, *et al.* 2021. *Nat Med.* 22:. [PubMed](#)
9. Thyagarajan B, *et al.* 2018. *J Immunol Methods.* 463:61. [PubMed](#)
10. Vierboom MPM, *et al.* 2021. *Cell Reports Medicine.* 2(1):100185. [PubMed](#)
11. Dijkman K, *et al.* 2021. *Cell Reports Medicine.* 2(1):100187. [PubMed](#)

RRID

AB_2562444 (BioLegend Cat. No. 305627)
AB_2563825 (BioLegend Cat. No. 305628)

Antigen Details

Structure	TNFR superfamily, type I transmembrane protein, 45 kD
Distribution	T cells and B cells, monocytes, neutrophils, fibroblasts, expression level upregulated on activated lymphocytes
Function	Mediates apoptosis
Ligand/Receptor	Fas ligand (CD178)
Cell Type	B cells, Fibroblasts, Lymphocytes, Monocytes, Neutrophils, T cells
Biology Area	Apoptosis/Tumor Suppressors/Cell Death, Cell Biology, Immunology, Neuroscience
Molecular Family	CD Molecules
Antigen References	1. Krammer P, <i>et al.</i> 1994. <i>Immunol. Rev.</i> 142:175. 2. Nagata S, <i>et al.</i> 1995. <i>Science</i> 267:1449.
Gene ID	355

Related Protocols


[Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

APC anti-human CD95 (Fas), Biotin anti-human CD95 (Fas), FITC anti-human CD95 (Fas), PE anti-human CD95 (Fas), PE/Cyanine5 anti-human CD95 (Fas), Purified anti-human CD95 (Fas), Alexa Fluor® 488 anti-human CD95 (Fas), Alexa Fluor® 647 anti-human CD95 (Fas), Brilliant Violet 421™ anti-human CD95 (Fas), Pacific Blue™ anti-human CD95 (Fas), PE/Cyanine7 anti-

human CD95 (Fas), Brilliant Violet 605™ anti-human CD95 (Fas), PerCP/Cyanine5.5 anti-human CD95 (Fas), Purified anti-human CD95 (Fas) (Maxpar® Ready), PE/Dazzle™ 594 anti-human CD95 (Fas), APC/Fire™ 750 anti-human CD95 (Fas), APC/Cyanine7 anti-human CD95 (Fas), Brilliant Violet 510™ anti-human CD95 (Fas), Brilliant Violet 711™ anti-human CD95 (Fas), Brilliant Violet 785™ anti-human CD95 (Fas), Brilliant Violet 650™ anti-human CD95 (Fas), Alexa Fluor® 700 anti-human CD95 (Fas), TotalSeq™-A0156 anti-human CD95 (Fas), TotalSeq™-C0156 anti-human CD95 (Fas), TotalSeq™-B0156 anti-human CD95 (Fas), Ultra-LEAF™ Purified anti-human CD95 (Fas), TotalSeq™-D0156 anti-human CD95 (Fas), PE/Fire™ 640 anti-human CD95 (Fas), KIRAVIA Blue 520™ anti-human CD95 (Fas), APC/Fire™ 810 anti-human CD95 (Fas) Antibody

Product Data

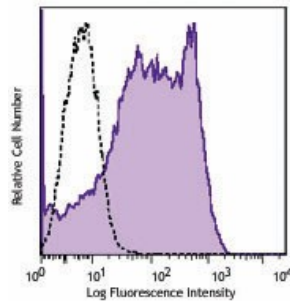

Human peripheral blood lymphocytes were stained with CD95 (clone DX2) Brilliant Violet 605™ (filled histogram) or mouse IgG1, κ Brilliant Violet 605™ (open histogram).

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587