

Brilliant Violet 785™ anti-mouse CD11c Antibody

Catalog# / Size	117335 / 125 μ L 117336 / 50 μ g
Clone	N418
Regulatory Status	RUO
Other Names	α X integrin, integrin α X chain, CR4, p150, ITGAX
Isotype	Armenian Hamster IgG
Description	CD11c is a 150 kD glycoprotein also known as α X integrin, CR4, and p150. CD11c forms a α X β 2 heterodimer with β 2 integrin (CD18). It is primarily expressed on dendritic cells, NK cells, a subset of intestinal intraepithelial lymphocytes (IEL), and some activated T cells. The α X β 2 integrin plays an important role in cell-cell contact by binding its ligands: iC3b, fibrinogen, and CD54.

Product Details

Verified Reactivity	Mouse
Antibody Type	Monoclonal
Host Species	Armenian Hamster
Immunogen	Mouse spleen dendritic cells
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA).
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 785™ under optimal conditions.
Concentration	μ g sizes: 0.2 mg/mL μ L sizes: lot-specific (to obtain lot-specific concentration, please enter the lot number in our Concentration and Expiration Lookup or Certificate of Analysis online tools.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	<p>Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining using the μg size, the suggested use of this reagent is \leq0.25 μg per million cells in 100 μl volume. For flow cytometric staining using the μl size, the suggested use of this reagent is 5 μl per million cells in 100 μl staining volume or 5 μl per 100 μl of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.</p> <p>Brilliant Violet 785™ excites at 405 nm and emits at 785 nm. The bandpass filter 780/60 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 785™ is a trademark of Sirigen Group Ltd.</p> <p>Learn more about Brilliant Violet™.</p> <p>This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.</p>
Excitation Laser	Violet Laser (405 nm)
Application Notes	Additional reported applications (for the relevant formats) include: immunoprecipitation ³ , immunohistochemical staining of acetone-fixed frozen sections ³ , immunofluorescence microscopy ⁵ .

⁹ (Alexa Fluor® 488 conjugated N418 was used for IHC in frozen sections¹⁰), and spatial biology (IBEX)^{22,23}.

Application References

(PubMed link indicates BioLegend citation)

1. Granucci F, et al. 1997. *J. Immunol.* 159:1794.
2. Stokes RW, et al. 1998. *J. Immunol.* 160:5514.
3. Metlay JP, et al. 1990. *J. Exp. Med.* 171:1753. (IHC, IP)
4. Ma XT, et al. 2006. *Cancer Research* 66:1169.
5. Chin RK, et al. 2006. *J. Immunol.* 177:290. (IF)
6. Cervantes-Barragan L, et al. 2007. *Blood* 109:1131. (FC) [PubMed](#)
7. Turnquist HR, et al. 2007. *J. Immunol.* 178:7018. (FC) [PubMed](#)
8. Benson MJ, et al. 2007. *J. Exp. Med.* doi:10.1084/jem.20070719. (FC) [PubMed](#)
9. You Y, et al. 2009. *J. Immunol.* 182:7343. (IF) [PubMed](#)
10. Roland CL, et al. 2009. *Mol. Cancer Res.* 8:1761. (IHC, FC) [PubMed](#)
11. Wikstrom M, et al. 2006. *J. Immunol.* 177:913. [PubMed](#)
12. Pericolini E, et al. 2008. *J. Leukocyte Biol.* 83:1286. [PubMed](#)
13. Randall LM, et al. 2008. *Infect. Immun.* 76:3312. [PubMed](#)
14. Fahlen-Yrild L, et al. 2009. *J. Immunol.* 183:5032. [PubMed](#)
15. Osterholzer JJ, et al. 2009. *J. Immunol.* 183:8044. [PubMed](#)
16. Bankoti J, et al. 2010. *Toxicol. Sci.* 115:422. (FC) [PubMed](#)
17. Eisenach PA, et al. 2010. *J Cell Sci.* 123:4182. [PubMed](#)
18. Leppin K, et al. 2014. *Invest. Ophthalmol. Vis. Sci.* 55:3603. [PubMed](#)
19. Sakai F, et al. 2014. *PLoS One.* 9:105370. [PubMed](#)
20. Gibbins JD, et al. 2014. *Blood.* 124:2953. [PubMed](#)
21. White CE, et al. 2015. *J Immunol.* 194:697. [PubMed](#)
22. Lu X, et al. 2015. *J Immunol.* 194:2011. [PubMed](#)
23. Radtke AJ, et al. 2020. *Proc Natl Acad Sci U S A.* 117:33455-65. (SB) [PubMed](#)
24. Radtke AJ, et al. 2022. *Nat Protoc.* 17:378-401. (SB) [PubMed](#)

Product Citations

1. Oliver AJ, et al. 2020. *Oncoimmunology.* 9:1802979. [PubMed](#)
2. Apte SH, et al. 2020. *Clin Transl Immunology.* 9:e1209. [PubMed](#)
3. Hu K, et al. 2022. *J Extracell Vesicles.* 11:e12199. [PubMed](#)
4. Han SJ et al. 2017. *Immunity.* 47(6):1154-1168. [PubMed](#)
5. Arce Vargas F et al. 2018. *Cancer cell.* 33(4):649-663. [PubMed](#)
6. Gruber T, et al. 2020. *JCI Insight.* 5:00. [PubMed](#)
7. Hays C, et al. 2020. *eLife.* 8:e48772. [PubMed](#)
8. Lechner AJ et al. 2017. *Cell stem cell.* 21(1):120-134. [PubMed](#)
9. Chow MT et al. 2019. *Immunity.* 50(6):1498-1512. [PubMed](#)
10. Tordesillas L, et al. 2018. *Nat Commun.* 9:5238. [PubMed](#)
11. Kremenovic M, et al. 2022. *J Immunother Cancer.* 10:. [PubMed](#)
12. Gómez-Díaz C, et al. 2021. *iScience.* 24:103241. [PubMed](#)
13. Zhang H, et al. 2020. *Cancer Cell.* 37(1):37-54.e9. [PubMed](#)
14. Devos M, et al. 2020. *J Exp Med.* 217:00:00. [PubMed](#)
15. Casadesús AV, et al. 2020. *Oncoimmunology.* 9:1770565. [PubMed](#)
16. Wu X, et al. 2021. *Elife.* 10:. [PubMed](#)
17. Park HB, et al. 2020. *Oncoimmunology.* 9:1772663. [PubMed](#)
18. Funk KE, et al. 2019. *J Neuroinflammation.* 16:22. [PubMed](#)
19. Watanabe K, et al. 2018. *JCI Insight.* 3. [PubMed](#)
20. Linehan JL et al. 2018. *Cell.* 172(4):784-796. [PubMed](#)
21. Fennell LM, et al. 2020. *EMBO J.* 39:e103303. [PubMed](#)
22. Schepanski S, et al. 2022. *Nat Commun.* 13:4571. [PubMed](#)
23. AM C, et al. 2015. *Proc Natl Acad Sci U S A.* 112 7557. [PubMed](#)
24. Ji DX, et al. 2019. *Nat Microbiol.* 1.644444444. [PubMed](#)
25. Tan X, et al. 2021. *Sci Adv.* 7: . [PubMed](#)
26. Bourque J, et al. 2021. *Heliyon.* 7:e08311. [PubMed](#)
27. Gilbert NM, et al. 2021. *Biol Open.* 10:. [PubMed](#)
28. Liu X, et al. 2020. *Cell Host Microbe.* 28(5):683-698.e6. [PubMed](#)
29. Oguri Y, et al. 2020. *Cell.* 182(3):563-577.e20. [PubMed](#)
30. Brown SM et al. 2017. *Endocrinology.* 158(10):3592-3604. [PubMed](#)
31. Mohamed E, et al. 2020. *Immunity.* 52(4):668-682.e7. [PubMed](#)
32. Andriessen EMMA, et al. 2021. *EMBO Mol Med.* 13:e11754. [PubMed](#)
33. Petkevicius K, et al. 2019. *Elife.* 8:e47990. [PubMed](#)
34. Baranwal G, et al. 2021. *Physiol Rep.* 9:e15094. [PubMed](#)
35. Hayashida E, et al. 2019. *J Neuroinflammation.* 0.789583333. [PubMed](#)
36. Fedele C, et al. 2021. *J Exp Med.* 218: . [PubMed](#)
37. E O'Koren, R Mathew, D Saban 2016. *Sci Rep.* 6:20636. [PubMed](#)
38. Jaeger N, et al. 2020. *Cell Rep.* 33:108331. [PubMed](#)
39. Gangoso E, et al. 2021. *Cell.* 184:2454. [PubMed](#)
40. Nascimento DC, et al. 2021. *Immunity.* . [PubMed](#)
41. McGinley AM, et al. 2020. *Immunity.* 52(2):342-356. [PubMed](#)

RRID

AB_11219204 (BioLegend Cat. No. 117335)
AB_2565268 (BioLegend Cat. No. 117336)

Antigen Details

Structure	Integrin α -chain, associates with integrin β_2 (CD18), 150 kD
Distribution	Dendritic cells, NK cells, intestinal intraepithelial lymphocytes (IEL), some activated T cells
Function	Cellular adhesion
Ligand/Receptor	iC3b, fibrinogen
Cell Type	Dendritic cells, Epithelial cells, NK cells, T cells, Tregs
Biology Area	Cell Adhesion, Cell Biology, Costimulatory Molecules, Immunology, Innate Immunity, Neuroscience, Neuroscience Cell Markers
Molecular Family	Adhesion Molecules, CD Molecules
Antigen References	1. Barclay A, <i>et al.</i> 1997. The Leukocyte Antigen Facts Book Academic Press. 2. Springer TA. 1994. <i>Cell</i> 76:301. 3. Lopez-Rodriguez C, <i>et al.</i> 1996. <i>J. Immunol.</i> 156:3780.
Gene ID	16411


Related Protocols

[Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

APC anti-mouse CD11c, Biotin anti-mouse CD11c, FITC anti-mouse CD11c, PE anti-mouse CD11c, Purified anti-mouse CD11c, Alexa Fluor® 488 anti-mouse CD11c, Alexa Fluor® 647 anti-mouse CD11c, PE/Cyanine5 anti-mouse CD11c, PE/Cyanine7 anti-mouse CD11c, Brilliant Violet 605™ anti-mouse CD11c, Alexa Fluor® 700 anti-mouse CD11c, Pacific Blue™ anti-mouse CD11c, APC/Cyanine7 anti-mouse CD11c, PerCP/Cyanine5.5 anti-mouse CD11c, PerCP anti-mouse CD11c, Brilliant Violet 421™ anti-mouse CD11c, Brilliant Violet 570™ anti-mouse CD11c, Brilliant Violet 785™ anti-mouse CD11c, Brilliant Violet 510™ anti-mouse CD11c, Brilliant Violet 650™ anti-mouse CD11c, Purified anti-mouse CD11c (Maxpar® Ready), Alexa Fluor® 594 anti-mouse CD11c, PE/Dazzle™ 594 anti-mouse CD11c, Brilliant Violet 711™ anti-mouse CD11c, APC/Fire™ 750 anti-mouse CD11c, TotalSeq™-A0106 anti-mouse CD11c, Brilliant Violet 750™ anti-mouse CD11c, TotalSeq™-B0106 anti-mouse CD11c, TotalSeq™-C0106 anti-mouse CD11c, KIRAVIA Blue 520™ anti-mouse CD11c, Spark Blue™ 550 anti-mouse CD11c, Spark NIR™ 685 anti-mouse CD11c, Spark UV™ 387 anti-mouse CD11c, Spark Red™ 718 anti-mouse CD11c

Product Data

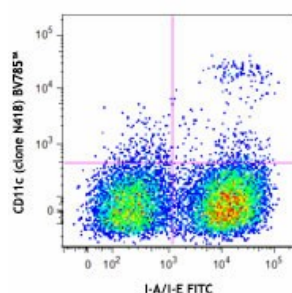

C57BL/6 mouse splenocytes were stained with I-A/I-E FITC and CD11c (clone N418) Brilliant Violet 785™.

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587

