

## Alexa Fluor® 488 anti-human CD56 (NCAM) Antibody

<b>Catalog# / Size</b>	318311 / 25 tests 318312 / 100 tests
<b>Clone</b>	HCD56
<b>Regulatory Status</b>	RUO
<b>Other Names</b>	Leu-19, NKH1
<b>Isotype</b>	Mouse IgG1, κ
<b>Description</b>	CD56 is a single transmembrane glycoprotein also known as NCAM (Neural Cell Adhesion Molecule), Leu-19, or NKH1. It is a member of the Ig superfamily. The 140 kD isoform is expressed on NK cells and NK-T cells. CD56 is also expressed in the brain (cerebellum and cortex) and at neuromuscular junctions. Certain large granular lymphocyte (LGL) leukemias, small-cell lung carcinomas, neuronal derived tumors, myelomas, and myeloid leukemias also express CD56. CD56 plays a role in homophilic and heterophilic adhesion via binding to itself or heparin sulfate.

### Product Details

<b>Verified Reactivity</b>	Human
<b>Reported Reactivity</b>	African Green, Baboon, Cynomolgus, Rhesus
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Mouse
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
<b>Preparation</b>	The antibody was purified by affinity chromatography and conjugated with Alexa Fluor® 488 under optimal conditions.
<b>Concentration</b>	Lot-specific (to obtain lot-specific concentration, please enter the lot number in our <a href="#">Concentration and Expiration Lookup</a> or <a href="#">Certificate of Analysis</a> online tools.)
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a> . For flow cytometric staining, the suggested use of this reagent is 5 µl per million cells in 100 µl staining volume or 5 µl per 100 µl of whole blood.  * Alexa Fluor® 488 has a maximum emission of 519 nm when it is excited at 488 nm.  Alexa Fluor® and Pacific Blue™ are trademarks of Life Technologies Corporation.  <a href="#">View full statement regarding label licenses</a>
<b>Excitation Laser</b>	Blue Laser (488 nm)
<b>Application Notes</b>	Clone HCD56 is not recommended for immunohistochemistry formalin-fixed paraffin-embedded tissue.
<b>Application References</b> (PubMed link indicates BioLegend citation)	<ol style="list-style-type: none"> <li>1. Kishimoto T, <i>et al.</i> Eds. 1997. Leucocyte Typing VI. Garland Publishing Inc. London.</li> <li>2. Correia DV, <i>et al.</i> 2011. <i>Blood</i> 118:992. (FC) <a href="#">PubMed</a></li> </ol>
<b>Product Citations</b>	<ol style="list-style-type: none"> <li>1. Karampatzakis A, <i>et al.</i> 2021. <i>Front Immunol.</i> 12:641521. <a href="#">PubMed</a></li> <li>2. Ângela Crespo, Jack Strominger, Tamara Tilburgs 2016. <i>Proc Natl Acad Sci U S A.</i> 113(52):15072-15077. <a href="#">PubMed</a></li> </ol>

- Ferry GM, *et al.* 2022. *Front Immunol.* 13:863155. [PubMed](#)
- Cohen Z, *et al.* 2020. *Biomed Rep.* 13:17. [PubMed](#)

**RRID** AB\_604094 (BioLegend Cat. No. 318311)  
AB\_604102 (BioLegend Cat. No. 318312)

## Antigen Details

<b>Structure</b>	Ig superfamily, single transmembrane or GPI-anchored glycoprotein
<b>Distribution</b>	NK cells, T subset, neural tissue, some LGL and myeloid leukemias
<b>Function</b>	Adhesion
<b>Ligand/Receptor</b>	Heparin sulfate
<b>Cell Type</b>	B cells, Leukemia, Mesenchymal Stem Cells, Neurons, NK cells, T cells
<b>Biology Area</b>	Cell Adhesion, Cell Biology, Costimulatory Molecules, Immunology, Innate Immunity, Neuroscience, Stem Cells, Synaptic Biology
<b>Molecular Family</b>	Adhesion Molecules, CD Molecules
<b>Antigen References</b>	1. Lanier L, <i>et al.</i> 1991. <i>J. Immunol.</i> 146:4421. 2. Hemperly J, <i>et al.</i> 1990. <i>J. Mol. Neurosci.</i> 2:71. 3. Cremer H, <i>et al.</i> 1994. <i>Nature</i> 367:455.
<b>Gene ID</b>	<a href="#">4684</a>


## Related Protocols

[Cell Surface Flow Cytometry Staining Protocol](#)

## Other Formats

Purified anti-human CD56 (NCAM), FITC anti-human CD56 (NCAM), PE anti-human CD56 (NCAM), PE/Cyanine5 anti-human CD56 (NCAM), APC anti-human CD56 (NCAM), Alexa Fluor® 488 anti-human CD56 (NCAM), Alexa Fluor® 647 anti-human CD56 (NCAM), Alexa Fluor® 700 anti-human CD56 (NCAM), PE/Cyanine7 anti-human CD56 (NCAM), Biotin anti-human CD56 (NCAM), PerCP/Cyanine5.5 anti-human CD56 (NCAM), Pacific Blue™ anti-human CD56 (NCAM), APC/Cyanine7 anti-human CD56 (NCAM), Brilliant Violet 421™ anti-human CD56 (NCAM), Brilliant Violet 570™ anti-human CD56 (NCAM), Brilliant Violet 605™ anti-human CD56 (NCAM), Brilliant Violet 711™ anti-human CD56 (NCAM), Brilliant Violet 510™ anti-human CD56 (NCAM), PerCP anti-human CD56 (NCAM), Brilliant Violet 650™ anti-human CD56 (NCAM), Purified anti-human CD56 (NCAM) (Maxpar® Ready), PE/Dazzle™ 594 anti-human CD56 (NCAM)

## Product Data

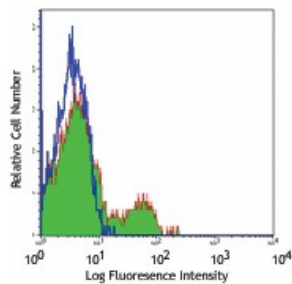

Human peripheral blood lymphocytes  
stained with HCD56 Alexa Fluor® 488

\*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, [www.biolegend.com/ordering#license](http://www.biolegend.com/ordering#license)). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587