

Brilliant Violet 650™ anti-human CD123 Antibody

Catalog# / Size	306019 / 25 tests 306020 / 100 tests
Clone	6H6
Regulatory Status	RUO
Other Names	IL-3R α , IL-3 Receptor alpha
Isotype	Mouse IgG1, κ
Description	CD123 is the 70 kD transmembrane α chain of the IL-3 receptor. Alone, CD123 binds IL-3 with low affinity; when CD123 associates with CDw131 (common β chain), it binds IL-3 with high affinity. CD123 does not transduce intracellular signals upon binding IL-3 and requires the β chain for this function. CD123 is expressed by myeloid precursors, macrophages, dendritic cells, mast cells, basophils, megakaryocytes, and some B cells.

Product Details

Verified Reactivity	Human
Reported Reactivity	Rhesus
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	Human IL-3R α transfected COS cells.
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA).
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 650™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration, please enter the lot number in our Concentration and Expiration Lookup or Certificate of Analysis online tools.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	<p>Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 μl per million cells in 100 μl staining volume or 5 μl per 100 μl of whole blood.</p> <p>Brilliant Violet 650™ excites at 405 nm and emits at 645 nm. The bandpass filter 660/20 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 650™ is a trademark of Sirigen Group Ltd.</p> <p>Learn more about Brilliant Violet™.</p> <p>This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.</p>
Excitation Laser	Violet Laser (405 nm)
Application Notes	Clone 6H6 does not inhibit IL-3 binding to low- or high-affinity IL-3Rs. Additional reported applications (for the relevant formats) include: Western blotting ¹ , immunoprecipitation ¹ , and immunohistochemical staining of acetone-fixed frozen sections ² and also paraformaldehyde fixed

paraffin embedded tissue⁷.

Application References

(PubMed link indicates BioLegend citation)

1. Sun Q, *et al.* 1996. *Blood* 87:83. (IP, WB)
2. Herling M, *et al.* 2003. *Blood* 101:5007. (IHC)
3. Charles N, *et al.* 2010. *Nat. Med.* 16:701. (FC) [PubMed](#)
4. Martin-Gayo E, *et al.* 2010. *Blood* 115:5366. [PubMed](#)
5. Chen SC, *et al.* 2010. *Arch Dermatol Res.* 302:113. [PubMed](#)
6. Liu Y, *et al.* 2012. *Food Chem Toxicol.* 50:1920. [PubMed](#)
7. Peduzzi E, *et al.* 2007. *J. Invest. Dermatol.* 127:638. (IHC)

Product Citations

1. de Mingo Pulido, *et al.* 2021. *Immunity.* 54(6):1154-1167.e7. [PubMed](#)
2. Herter JM, *et al.* 2022. *Strahlenther Onkol.* Online ahead of print. [PubMed](#)
3. Cillo AR, *et al.* 2020. *Immunity.* 52(1):183-199.e9. [PubMed](#)
4. Lu Q, *et al.* 2021. *Immunity.* 54(6):1304-1319.e9. [PubMed](#)
5. Heger L, *et al.* 2018. *Front Immunol.* 9:744. [PubMed](#)
6. Swanson E, *et al.* 2021. *eLife.* 10:00. [PubMed](#)
7. Mishra A, *et al.* 2018. *J Allergy Clin Immunol.* 142:1066. [PubMed](#)
8. Savage AK, *et al.* 2021. *iScience.* 24(5):102404. [PubMed](#)
9. Genge PC, *et al.* 2021. *STAR Protoc.* 2:100900. [PubMed](#)
10. Abd Hamid M *et al.* 2019. *Cancer Immunol Res.* 7(8):1293-1306. [PubMed](#)

RRID

AB_11218792 (BioLegend Cat. No. 306019)
AB_2563827 (BioLegend Cat. No. 306020)

Antigen Details

Structure	Ig superfamily, type I transmembrane glycoprotein, associates with CDw131, 70 kD
Distribution	Myeloid precursors, basophils, mast cells, macrophages, dendritic cells, megakaryocytes, subset of lymphocytes
Function	Hematopoietic cell proliferation, differentiation
Ligand/Receptor	IL-3
Cell Type	Basophils, Dendritic cells, Hematopoietic stem and progenitors, Lymphocytes, Macrophages, Mast cells, Megakaryocytes
Biology Area	Immunology
Molecular Family	CD Molecules, Cytokine/Chemokine Receptors
Antigen References	1. Miyajima A, <i>et al.</i> 1993. <i>Blood</i> 82:1960.
Gene ID	3563



Related Protocols

[Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

Biotin anti-human CD123, PE anti-human CD123, Purified anti-human CD123, PE/Cyanine5 anti-human CD123, PE/Cyanine7 anti-human CD123, APC anti-human CD123, FITC anti-human CD123, PerCP/Cyanine5.5 anti-human CD123, Brilliant Violet 421™ anti-human CD123, Brilliant Violet 650™ anti-human CD123, Brilliant Violet 510™ anti-human CD123, Alexa Fluor® 647 anti-human CD123, Brilliant Violet 605™ anti-human CD123, Purified anti-human CD123 (Maxpar® Ready), Brilliant Violet 711™ anti-human CD123, Brilliant Violet 785™ anti-human CD123, PE/Dazzle™ 594 anti-human CD123, Alexa Fluor® 488 anti-human CD123, TotalSeq™-A0064 anti-human CD123, Alexa Fluor® 700 anti-human CD123, APC/Fire™ 750 anti-human CD123, Pacific Blue™ anti-human CD123, TotalSeq™-C0064 anti-human CD123, TotalSeq™-B0064 anti-human CD123, TotalSeq™-D0064 anti-human CD123, PerCP anti-human CD123, GMP PE/Cyanine7 anti-human CD123, APC/Fire™ 810 anti-human CD123 Antibody

Product Data

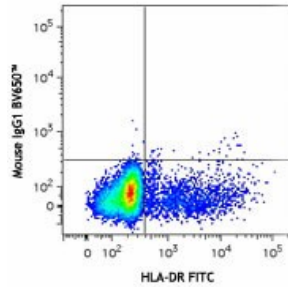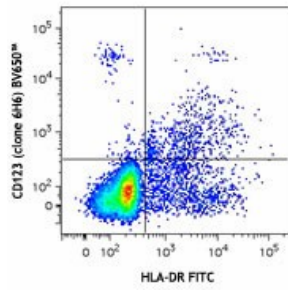

Human peripheral blood leukocytes were stained with HLA-DR FITC and CD123 (clone 6H6) Brilliant Violet 650™ (top) or mouse IgG1 Brilliant Violet 650™ isotype control (bottom).

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
 Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587