

Brilliant Violet 785™ anti-human CD19 Antibody

Catalog# / Size	302239 / 25 tests 302240 / 100 tests
Clone	HIB19
Regulatory Status	RUO
Workshop	V CD19.11
Other Names	B4
Isotype	Mouse IgG1, κ
Description	CD19 is a 95 kD type I transmembrane glycoprotein also known as B4. It is a member of the immunoglobulin superfamily expressed on B-cells (from pro-B to blastoid B cells, absent on plasma cells) and follicular dendritic cells. CD19 is involved in B cell development, activation, and differentiation. CD19 forms a complex with CD21 (CR2) and CD81 (TAPA-1), and functions as a BCR co-receptor.

Product Details

Verified Reactivity	Human
Reported Reactivity	Chimpanzee
Antibody Type	Monoclonal
Host Species	Mouse
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA).
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 785™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration, please enter the lot number in our Concentration and Expiration Lookup or Certificate of Analysis online tools.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	<p>Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 µl per million cells in 100 µl staining volume or 5 µl per 100 µl of whole blood.</p> <p>Brilliant Violet 785™ excites at 405 nm and emits at 785 nm. The bandpass filter 780/60 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 785™ is a trademark of Sirigen Group Ltd.</p> <p>Learn more about Brilliant Violet™.</p> <p>This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.</p>
Excitation Laser	Violet Laser (405 nm)
Application Notes	Additional reported applications (for the relevant formats) include: immunohistochemical staining of acetone-fixed frozen tissue sections ⁸ and blocking of B cell proliferation. Clone HIB19 is not recommended for formalin-fixed paraffin-embedded sections. The Ultra-LEAF™ purified antibody

(Endotoxin < 0.01 EU/μg, Azide-Free, 0.2 μm filtered) is recommended for functional assays (Cat. No. 302267 & 302268).

Clone H1B19 partially blocks anti-human CD19 clones 4G7 and SJ25C1 staining based on in-house testing

Application References

(PubMed link indicates BioLegend citation)

1. Schlossman S, *et al.* 1995. Leucocyte Typing V. Oxford University Press. New York.
2. Knapp W, *et al.* 1989. Leucocyte Typing IV. Oxford University Press. New York.
3. Bradbury L, *et al.* 1993. *J. Immunol.* 151:2915.
4. Joseph A, *et al.* 2010. *J. Virol.* 84:6645. [PubMed](#)
5. Wang X, *et al.* 2010. *Haematologica.* 95:884. (FC) [PubMed](#)
6. Walker JD, *et al.* 2009. *J. Immunol.* 182:1548. (Block) [PubMed](#)
7. Yoshino N, *et al.* 2000. *Exp. Anim. (Tokyo)* 49:97. (FC)
8. Hansen A, *et al.* 2002. *Arthritis Rheum.* 46:2160. (IHC)
9. Stoeckius M, *et al.* 2017. *Nat. Methods.* 14:865. (PG)
10. Peterson VM, *et al.* 2017. *Nat. Biotechnol.* 35:936. (PG)

Product Citations

1. Argüello RJ, *et al.* 2020. *Cell Metab.* 32:1063. [PubMed](#)
2. Tran TM, *et al.* 2020. *Immunity.* 51(4):750-765. [PubMed](#)
3. Ye Y, *et al.* 2022. *Nat Commun.* 13:6458. [PubMed](#)
4. Arce Vargas F *et al.* 2018. *Cancer cell.* 33(4):649-663. [PubMed](#)
5. Meehan S, *et al.* 2019. *Commun Biol.* 2:229. [PubMed](#)
6. Fajgenbaum DC, *et al.* 2019. *J Clin Invest.* 130:4451. [PubMed](#)
7. Lindner JM, *et al.* 2019. *Immunity.* 50:668. [PubMed](#)
8. Tauriainen J, *et al.* 2017. *Sci Rep.* 7:40354. [PubMed](#)
9. Antonucci L, *et al.* 2020. *Journal of Immunology Research.* 3257.929861. [PubMed](#)
10. Laing AG, *et al.* 2020. *Nat Med.* 26:1623. [PubMed](#)
11. Lucas C, *et al.* 2021. *Nature.* Online ahead of print. [PubMed](#)
12. Onodera T, *et al.* 2019. *J Immunol.* 203:3282. [PubMed](#)
13. Costa A *et al.* 2018. *Cancer cell.* 33(3):463-479. [PubMed](#)
14. Maat H, *et al.* 2021. *iScience.* 24(5):102435. [PubMed](#)
15. Vardam-Kaur T, *et al.* 2021. *Oncotarget.* 12:2051. [PubMed](#)
16. Brown CC, *et al.* 2020. *Cell.* 179(4):846-863.e24. [PubMed](#)
17. Vikkurthi R, *et al.* 2022. *Nat Microbiol.* 7:974. [PubMed](#)
18. Combes AJ, *et al.* 2021. *Nature.* 591:124. [PubMed](#)
19. Adachi Y, *et al.* 2019. *Nat Commun.* 10:3883. [PubMed](#)
20. Ruella M, *et al.* 2018. *Nat Med.* 24:1499. [PubMed](#)
21. Chung SH, *et al.* 2021. *Hum Gene Ther.* 32:682. [PubMed](#)
22. Kim ST, *et al.* 2021. *J Immunother Cancer.* 9. [PubMed](#)
23. Garfall AL, *et al.* 2018. *JCI Insight.* 3. [PubMed](#)
24. Wang J, *et al.* 2020. *Cell.* 183(7):1867-1883.e26. [PubMed](#)
25. Kroloff M, *et al.* 2022. *J Immunother Cancer.* 10. [PubMed](#)
26. Friedensohn S, *et al.* 2018. *Front Immunol.* 9:1401. [PubMed](#)

RRID

AB_11218596 (BioLegend Cat. No. 302239)
AB_2563442 (BioLegend Cat. No. 302240)

Antigen Details


Structure	Ig superfamily, type I transmembrane glycoprotein, 95 kD
Distribution	B lineage (except plasma cells), follicular dendritic cells
Function	B cell activation and differentiation
Ligand/Receptor	Forms complex with CD21 (CR2) and CD81 (TAPA-1), BCR coreceptor
Cell Type	B cells, Dendritic cells
Biology Area	Costimulatory Molecules, Immunology
Molecular Family	CD Molecules
Antigen References	1. Tedder T, <i>et al.</i> 1994. <i>Immunol. Today</i> 15:437. 2. Bradbury L, <i>et al.</i> 1993. <i>J. Immunol.</i> 151:2915.
Gene ID	930

Related Protocols

Other Formats

APC anti-human CD19, Biotin anti-human CD19, FITC anti-human CD19, PE anti-human CD19, PE/Cyanine5 anti-human CD19, Purified anti-human CD19, APC/Cyanine7 anti-human CD19, PE/Cyanine7 anti-human CD19, Alexa Fluor® 488 anti-human CD19, Alexa Fluor® 647 anti-human CD19, Pacific Blue™ anti-human CD19, Alexa Fluor® 700 anti-human CD19, PerCP anti-human CD19, PerCP/Cyanine5.5 anti-human CD19, Brilliant Violet 421™ anti-human CD19, Brilliant Violet 570™ anti-human CD19, Brilliant Violet 650™ anti-human CD19, Brilliant Violet 785™ anti-human CD19, Brilliant Violet 510™ anti-human CD19, Brilliant Violet 605™ anti-human CD19, Brilliant Violet 711™ anti-human CD19, Purified anti-human CD19 (Maxpar® Ready), Alexa Fluor® 594 anti-human CD19, PE/Dazzle™ 594 anti-human CD19, APC/Fire™ 750 anti-human CD19, TotalSeq™-A0050 anti-human CD19, Brilliant Violet 750™ anti-human CD19, TotalSeq™-B0050 anti-human CD19, TotalSeq™-C0050 anti-human CD19, Spark NIR™ 685 anti-human CD19, Ultra-LEAF™ Purified anti-human CD19, APC/Fire™ 810 anti-human CD19, PE/Fire™ 640 anti-human CD19, PE/Fire™ 700 anti-human CD19, TotalSeq™-D0050 anti-human CD19, Spark YG™ 593 anti-human CD19, GMP Pacific Blue™ anti-human CD19, Spark Violet™ 423 anti-human CD19, GMP PE anti-human CD19, GMP APC anti-human CD19, KIRAVIA Blue 520™ anti-human CD19, GMP PerCP/Cyanine5.5 anti-human CD19, GMP PE/Cyanine7 anti-human CD19, Spark Violet™ 500 anti-human CD19

Product Data

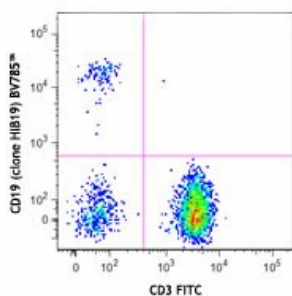

Human peripheral blood lymphocytes were stained with CD3 FITC and CD19 (clone H1B19) Brilliant Violet 785™.

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587