

PE anti-human IFN- γ Antibody

| | |
|--------------------------|---|
| Catalog# / Size | 506506 / 25 tests 506507 / 100 tests |
| Clone | B27 |
| Regulatory Status | RUO |
| Other Names | Interferon- γ , Immune interferon, Type II interferon, T cell interferon, Macrophage-activating factor (MAF), IFN-g, IFN-gamma |
| Isotype | Mouse IgG1, κ |
| Description | Interferon- γ is a potent multifunctional cytokine which is secreted primarily by activated NK cells and T cells. Originally characterized based on anti-viral activities, IFN- γ also exerts anti-proliferative, immunoregulatory, and proinflammatory activities. IFN- γ can upregulate MHC class I and II antigen expression by antigen-presenting cells. The B27 antibody reacts with the human interferon- γ . The B27 antibody can neutralize the bioactivity of natural or recombinant IFN- γ . |

Product Details

| | |
|---|---|
| Verified Reactivity | Human |
| Reported Reactivity | Chimpanzee, Baboon, Cynomolgus, Rhesus, African Green, Pigtailed Macaque, Sooty Mangabey |
| Antibody Type | Monoclonal |
| Host Species | Mouse |
| Immunogen | <i>E. coli</i> -expressed recombinant human IFN- γ |
| Formulation | Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA) |
| Preparation | The antibody was purified by affinity chromatography, and conjugated with PE under optimal conditions. |
| Storage & Handling | The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze. |
| Application | ICFC - Quality tested |
| Recommended Usage | Each lot of this antibody is quality control tested by intracellular immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 μ l per million cells in 100 μ l staining volume or 5 μ l per 100 μ l of whole blood. |
| Excitation Laser | Blue Laser (488 nm) Green Laser (532 nm)/Yellow-Green Laser (561 nm) |
| Application Notes | Flow Cytometry²: The fluorochrome-labeled B27 antibody is useful for intracellular immunofluorescent staining and flow cytometric analysis to identify IFN- γ -producing cells within mixed cell populations. For intracellular cytokine staining protocol, please visit www.biolegend.com and click on the support section. Neutralization^{1,3,6,7}: The Ultra-LEAF™ Purified antibody (Endotoxin <0.01 EU/ μ g, Azide-Free, 0.2 μ m filtered) is recommended for neutralization of human IFN- γ bioactivity (Cat. No. 506531). |
| Application References | <ol style="list-style-type: none"> 1. Favre C, <i>et al.</i> 1989. <i>Molec. Immunol.</i> 26:17. (Neut) 2. Kaur A, <i>et al.</i> 2002. <i>J Virol.</i> 76:3646. 3. Abrams J, <i>et al.</i> 1992. <i>Immunol. Rev.</i> 127:5. (Neut) 4. Andersson U, <i>et al.</i> 1999. <i>Detection and quantification of gene expression</i>. New York:Springer-Verlag. 5. Rout N, <i>et al.</i> 2010. <i>PLoS One</i> 5:e9787. (FC) 6. Acosta-Rodriguez EV, <i>et al.</i> 2007. <i>Nat. Immunol.</i> 9:942-9. (Neut) 7. Gangur V, <i>et al.</i> 1998. <i>FASEB J.</i> 12:705-13. (Neut) |
| (PubMed link indicates BioLegend citation) | |

Product Citations

1. Mo Y, *et al.* 2022. *Front Immunol.* 12:799896. [PubMed](#)
2. Rappazzo CG, *et al.* 2020. *bioRxiv.* . [PubMed](#)
3. Roberts A, *et al.* 2021. *Sci Rep.* 11:4030. [PubMed](#)
4. Bilich T, *et al.* 2021. *Cancer Discov.* 11:1982. [PubMed](#)
5. Ali A, *et al.* 2020. *J Infect Dis.* 222:853. [PubMed](#)
6. Xun X, *et al.* 2021. *Am J Transl Res.* 13:4360. [PubMed](#)
7. Peng Q, *et al.* 2022. *EBioMedicine.* 77:103904. [PubMed](#)
8. Bauer J, *et al.* 2022. *Nat Commun.* 13:6401. [PubMed](#)
9. Das J, *et al.* 2021. *Med (N Y).* 2:1269. [PubMed](#)
10. Lübke M, *et al.* 2020. *J Exp Med.* 217:00:00. [PubMed](#)
11. Qiu C, *et al.* 2022. *Front Immunol.* 12:764949. [PubMed](#)
12. Rice TF, *et al.* 2021. *EBioMedicine.* 72:103612. [PubMed](#)
13. Heitmann JS, *et al.* 2021. *Nature.* Online ahead of print. [PubMed](#)
14. Minassian AM, *et al.* 2021. *Med (N Y).* 2:701. [PubMed](#)
15. Bilich T, *et al.* 2021. *Science Translational Medicine.* 13(590):. [PubMed](#)
16. Punik J, *et al.* 2021. *Cell Reports.* 35(13):109320. [PubMed](#)
17. Such L, *et al.* 2020. *J Clin Invest.* 4266:130. [PubMed](#)
18. Ni J, *et al.* 2020. *Immunity.* 52(6):1075-1087.e8. [PubMed](#)
19. Grace PS, *et al.* 2021. *Front Immunol.* 12:679973. [PubMed](#)
20. Xu C, *et al.* 2020. *Methods Mol Biol.* 2108:159. [PubMed](#)
21. Imura Y, *et al.* 2020. *JCI Insight.* :5. [PubMed](#)
22. Duque JSR, *et al.* 2022. *Nat Commun.* 13:3700. [PubMed](#)
23. Reuschl AK, *et al.* 2022. *Cell Rep.* 39:110650. [PubMed](#)
24. Atyeo C, *et al.* 2022. *Nat Commun.* 13:3571. [PubMed](#)
25. Crump NT, *et al.* 2021. *Cell Reports.* 35(6):109101. [PubMed](#)
26. Zhong M, *et al.* 2022. *Cell Death Dis.* 13:671. [PubMed](#)
27. Luo Z, *et al.* 2021. *Int J Mol Sci.* 22:. [PubMed](#)
28. Niu M, *et al.* 2021. *Cell Reports.* 36(8):109611. [PubMed](#)
29. Jeong HO, *et al.* 2022. *iScience.* 25:105358. [PubMed](#)
30. Li L, *et al.* 2022. *JCI Insight.* 7:. [PubMed](#)
31. Cimen Bozkus C, *et al.* 2021. *STAR Protoc.* 2:100758. [PubMed](#)

RRID AB_315439 (BioLegend Cat. No. 506506)
 AB_315440 (BioLegend Cat. No. 506507)

Antigen Details

| | |
|---------------------------|--|
| Structure | Cytokine; dimer; 20-25 kD (Mammalian) |
| Bioactivity | Antiviral/antiparasitic activities; inhibits proliferation; enhances MHC class I and II expression on APC |
| Cell Sources | CD8 ⁺ and CD4 ⁺ T cells, NK cells |
| Cell Targets | T cells, B cells, macrophages, NK cells, endothelial cells, fibroblasts |
| Receptors | IFN-γRα (CDw119) dimerized with IFN-γRβ (AF-1) |
| Cell Type | Tregs |
| Biology Area | Cell Biology, Immunology, Neuroinflammation, Neuroscience |
| Molecular Family | Cytokines/Chemokines |
| Antigen References | <ol style="list-style-type: none"> 1. Fitzgerald K, <i>et al.</i> Eds. 2001. <i>The Cytokine FactsBook.</i> Academic Press San Diego. 2. De Maeyer E, <i>et al.</i> 1992. <i>Curr. Opin. Immunol.</i> 4:321. 3. Farrar M, <i>et al.</i> 1993. <i>Annu. Rev. Immunol.</i> 11:571. 4. Gray P, <i>et al.</i> 1987. <i>Lymphokines</i> 13:151. |
| Regulation | Upregulated by IL-2, FGF-basic, EGF; downregulated by vitamin D3 or DMN; labile at pH2 |
| Gene ID | 3458 |

Related Protocols


[Surface and Intracellular Cytokine Staining for Flow Cytometry - Video](#)

[Intracellular Flow Cytometry Staining Protocol](#)

Other Formats

APC anti-human IFN- γ , FITC anti-human IFN- γ , PE anti-human IFN- γ , Purified anti-human IFN- γ , Alexa Fluor® 700 anti-human IFN- γ , PE/Cyanine7 anti-human IFN- γ , Purified anti-human IFN- γ (Maxpar® Ready), APC/Cyanine7 anti-human IFN- γ , Pacific Blue™ anti-human IFN- γ , PerCP/Cyanine5.5 anti-human IFN- γ , PE/Dazzle™ 594 anti-human IFN- γ , Ultra-LEAF™ Purified anti-human IFN- γ , Brilliant Violet 605™ anti-human IFN- γ , Brilliant Violet 510™ anti-human IFN- γ , Brilliant Violet 421™ anti-human IFN- γ

Product Data

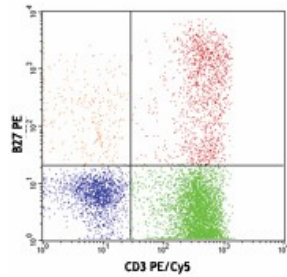

PMA+ionomycin-stimulated (6 hours) human peripheral blood lymphocytes were surface stained with CD3 PE/Cy5 and then intracellularly stained with B27 PE

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587