

## APC/Cyanine7 anti-human CD62L Antibody

<b>Catalog# / Size</b>	304813 / 25 tests 304814 / 100 tests
<b>Clone</b>	DREG-56
<b>Regulatory Status</b>	RUO
<b>Workshop</b>	V S056
<b>Other Names</b>	L-selectin, LECAM-1, LAM-1, Leu-8, TQ-1
<b>Isotype</b>	Mouse IgG1, κ
<b>Description</b>	CD62L is a 74-95 kD single chain type I glycoprotein referred to as L-selectin or LECAM-1. It is expressed on most peripheral blood B cells, subsets of T and NK cells, monocytes, granulocytes, and certain hematopoietic malignant cells. CD62L binds to carbohydrates present on certain glycoforms of CD34, glycam-1, and MAdCAM-1 and with a low affinity to anionic oligosaccharide sequences related to sialylated Lewis X (sLex, CD15s) through its C-type lectin domain. CD62L is important for the homing of naïve lymphocytes to high endothelial venules in peripheral lymph nodes and Peyer's patches. It also plays a role in leukocyte rolling on activated endothelial cells.

### Product Details

---

<b>Verified Reactivity</b>	Human
<b>Reported Reactivity</b>	Chimpanzee, Cow
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Mouse
<b>Immunogen</b>	Concentrated supernatant from PMA-activated human peripheral blood leukocytes
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
<b>Preparation</b>	The antibody was purified by affinity chromatography, and conjugated with APC/Cyanine7 under optimal conditions.
<b>Concentration</b>	Lot-specific (to obtain lot-specific concentration, please enter the lot number in our <a href="#">Concentration and Expiration Lookup</a> or <a href="#">Certificate of Analysis</a> online tools.)
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a> . For flow cytometric staining, the suggested use of this reagent is 5 µl per million cells in 100 µl staining volume or 5 µl per 100 µl of whole blood.
<b>Excitation Laser</b>	Red Laser (633 nm)
<b>Application Notes</b>	Additional reported applications (for the relevant formats) include: Western blotting <sup>2,3,9</sup> and <i>in vitro</i> blocking of lymphocytes binding to high endothelial venules (HEV) <sup>2</sup> . The Ultra-LEAF™ purified antibody (Endotoxin < 0.01 EU/µg, Azide-Free, 0.2 µm filtered) is recommended for functional assays (Cat. Nos. 304853-304858).
<b>Additional Product Notes</b>	BioLegend is in the process of converting the name APC/Cy7 to APC/Cyanine7. The dye molecule remains the same, so you should expect the same quality and performance from our APC/Cyanine7 products. Please contact <a href="#">Technical Service</a> if you have any questions.
<b>Application References</b>	<ol style="list-style-type: none"> <li>Schlossman S, et al. Eds. 1995. Leucocyte Typing V. Oxford University Press. New York.</li> <li>Kishimoto TK, et al. 1990. <i>Proc. Natl. Acad. Sci. USA</i> 87:2244. (WB, Block)</li> <li>Jutila M, et al. 2002. <i>J. Immunol.</i> 169:1768. (WB)</li> </ol>
<b>(PubMed link indicates BioLegend citation)</b>	

4. Tamassia N, *et al.* 2008. *J. Immunol.* 181:6563. (FC) [PubMed](#)
5. Kmiecik M, *et al.* 2009. *J. Transl. Med.* 7:89. (FC) [PubMed](#)
6. Thakral D, *et al.* 2008. *J. Immunol.* 180:7431. (FC) [PubMed](#)
7. Charles N, *et al.* 2010. *Nat. Med.* 16:701. (FC) [PubMed](#)
8. Yoshino N, *et al.* 2000. *Exp. Anim. (Tokyo)* 49:97. (FC)
9. Koenig JM, *et al.* 1996. *Pediatr. Res.* 39:616. (WB)
10. Shi C, *et al.* 2011. *J. Immunol.* 187:5293. (FC) [PubMed](#)
11. Burges M, *et al.* 2013. *Clin Cancer Res.* 19:5675. [PubMed](#)
12. Cash JL, *et al.* 2013. *EMBO Rep.* 14:999. (FC) [PubMed](#)

#### Product Citations

1. Ruck T, *et al.* 2015. *Math Med Biol.* 10.1093/imammb/dqv038. [PubMed](#)
2. Liu L, *et al.* 2016. *Nat Biotechnol.* 10.1038/nbt.3461. [PubMed](#)
3. Ajith A, *et al.* 2021. *Front Immunol.* 12:687715. [PubMed](#)
4. Carre C, *et al.* 2021. *iScience.* 24:102970. [PubMed](#)
5. Wei&Bszig;müller S, *et al.* 2016. *PLoS One.* 11: 0149093. [PubMed](#)
6. Björklund &, *et al.* 2016. *Nat Immunol.* 17:451-460. [PubMed](#)
7. Caduff N, *et al.* 2021. *Cell Reports.* 35(5):109056. [PubMed](#)
8. Piekarska K, *et al.* 2022. *Nat Commun.* 13:856. [PubMed](#)
9. Amarnath S, *et al.* 2010. *PLoS One.* 5:e1000302. [PubMed](#)
10. Hunter S, *et al.* 2018. *J Hepatol.* 69:654. [PubMed](#)
11. Dettling S, *et al.* 2018. *Clin Cancer Res.* 24(12). [PubMed](#)
12. Engelmann C, *et al.* 2016. *Sci Rep.* 6:37926. [PubMed](#)
13. Koh WH, *et al.* 2020. *STAR Protoc.* 1:100203. [PubMed](#)
14. Kim ML, *et al.* 2021. *iScience.* 24:103509. [PubMed](#)
15. Choe JH, *et al.* 2021. *Sci Transl Med.* 13:. [PubMed](#)

#### RRID

AB\_493583 (BioLegend Cat. No. 304813)  
 AB\_493582 (BioLegend Cat. No. 304814)

## Antigen Details

---

<b>Structure</b>	Selectin, single chain glycoprotein, 74-95 kD
<b>Distribution</b>	Majority of B cells, naïve T cells, subset of memory T and NK cells, monocytes, granulocytes, thymocytes
<b>Function</b>	Leukocyte homing, leukocyte tethering, rolling
<b>Ligand/Receptor</b>	CD34, GlyCAM, MAdCAM-1
<b>Cell Type</b>	B cells, Granulocytes, Monocytes, Neutrophils, NK cells, T cells, Thymocytes, Tregs
<b>Biology Area</b>	Cell Adhesion, Cell Biology, Costimulatory Molecules, Immunology, Innate Immunity
<b>Molecular Family</b>	Adhesion Molecules, CD Molecules
<b>Antigen References</b>	<ol style="list-style-type: none"> <li>1. Kishimoto T, <i>et al.</i> 1990. <i>P. Natl. Acad. Sci. USA</i> 87:2244.</li> <li>2. Kishimoto T, <i>et al.</i> 1991. <i>Blood</i> 78:805.</li> </ol>

#### Gene ID

[6402](#)

## Related Protocols

---

[Cell Surface Flow Cytometry Staining Protocol](#)


## Other Formats

---

APC anti-human CD62L, FITC anti-human CD62L, PE anti-human CD62L, PE/Cyanine5 anti-human CD62L, Purified anti-human CD62L, APC/Cyanine7 anti-human CD62L, Alexa Fluor® 488 anti-human CD62L, Alexa Fluor® 647 anti-human CD62L, Alexa Fluor® 700 anti-human CD62L, PE/Cyanine7 anti-human CD62L, PerCP/Cyanine5.5 anti-human CD62L, Pacific Blue™ anti-human CD62L, Brilliant Violet 421™ anti-human CD62L, Brilliant Violet 785™ anti-human CD62L, Brilliant Violet 650™ anti-human CD62L, PE/Dazzle™ 594 anti-human CD62L, Brilliant Violet 605™ anti-human CD62L, Purified anti-human CD62L (Maxpar® Ready), APC/Fire™ 750 anti-human CD62L, Brilliant Violet 510™ anti-human CD62L, TotalSeq™-A0147 anti-human CD62L, TotalSeq™-B0147 anti-human CD62L, TotalSeq™-C0147 anti-human CD62L, Ultra-LEAF™ Purified anti-human CD62L, Brilliant Violet 711™ anti-human CD62L, Spark NIR™ 685 anti-human CD62L, TotalSeq™-D0147 anti-human CD62L, APC/Fire™ 810 anti-human

## Product Data

---

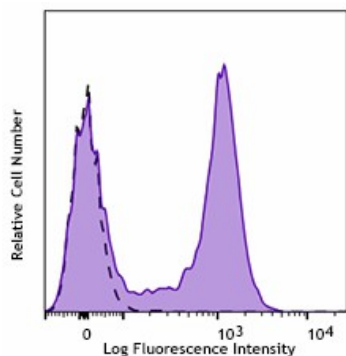

Human peripheral blood lymphocytes were stained with CD62L (clone DREG-56) APC/Cyanine7 (filled histogram) or Mouse IgG1, ? APC/Cyanine7 isotype control (open histogram).

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

\*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, [www.biolegend.com/ordering#license](http://www.biolegend.com/ordering#license)). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587