

Alexa Fluor[®] 594 anti-mouse CD11c Antibody

Catalog# / Size	117346 / 100 µg
Clone	N418
Regulatory Status	RUO
Other Names	αX integrin, integrin αX chain, CR4, p150, ITGAX
Isotype	Armenian Hamster IgG
Description	CD11c is a 150 kD glycoprotein also known as α _X integrin, CR4, and p150. CD11c forms a α _X β ₂ heterodimer with β ₂ integrin (CD18). It is primarily expressed on dendritic cells, NK cells, a subset of intestinal intraepithelial lymphocytes (IEL), and some activated T cells. The α _X β ₂ integrin plays an important role in cell-cell contact by binding its ligands: iC3b, fibrinogen, and CD54.

Product Details

Verified Reactivity	Mouse
Antibody Type	Monoclonal
Host Species	Armenian Hamster
Immunogen	Mouse spleen dendritic cells
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide.
Preparation	The antibody was purified by affinity chromatography and conjugated with Alexa Fluor [®] 594 under optimal conditions.
Concentration	0.5 mg/mL
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	IHC-F - Quality tested 3D IHC- Verified
Recommended Usage	<p>Each lot of this antibody is quality control tested by immunohistochemical staining on frozen tissue sections. For immunohistochemistry, a concentration range of 2.0 - 5.0 µg/mL is suggested. For 3D immunohistochemistry on formalin-fixed tissues, a concentration of 5.0 µg/mL is suggested. It is recommended that the reagent be titrated for optimal performance for each application.</p> <p>* Alexa Fluor[®] 594 has an excitation maximum of 590 nm, and a maximum emission of 617 nm.</p> <p>Alexa Fluor[®] and Pacific Blue™ are trademarks of Life Technologies Corporation.</p> <p>View full statement regarding label licenses</p>
Excitation Laser	Green Laser (532 nm)/Yellow-Green Laser (561 nm)
Application Notes	Additional reported applications (for the relevant formats) include: immunoprecipitation ³ , immunohistochemical staining of acetone-fixed frozen sections ³ , immunofluorescence microscopy ^{5,9} (Alexa Fluor [®] 488 conjugated N418 was used for IHC in frozen sections ¹⁰), and spatial biology (IBEX) ^{22,23} .
Application References	<ol style="list-style-type: none"> Granucci F, <i>et al.</i> 1997. <i>J. Immunol.</i> 159:1794. Stokes RW, <i>et al.</i> 1998. <i>J. Immunol.</i> 160:5514. Metlay JP, <i>et al.</i> 1990. <i>J. Exp. Med.</i> 171:1753. (IHC, IP) Ma XT, <i>et al.</i> 2006. <i>Cancer Research</i> 66:1169. Chin RK, <i>et al.</i> 2006. <i>J. Immunol.</i> 177:290. (IF) Cervantes-Barragan L, <i>et al.</i> 2007. <i>Blood</i> 109:1131. (FC) PubMed Turnquist HR, <i>et al.</i> 2007. <i>J. Immunol.</i> 178:7018. (FC) PubMed Benson MJ, <i>et al.</i> 2007. <i>J. Exp. Med.</i> doi:10.1084/jem.20070719. (FC) PubMed
(PubMed link indicates BioLegend citation)	

9. You Y, *et al.* 2009. *J. Immunol.* 182:7343. (IF) [PubMed](#)
10. Roland CL, *et al.* 2009. *Mol. Cancer Res.* 8:1761. (IHC, FC) [PubMed](#)
11. Wikstrom M, *et al.* 2006. *J. Immunol.* 177:913. [PubMed](#)
12. Pericolini E, *et al.* 2008. *J. Leukocyte Biol.* 83:1286. [PubMed](#)
13. Randall LM, *et al.* 2008. *Infect. Immun.* 76:3312. [PubMed](#)
14. Fahlen-Yrild L, *et al.* 2009. *J. Immunol.* 183:5032. [PubMed](#)
15. Osterholzer JJ, *et al.* 2009. *J. Immunol.* 183:8044. [PubMed](#)
16. Bankoti J, *et al.* 2010. *Toxicol. Sci.* 115:422. (FC) [PubMed](#)
17. Eisenach PA, *et al.* 2010. *J Cell Sci.* 123:4182. [PubMed](#)
18. Leppin K, *et al.* 2014. *Invest. Ophthalmol. Vis. Sci.* 55:3603. [PubMed](#)
19. Sakai F, *et al.* 2014. *PLoS One.* 9:105370. [PubMed](#)
20. Gibbins JD, *et al.* 2014. *Blood.* 124:2953. [PubMed](#)
21. White CE, *et al.* 2015. *J Immunol.* 194:697. [PubMed](#)
22. Lu X, *et al.* 2015. *J Immunol.* 194:2011. [PubMed](#)
23. Radtke AJ, *et al.* 2020. *Proc Natl Acad Sci U S A.* 117:33455-65. (SB) [PubMed](#)
24. Radtke AJ, *et al.* 2022. *Nat Protoc.* 17:378-401. (SB) [PubMed](#)

Product Citations

1. Baptista AP *et al.* 2019. *Immunity.* 50(5):1188-1201 . [PubMed](#)
2. Bittner-Eddy P, *et al.* 2016. *J Immunol.* 197: 1435 - 1446. [PubMed](#)
3. Haugh KA, *et al.* 2021. *eLife.* 0.416666666666667. [PubMed](#)
4. Sokol CL *et al.* 2018. *Immunity.* 49(3):449-463 . [PubMed](#)
5. Samarchith P Kurup *et al.* 2019. *Cell host & microbe.* 25(4):565-577 . [PubMed](#)
6. Adeegbe DO, *et al.* 2017. *Cancer Discov.* 0.883333333. [PubMed](#)
7. De Pascalis R, *et al.* 2020. *PLoS One.* 15:e0237034. [PubMed](#)
8. Hofmann J, *et al.* 2021. *Front Immunol.* 11:599495. [PubMed](#)

RRID

AB_2563323 (BioLegend Cat. No. 117346)

Antigen Details

Structure	Integrin α -chain, associates with integrin β 2 (CD18), 150 kD
Distribution	Dendritic cells, NK cells, intestinal intraepithelial lymphocytes (IEL), some activated T cells
Function	Cellular adhesion
Ligand/Receptor	iC3b, fibrinogen
Cell Type	Dendritic cells, Epithelial cells, NK cells, T cells, Tregs
Biology Area	Cell Adhesion, Cell Biology, Costimulatory Molecules, Immunology, Innate Immunity, Neuroscience, Neuroscience Cell Markers
Molecular Family	Adhesion Molecules, CD Molecules
Antigen References	<ol style="list-style-type: none"> 1. Barclay A, <i>et al.</i> 1997. <i>The Leukocyte Antigen Facts Book</i> Academic Press. 2. Springer TA. 1994. <i>Cell</i> 76:301. 3. Lopez-Rodriguez C, <i>et al.</i> 1996. <i>J. Immunol.</i> 156:3780.
Gene ID	16411

Related Protocols

[Immunohistochemistry Protocol for Frozen Sections](#)


[Ce3D™ Tissue Clearing Kit](#)

Other Formats


APC anti-mouse CD11c, Biotin anti-mouse CD11c, FITC anti-mouse CD11c, PE anti-mouse CD11c, Purified anti-mouse CD11c, Alexa Fluor® 488 anti-mouse CD11c, Alexa Fluor® 647 anti-mouse CD11c, PE/Cyanine5 anti-mouse CD11c, PE/Cyanine7 anti-mouse CD11c, Brilliant Violet 605™ anti-mouse CD11c, Alexa Fluor® 700 anti-mouse CD11c, Pacific Blue™ anti-mouse CD11c, APC/Cyanine7 anti-mouse CD11c, PerCP/Cyanine5.5 anti-mouse CD11c, PerCP anti-mouse CD11c, Brilliant Violet 421™ anti-mouse CD11c, Brilliant Violet 570™ anti-mouse CD11c, Brilliant Violet 785™ anti-mouse CD11c, Brilliant Violet 510™ anti-mouse CD11c, Brilliant Violet 650™ anti-mouse CD11c, Purified anti-mouse CD11c (Maxpar® Ready), Alexa Fluor® 594 anti-mouse CD11c, PE/Dazzle™ 594 anti-mouse CD11c, Brilliant Violet 711™ anti-mouse CD11c, APC/Fire™ 750 anti-mouse CD11c, TotalSeq™-A0106 anti-mouse CD11c, Brilliant Violet 750™ anti-mouse CD11c, TotalSeq™-B0106 anti-mouse CD11c, TotalSeq™-C0106 anti-

mouse CD11c, KIRAVIA Blue 520™ anti-mouse CD11c, Spark Blue™ 550 anti-mouse CD11c, Spark NIR™ 685 anti-mouse CD11c, Spark UV™ 387 anti-mouse CD11c, Spark Red™ 718 anti-mouse CD11c

Product Data

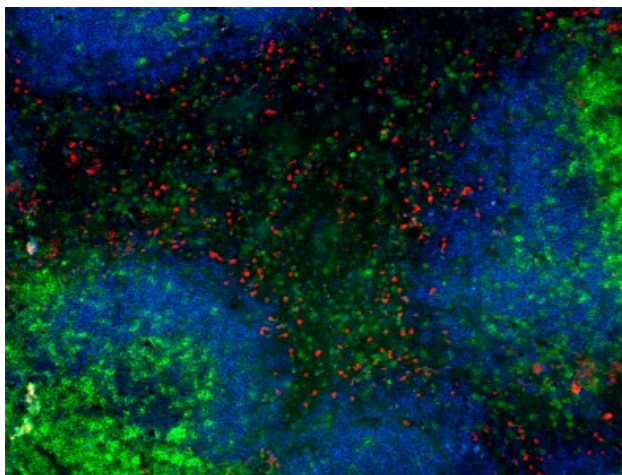

C57BL/6 mouse frozen spleen section was fixed with 4% paraformaldehyde (PFA) for 10 minutes at room temperature and blocked with 5% FBS plus 5% rat serum for 1 hour at room temperature. Then the section was stained with 5 µg/ml of CD11c (clone N418) Alexa Fluor® 594 (red), 2.5 µg/ml of CD3 (clone 145-2C11) Alexa Fluor® 647 (green), and 2.5 µg/ml of B220 (clone RA3-6B2) Alexa Fluor® 488 (blue) overnight at 4°C. The image was captured by 10X objective.

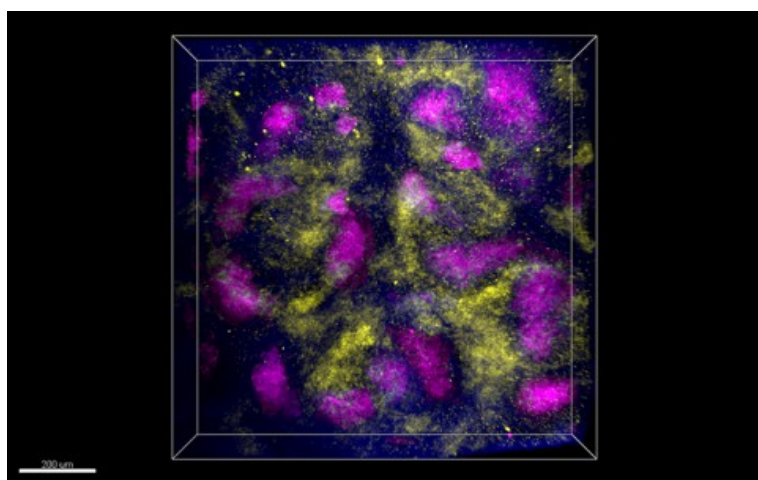

Paraformaldehyde-fixed (4%), 500 µm-thick mouse spleen section was processed according to the Ce3DTM Tissue Clearing Kit protocol (cat. no. 427701). The section was costained with anti-mouse CD11c Antibody (clone N418) Alexa Fluor® 594 at 5 µg/mL (yellow), and anti-mouse CD21/CD35 (CR2/CR1) Antibody (Clone 7E9) Alexa Fluor® 647 at 5 µg/mL (magenta) and counterstained with DAPI (blue). The section was then optically cleared and mounted in a sample chamber. The image was captured with a 10X objective using Zeiss 780 confocal microscope and processed by Imaris image analysis software.

[Watch the video.](#)

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587