

Alexa Fluor[®] 700 anti-human CD27 Antibody

Catalog# / Size	302813 / 25 µg 302814 / 100 µg
Clone	O323
Regulatory Status	RUO
Workshop	IV T-186
Other Names	S152, T14, TNFRSF7
Isotype	Mouse IgG1, κ
Description	CD27 is a 50-55 kD type I membrane protein also known as S152 and T14. It is a lymphocyte-specific member of the TNF-receptor superfamily. CD27 is expressed on medullary thymocytes, virtually all mature T cells, some B cells, and NK cells. CD27 binds to CD70 and plays an important role in costimulation of T cell activation, and regulation of B cell differentiation and proliferation. The cytoplasmic domains of CD27 have also been shown to interact with TRAF2 and TRAF5 to elicit NF-κB and SAPK/JNK activation.

Product Details

Verified Reactivity	Human, Cynomolgus, Rhesus
Reported Reactivity	African Green, Baboon, Squirrel Monkey
Antibody Type	Monoclonal
Host Species	Mouse
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide.
Preparation	The antibody was purified by affinity chromatography and conjugated with Alexa Fluor [®] 700 under optimal conditions.
Concentration	0.5 mg/ml
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	<p>Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. The suggested use of this reagent is ≤1.0 µg per million cells in 100 µl volume. It is highly recommended that the reagent be titrated for optimal performance for each application.</p> <p>* Alexa Fluor[®] 700 has a maximum emission of 719 nm when it is excited at 633 nm / 635 nm. Prior to using Alexa Fluor[®] 700 conjugate for flow cytometric analysis, please verify your flow cytometer's capability of exciting and detecting the fluorochrome.</p> <p>Alexa Fluor[®] and Pacific Blue™ are trademarks of Life Technologies Corporation.</p> <p>View full statement regarding label licenses</p>
Excitation Laser	Red Laser (633 nm)
Application References	1. Knapp W, <i>et al.</i> Eds. 1989. Leucocyte Typing IV. Oxford University Press. New York. 2. Correia DV, <i>et al.</i> 2011. <i>Blood</i> 118:992. (FC) PubMed
(PubMed link indicates BioLegend citation)	
Product Citations	1. Ravens S, <i>et al.</i> 2020. Proc Natl Acad Sci U S A. 117:18649. PubMed 2. Sauce D, <i>et al.</i> 2012. J Immunol. 189:5541. PubMed 3. Esaulova E, <i>et al.</i> 2020. Cell Host Microbe. 29(2):165-178.e8. PubMed 4. Eller MA, <i>et al.</i> 2016. J Virol. 90: 4005-4016. PubMed

5. Fali T, *et al.* 2019. *J Immunol.* 202:1962. [PubMed](#)
6. Onodera T, *et al.* 2019. *J Immunol.* 203:3282. [PubMed](#)
7. Sayin I, *et al.* 2018. *J Exp Med.* 7:40286. [PubMed](#)
8. Buggert M, *et al.* 2020. *Cell.* 183(7):1946-1961.e15. [PubMed](#)
9. Wang X, *et al.* 2019. *J Immunol.* 203:2171. [PubMed](#)
10. Cherukuri A, *et al.* 2014. *J Am Soc Nephrol.* 7:1575. [PubMed](#)
11. Rouers A, *et al.* 2021. *iScience.* 24:102482. [PubMed](#)
12. Adachi Y, *et al.* 2019. *Nat Commun.* 10:3883. [PubMed](#)
13. NULL, *et al.* 2022. *Cell.* 185:916. [PubMed](#)
14. Almeida J, *et al.* 2007. *J Exp Med.* 204:2473. [PubMed](#)
15. Papagno L, *et al.* 2020. *Sci Rep.* 10:11620. [PubMed](#)
16. Sauce D, *et al.* 2011. *Blood.* 117:5142. [PubMed](#)

RRID AB_493756 (BioLegend Cat. No. 302813)
 AB_493757 (BioLegend Cat. No. 302814)

Antigen Details

Structure	TNF-R superfamily, type I transmembrane glycoprotein, 50-55 kD
Distribution	Medullary thymocytes, T and B cell subsets, NK cells
Function	T cell costimulation
Ligand/Receptor	CD70
Cell Type	B cells, NK cells, T cells, Thymocytes, Tregs
Biology Area	Costimulatory Molecules, Immunology
Molecular Family	CD Molecules
Antigen References	1. Hintzen R, <i>et al.</i> 1994. <i>Immunol. Today</i> 15:307. 2. Agematsu K, <i>et al.</i> 1995. <i>J. Immunol.</i> 154:3627.
Gene ID	939


Related Protocols

[Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

APC anti-human CD27, Biotin anti-human CD27, FITC anti-human CD27, PE anti-human CD27, Purified anti-human CD27, Alexa Fluor® 647 anti-human CD27, Alexa Fluor® 700 anti-human CD27, APC/Cyanine7 anti-human CD27, PerCP anti-human CD27, PerCP/Cyanine5.5 anti-human CD27, Pacific Blue™ anti-human CD27, Brilliant Violet 421™ anti-human CD27, Brilliant Violet 570™ anti-human CD27, Brilliant Violet 650™ anti-human CD27, Brilliant Violet 605™ anti-human CD27, Brilliant Violet 711™ anti-human CD27, Brilliant Violet 785™ anti-human CD27, Brilliant Violet 510™ anti-human CD27, PE/Cyanine7 anti-human CD27, Purified anti-human CD27 (Maxpar® Ready), PE/Dazzle™ 594 anti-human CD27, APC/Fire™ 750 anti-human CD27, TotalSeq™-A0154 anti-human CD27, Brilliant Violet 750™ anti-human CD27, TotalSeq™-B0154 anti-human CD27, TotalSeq™-C0154 anti-human CD27, Spark NIR™ 685 anti-human CD27, PE/Fire™ 810 anti-human CD27, TotalSeq™-D0154 anti-human CD27, APC/Fire™ 810 anti-human CD27, PE/Cyanine5 anti-human CD27 Antibody, Spark UV™ 387 anti-human CD27

Product Data

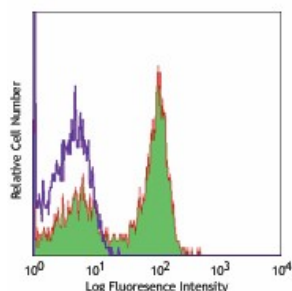

Human peripheral blood lymphocytes stained with O323 Alexa Fluor® 700

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587