

## Alexa Fluor<sup>®</sup> 700 anti-mouse CD86 Antibody

<b>Catalog# / Size</b>	105023 / 25 µg 105024 / 100 µg
<b>Clone</b>	GL-1
<b>Regulatory Status</b>	RUO
<b>Other Names</b>	B7-2, B70, Ly-58
<b>Isotype</b>	Rat IgG2a, κ
<b>Description</b>	CD86 is an 80 kD immunoglobulin superfamily member also known as B7-2, B70, and Ly-58. CD86 is expressed on activated B and T cells, macrophages, dendritic cells, and astrocytes. CD86, along with CD80, is a ligand of CD28 and CD152 (CTLA-4). CD86 is expressed earlier in the immune response than CD80. CD86 has also been shown to be involved in immunoglobulin class-switching and triggering of NK cell-mediated cytotoxicity. CD86 binds to CD28 to transduce co-stimulatory signals for T cell activation, proliferation, and cytokine production. CD86 can also bind to CD152, also known as CTLA-4, to deliver an inhibitory signal to T cells.

### Product Details

<b>Verified Reactivity</b>	Mouse
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Rat
<b>Immunogen</b>	LPS-activated CBA/Ca mouse splenic B cells
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide.
<b>Preparation</b>	The antibody was purified by affinity chromatography and conjugated with Alexa Fluor <sup>®</sup> 700 under optimal conditions.
<b>Concentration</b>	0.5 mg/ml
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	<p>Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a>. The suggested use of this reagent is ≤1.0 µg per million cells in 100 µl volume. It is highly recommended that the reagent be titrated for optimal performance for each application.</p> <p>* Alexa Fluor<sup>®</sup> 700 has a maximum emission of 719 nm when it is excited at 633 nm / 635 nm. Prior to using Alexa Fluor<sup>®</sup> 700 conjugate for flow cytometric analysis, please verify your flow cytometer's capability of exciting and detecting the fluorochrome.</p> <p>Alexa Fluor<sup>®</sup> and Pacific Blue™ are trademarks of Life Technologies Corporation.</p> <p><a href="#">View full statement regarding label licenses</a></p>
<b>Excitation Laser</b>	Red Laser (633 nm)
<b>Application Notes</b>	The GL-1 antibody can block the mixed lymphocyte reaction <i>in vitro</i> and has been shown to inhibit the priming of cytotoxic T lymphocytes <i>in vivo</i> (along with antibodies against B7-1). Additional reported applications (for the relevant formats) include: immunoprecipitation <sup>1</sup> , immunohistochemical staining of acetone-fixed frozen sections <sup>2,6</sup> , immunofluorescence microscopy, and <i>in vivo</i> and <i>in vitro</i> blocking of T cell responses <sup>1-6</sup> . GL-1 is not suitable for immunohistochemical staining of formalin-fixed paraffin sections. The Ultra-LEAF™ purified antibody (Endotoxin < 0.01 EU/µg, Azide-Free, 0.2 µm filtered) is recommended for functional assays (Cat. No. 105051-105056).
<b>Application References</b>	

**(PubMed link indicates BioLegend citation)**

1. Hathcock KS, *et al.* 1993. *Science* 262:905. (Block, IP)
2. Inaba KM, *et al.* 1994. *J. Exp. Med.* 180:1849. (Block, IHC)
3. Hathcock KS, *et al.* 1994. *J. Exp. Med.* 180:631. (Block)
4. Krummel MF, *et al.* 1995. *J. Exp. Med.* 182:459. (Block)
5. Liu Y, *et al.* 1997. *J. Exp. Med.* 185:251. (Block)
6. Herold KC, *et al.* 1997. *J. Immunol.* 158:984. (Block, IHC)
7. Shih FF, *et al.* 2006. *J. Immunol.* 176:3438. (FC)
8. Lawson BR, *et al.* 2007. *J. Immunol.* 178:5366.
9. Turnquist HR, *et al.* 2007. *J. Immunol.* 178:7018.
10. Klinger MB, *et al.* 2007. *Am. J. Physiol. Regul. Integr. Comp. Physiol.* 293:R677. [PubMed](#)
11. de Verteuil DA, *et al.* 2014. *J Immunol.* 193:1121. [PubMed](#)

**Product Citations**

1. Dando SJ, *et al.* 2021. *Invest Ophthalmol Vis Sci.* 62:10. [PubMed](#)
2. Enriquez AB, *et al.* 2022. *iScience.* 25:104305. [PubMed](#)
3. Swain S, *et al.* 2011. *Infect Immun.* 79:1905. [PubMed](#)
4. Enriquez AB, *et al.* 2022. *iScience.* 25:104305. [PubMed](#)
5. Gubser C, *et al.* 2016. *Sci Rep.* 6: 25758. [PubMed](#)
6. Chen L, *et al.* 2020. *J Exp Med.* 217:00:00. [PubMed](#)
7. Leylek R, *et al.* 2019. *Cell Rep.* 29:3736. [PubMed](#)
8. Kumar MP, *et al.* 2018. *Cell Rep.* 25:1458. [PubMed](#)
9. Roberta Zappasodi *et al.* 2018. *Cancer cell.* 33(6):1017-1032 . [PubMed](#)
10. Campbell C *et al.* 2018. *Immunity.* 48(6):1245-1257 . [PubMed](#)
11. Chung L *et al.* 2018. *Cell host & microbe.* 23(2):203-214 . [PubMed](#)
12. Walker MT, *et al.* 2021. *Front. Allergy.* 2:. [PubMed](#)
13. Molina MS, *et al.* 2021. *Front Immunol.* 12:699128. [PubMed](#)
14. Clemente-Casares X, *et al.* 2017. *Immunity.* 47:974. [PubMed](#)
15. Sadtler K, *et al.* 2017. *Tissue Eng Part A.* 23:1044. [PubMed](#)

**RRID**

AB\_493720 (BioLegend Cat. No. 105023)  
AB\_493721 (BioLegend Cat. No. 105024)

**Antigen Details**

---

<b>Structure</b>	Ig superfamily, 80 kD
<b>Distribution</b>	B cells and T cells (upregulated upon activation), macrophages, dendritic cells, and astrocytes
<b>Function</b>	T cell costimulation, Ig class-switching, NK cell cytotoxicity
<b>Ligand/Receptor</b>	CD28, CD152 (CTLA-4)
<b>Cell Type</b>	Astrocytes, B cells, Dendritic cells, Macrophages, T cells, Tregs
<b>Biology Area</b>	Cell Biology, Costimulatory Molecules, Immunology, Neuroscience, Neuroscience Cell Markers
<b>Molecular Family</b>	CD Molecules, Immune Checkpoint Receptors
<b>Antigen References</b>	<ol style="list-style-type: none"><li>1. Barclay A, <i>et al.</i> 1997. <i>The Leukocyte Antigen FactsBook</i> Academic Press.</li><li>2. Hathcock KS, <i>et al.</i> 1993. <i>Science</i> 262:905.</li><li>3. Freeman GJ, <i>et al.</i> 1993. <i>Science</i> 262:907.</li><li>4. Carreno BM, <i>et al.</i> 2002. <i>Annu. Rev. Immunol.</i> 20:29.</li></ol>
<b>Gene ID</b>	<a href="#">12524</a>

**Related Protocols**

---

[Cell Surface Flow Cytometry Staining Protocol](#)


**Other Formats**

---

Brilliant Violet 650™ anti-mouse CD86, Biotin anti-mouse CD86, FITC anti-mouse CD86, PE anti-mouse CD86, Purified anti-mouse CD86, Brilliant Violet 605™ anti-mouse CD86, APC anti-mouse CD86, PE/Cyanine7 anti-mouse CD86, Alexa Fluor® 488 anti-mouse CD86, Alexa Fluor® 647 anti-mouse CD86, Pacific Blue™ anti-mouse CD86, PE/Cyanine5 anti-mouse CD86, Alexa Fluor® 700 anti-mouse CD86, PerCP/Cyanine5.5 anti-mouse CD86, PerCP anti-mouse CD86, APC/Cyanine7 anti-mouse CD86, Brilliant Violet 421™ anti-mouse CD86, Brilliant Violet 510™ anti-mouse CD86, PE/Dazzle™ 594 anti-mouse CD86, Brilliant Violet 785™ anti-mouse CD86, APC/Fire™ 750 anti-mouse CD86, TotalSeq™-A0200 anti-mouse CD86, TotalSeq™-B0200 anti-mouse CD86, Ultra-LEAF™ Purified anti-mouse CD86, TotalSeq™-C0200 anti-mouse CD86

## Product Data

---

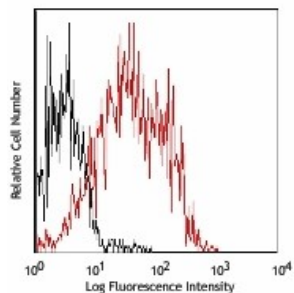

LPS-stimulated (3 days) C57BL/6 mouse splenocytes stained with GL-1 Alexa Fluor® 700

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

\*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, [www.biolegend.com/ordering#license](http://www.biolegend.com/ordering#license)). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587