

Brilliant Violet 510™ anti-mouse CD11c Antibody

Catalog# / Size	117337 / 125 µL 117353 / 50 µg 117338 / 500 µL
Clone	N418
Regulatory Status	RUO
Other Names	αX integrin, integrin αX chain, CR4, p150, ITGAX
Isotype	Armenian Hamster IgG
Description	CD11c is a 150 kD glycoprotein also known as αX integrin, CR4, and p150. CD11c forms a αXβ2 heterodimer with β2 integrin (CD18). It is primarily expressed on dendritic cells, NK cells, a subset of intestinal intraepithelial lymphocytes (IEL), and some activated T cells. The αXβ2 integrin plays an important role in cell-cell contact by binding its ligands: iC3b, fibrinogen, and CD54.

Product Details

Verified Reactivity	Mouse
Antibody Type	Monoclonal
Host Species	Armenian Hamster
Immunogen	Mouse spleen dendritic cells
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA).
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 510™ under optimal conditions.
Concentration	µg size: 0.2 mg/mL µL size: lot-specific (to obtain lot-specific concentration, please enter the lot number in our Concentration and Expiration Lookup or Certificate of Analysis online tools.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	<p>Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining using the µL size, the suggested use of this reagent is 5 µl per million cells in 100 µl staining volume or 5 µl per 100 µl of whole blood. For flow cytometric staining using the µg size, the suggested use of this reagent is ≤0.25 µg per million cells in 100 µl volume. It is recommended that the reagent be titrated for optimal performance for each application.</p> <p>Brilliant Violet 510™ excites at 405 nm and emits at 510 nm. The bandpass filter 510/50 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 510™ is a trademark of Sirigen Group Ltd.</p> <p>Learn more about Brilliant Violet™.</p> <p>This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.</p>
Excitation Laser	Violet Laser (405 nm)
Application Notes	Additional reported applications (for the relevant formats) include: immunoprecipitation ³ ,

immunohistochemical staining of acetone-fixed frozen sections³, immunofluorescence microscopy⁵,
⁹ (Alexa Fluor® 488 conjugated N418 was used for IHC in frozen sections¹⁰), and spatial biology
(IBEX)^{22,23}.

Application References

(PubMed link indicates
BioLegend citation)

1. Granucci F, *et al.* 1997. *J. Immunol.* 159:1794.
2. Stokes RW, *et al.* 1998. *J. Immunol.* 160:5514.
3. Metlay JP, *et al.* 1990. *J. Exp. Med.* 171:1753. (IHC, IP)
4. Ma XT, *et al.* 2006. *Cancer Research* 66:1169.
5. Chin RK, *et al.* 2006. *J. Immunol.* 177:290. (IF)
6. Cervantes-Barragan L, *et al.* 2007. *Blood* 109:1131. (FC) [PubMed](#)
7. Turnquist HR, *et al.* 2007. *J. Immunol.* 178:7018. (FC) [PubMed](#)
8. Benson MJ, *et al.* 2007. *J. Exp. Med.* doi:10.1084/jem.20070719. (FC) [PubMed](#)
9. You Y, *et al.* 2009. *J. Immunol.* 182:7343. (IF) [PubMed](#)
10. Roland CL, *et al.* 2009. *Mol. Cancer Res.* 8:1761. (IHC, FC) [PubMed](#)
11. Wikstrom M, *et al.* 2006. *J. Immunol.* 177:913. [PubMed](#)
12. Pericolini E, *et al.* 2008. *J. Leukocyte Biol.* 83:1286. [PubMed](#)
13. Randall LM, *et al.* 2008. *Infect. Immun.* 76:3312. [PubMed](#)
14. Fahlen-Yrild L, *et al.* 2009. *J. Immunol.* 183:5032. [PubMed](#)
15. Osterholzer JJ, *et al.* 2009. *J. Immunol.* 183:8044. [PubMed](#)
16. Bankoti J, *et al.* 2010. *Toxicol. Sci.* 115:422. (FC) [PubMed](#)
17. Eisenach PA, *et al.* 2010. *J Cell Sci.* 123:4182. [PubMed](#)
18. Leppin K, *et al.* 2014. *Invest. Ophthalmol. Vis. Sci.* 55:3603. [PubMed](#)
19. Sakai F, *et al.* 2014. *PLoS One.* 9:105370. [PubMed](#)
20. Gibbins JD, *et al.* 2014. *Blood.* 124:2953. [PubMed](#)
21. White CE, *et al.* 2015. *J Immunol.* 194:697. [PubMed](#)
22. Lu X, *et al.* 2015. *J Immunol.* 194:2011. [PubMed](#)
23. Radtke AJ, *et al.* 2020. *Proc Natl Acad Sci U S A.* 117:33455-65. (SB) [PubMed](#)
24. Radtke AJ, *et al.* 2022. *Nat Protoc.* 17:378-401. (SB) [PubMed](#)

Product Citations

1. Marchelletta RR, *et al.* 2021. *J Clin Invest.* 131:. [PubMed](#)
2. Yao RQ, *et al.* 2022. *Theranostics.* 12:4606. [PubMed](#)
3. Ramakrishna C, *et al.* 2019. *Nat Commun.* 10:2153. [PubMed](#)
4. Ferretti A, *et al.* 2016. *J Immunol.* 197: 470 - 479. [PubMed](#)
5. Barry KC, *et al.* 2018. *Nat Med.* 24:1178. [PubMed](#)
6. Larsen SE, *et al.* 2021. *Sci Rep.* 11:9040. [PubMed](#)
7. Lou F, *et al.* 2020. *STAR Protoc.* 1:100115. [PubMed](#)
8. Zhu B, *et al.* 2021. *Immunity.* 54(6):1200-1218.e9. [PubMed](#)
9. Walls AC, *et al.* 2021. *Cell.* 184:5432. [PubMed](#)
10. Ma M, *et al.* 2021. *Front Immunol.* 12:718779. [PubMed](#)
11. LaMothe RA *et al.* 2018. *Front Immunol.* 0.570138889 . [PubMed](#)
12. Perner C, *et al.* 2020. *Immunity.* 53(5):1063-1077.e7. [PubMed](#)
13. Polavaram NS, *et al.* 2021. *Bone Res.* 9:24. [PubMed](#)
14. Colliou N, *et al.* 2018. *Gut Microbes.* 9:279. [PubMed](#)
15. Chen D, *et al.* 2020. *Cancer Immunol Res.* 8:883. [PubMed](#)
16. Horn KJ, *et al.* 2022. *Nat Commun.* 13:3321. [PubMed](#)
17. Rohner L, *et al.* 2020. *Sci Rep.* 0.570833333. [PubMed](#)
18. Liu Z, *et al.* 2020. *Cell.* 178(6):1509-1525.e19. [PubMed](#)
19. Bitschar K, *et al.* 2019. *Nat Commun.* 10:2730. [PubMed](#)
20. Yang YS, *et al.* 2022. *Nat Commun.* 13:6175. [PubMed](#)
21. Tähtinen S, *et al.* 2022. *Nat Immunol.* 23:532. [PubMed](#)
22. Hoyer FF, *et al.* 2020. *Immunity.* 51(5):899-914.e7. [PubMed](#)
23. Dai X, *et al.* 2021. *Molecular Cell.* 81(11):2317-2331.e6. [PubMed](#)
24. Barsoumian HB, *et al.* 2022. *Cancers (Basel).* 14:. [PubMed](#)
25. Chen KE, *et al.* 2021. *J Neuroinflammation.* 18:140. [PubMed](#)
26. Pattabiraman G, *et al.* 2021. *Am J Physiol Renal Physiol.* . [PubMed](#)
27. Liu Z, *et al.* 2020. *STAR Protoc.* 1:100029. [PubMed](#)
28. Kong IY, *et al.* 2020. *Cell Rep.* 33:108290. [PubMed](#)
29. Tran S, *et al.* 2020. *Immunity.* 53(3):627-640.e5. [PubMed](#)
30. Jiang W, *et al.* 2021. *Cell Rep.* 37:110112. [PubMed](#)
31. Yu YR, *et al.* 2020. *STAR Protocols.* 11(3):e0150606. [PubMed](#)
32. Yáñez A *et al.* 2017. *Immunity.* 47(5):890-902 . [PubMed](#)
33. Spanier J, *et al.* 2016. *Nat Commun.* 7:11804. [PubMed](#)

RRID

AB_2562010 (BioLegend Cat. No. 117337)
AB_2686978 (BioLegend Cat. No. 117353)
AB_2562016 (BioLegend Cat. No. 117338)

Antigen Details


Structure	Integrin α -chain, associates with integrin β_2 (CD18), 150 kD
Distribution	Dendritic cells, NK cells, intestinal intraepithelial lymphocytes (IEL), some activated T cells
Function	Cellular adhesion

Ligand/Receptor	iC3b, fibrinogen
Cell Type	Dendritic cells, Epithelial cells, NK cells, T cells, Tregs
Biology Area	Cell Adhesion, Cell Biology, Costimulatory Molecules, Immunology, Innate Immunity, Neuroscience, Neuroscience Cell Markers
Molecular Family	Adhesion Molecules, CD Molecules
Antigen References	1. Barclay A, <i>et al.</i> 1997. The Leukocyte Antigen Facts Book Academic Press. 2. Springer TA. 1994. <i>Cell</i> 76:301. 3. Lopez-Rodriguez C, <i>et al.</i> 1996. <i>J. Immunol.</i> 156:3780.
Gene ID	16411

Other Formats

APC anti-mouse CD11c, Biotin anti-mouse CD11c, FITC anti-mouse CD11c, PE anti-mouse CD11c, Purified anti-mouse CD11c, Alexa Fluor® 488 anti-mouse CD11c, Alexa Fluor® 647 anti-mouse CD11c, PE/Cyanine5 anti-mouse CD11c, PE/Cyanine7 anti-mouse CD11c, Brilliant Violet 605™ anti-mouse CD11c, Alexa Fluor® 700 anti-mouse CD11c, Pacific Blue™ anti-mouse CD11c, APC/Cyanine7 anti-mouse CD11c, PerCP/Cyanine5.5 anti-mouse CD11c, PerCP anti-mouse CD11c, Brilliant Violet 421™ anti-mouse CD11c, Brilliant Violet 570™ anti-mouse CD11c, Brilliant Violet 785™ anti-mouse CD11c, Brilliant Violet 510™ anti-mouse CD11c, Brilliant Violet 650™ anti-mouse CD11c, Purified anti-mouse CD11c (Maxpar® Ready), Alexa Fluor® 594 anti-mouse CD11c, PE/Dazzle™ 594 anti-mouse CD11c, Brilliant Violet 711™ anti-mouse CD11c, APC/Fire™ 750 anti-mouse CD11c, TotalSeq™-A0106 anti-mouse CD11c, Brilliant Violet 750™ anti-mouse CD11c, TotalSeq™-B0106 anti-mouse CD11c, TotalSeq™-C0106 anti-mouse CD11c, KIRAVIA Blue 520™ anti-mouse CD11c, Spark Blue™ 550 anti-mouse CD11c, Spark NIR™ 685 anti-mouse CD11c, Spark UV™ 387 anti-mouse CD11c, Spark Red™ 718 anti-mouse CD11c

Product Data

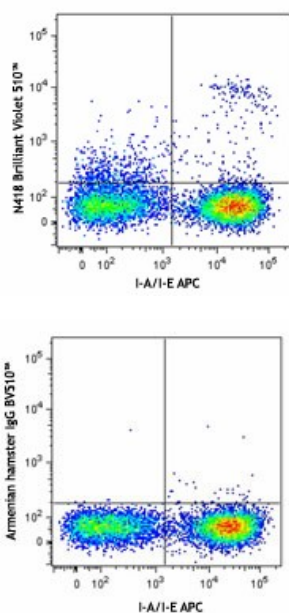

C57BL/6 mouse splenocytes were stained with mouse I-A/I-E APC and CD11c (clone N418) Brilliant Violet 510™ (top) or Armenian hamster IgG Brilliant Violet 510™ isotype control (bottom).

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587