

Brilliant Violet 605™ anti-human CD117 (c-kit) Antibody

Catalog# / Size	313217 / 25 tests 313218 / 100 tests
Clone	104D2
Regulatory Status	RUO
Other Names	Stem cell factor receptor, c-kit, mast cell growth factor receptor, steel factor receptor
Isotype	Mouse IgG1, κ
Description	CD117 is a 145 kD protein tyrosine kinase also known as c-Kit. It is a receptor for stem cell factor or c-Kit ligand. CD117 is expressed on pluripotent hematopoietic progenitor cells (approximately 1-4% bone marrow cells), mast cells, and acute myeloid leukemia cells (AML). CD117 binding of c-Kit ligand induces phosphorylation of CD117 and stimulates proliferation and survival of primitive hematopoietic stem cells as well as erythroid-committed and granulomonocytic committed cells.

Product Details

Verified Reactivity	Human
Reported Reactivity	Cynomolgus, Cow
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	MOLM-1 megakaryocytic cell line
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA).
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 605™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration, please enter the lot number in our Concentration and Expiration Lookup or Certificate of Analysis online tools.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	<p>Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 µl per million cells in 100 µl staining volume or 5 µl per 100 µl of whole blood.</p> <p>Brilliant Violet 605™ excites at 405 nm and emits at 603 nm. The bandpass filter 610/20 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 605™ is a trademark of Sirigen Group Ltd.</p> <p>Learn more about Brilliant Violet™.</p> <p>This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.</p>
Excitation Laser	Violet Laser (405 nm)
Application Notes	The 104D2 antibody does not block binding of c-Kit ligand. Additional reported applications (for the relevant formats) include: immunoprecipitation ¹ , immunofluorescence microscopy ¹ , and spatial

biology (IBEX)^{4,5}.

Application References

(PubMed link indicates BioLegend citation)

1. Broudy VC, *et al.* 1999. *Blood* 94:1979. (IF, IP)
2. Yoshino N, *et al.* 2000. *Exp. Anim. (Tokyo)* 49:97. (FC)
3. Nagano M, *et al.* 2007. *Blood* 110:151. (FC) [PubMed](#)
4. Radtke AJ, *et al.* 2020. *Proc Natl Acad Sci U S A.* 117:33455-65. (SB) [PubMed](#)
5. Radtke AJ, *et al.* 2022. *Nat Protoc.* 17:378-401. (SB) [PubMed](#)

Product Citations

1. Mashiko S, *et al.* 2015. *J Allergy Clin Immunol.* 136: 351-359. [PubMed](#)
2. Vijayakumar B, *et al.* 2022. *Immunity.* . [PubMed](#)
3. Bennion BG, *et al.* 2020. *Cell Reports.* 31(11):107771. [PubMed](#)
4. Caduff N, *et al.* 2021. *Cell Reports.* 35(5):109056. [PubMed](#)
5. Izadi D, *et al.* 2019. *Sci Adv.* 5:eaay0370. [PubMed](#)
6. Yang R, *et al.* 2020. *Cell.* 183(7):1826-1847.e31. [PubMed](#)
7. Alhaj Hussen K, *et al.* 2017. *Immunity.* 47:680. [PubMed](#)
8. Surace L, *et al.* 2021. *Nat Immunol.* 22:1367. [PubMed](#)
9. Schuldt N, *et al.* 2015. *PLoS One.* 10: 0145762. [PubMed](#)
10. Pang WW, *et al.* 2019. *Blood.* 133:2069. [PubMed](#)
11. Riding AM, *et al.* 2022. *iScience.* 25:104660. [PubMed](#)
12. Hernandez DC, *et al.* 2021. *Immunity.* .: [PubMed](#)

RRID

AB_2562024 (BioLegend Cat. No. 313217)
AB_2562025 (BioLegend Cat. No. 313218)

Antigen Details

Structure	Growth factor receptor with tyrosine kinase activity, subclass III, approximately 145 kD
Distribution	Pluripotent hematopoietic progenitor cells (approximately 1-4% bone marrow cells), mast cells, acute myeloid leukemic cells (AML)
Function	Growth factor receptor for stem cell factor. Induces proliferation and survival of primitive hematopoietic progenitors. Potent inducer of proliferation in erythroid-committed progenitor cells. Defects in CD117 have been linked to severe anemia and a decreased number of hematopoietic progenitor cells.
Ligand/Receptor	c-Kit ligand
Modification	Multiple phosphorylation sites
Cell Type	Embryonic Stem Cells, Hematopoietic stem and progenitors, Leukemia, Mast cells, Mesenchymal Stem Cells
Biology Area	Immunology, Stem Cells
Molecular Family	CD Molecules
Antigen References	1. Giebel LB, <i>et al.</i> 1992. <i>Oncogene</i> 7:2207. 2. Furitsu T, <i>et al.</i> 1993. <i>J. Clin. Invest.</i> 92:1736.
Gene ID	3815

Related Protocols


[Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

Purified anti-human CD117 (c-kit), PE anti-human CD117 (c-kit), APC anti-human CD117 (c-kit), Biotin anti-human CD117 (c-kit), PE/Cyanine5 anti-human CD117 (c-kit), PE/Cyanine7 anti-human CD117 (c-kit), PerCP/Cyanine5.5 anti-human CD117 (c-kit), Brilliant Violet 421™ anti-human CD117 (c-kit), Brilliant Violet 605™ anti-human CD117 (c-kit), Brilliant Violet 510™ anti-human CD117 (c-kit), Brilliant Violet 650™ anti-human CD117 (c-kit), Purified anti-human CD117 (c-kit) (Maxpar® Ready), PE/Dazzle™ 594 anti-human CD117 (c-kit), APC/Cyanine7 anti-human CD117 (c-kit), Brilliant Violet 711™ anti-human CD117 (c-kit), FITC anti-human CD117 (c-kit), Alexa Fluor® 488 anti-human CD117 (c-kit), Alexa Fluor® 647 anti-human CD117 (c-kit), APC/Fire™ 750 anti-human CD117 (c-kit), Brilliant Violet 785™ anti-human CD117 (c-kit), TotalSeq™-A0061 anti-human CD117 (c-kit), TotalSeq™-C0061 anti-human CD117 (c-kit), TotalSeq™-B0061 anti-human CD117 (c-kit), Alexa Fluor® 700 anti-human CD117 (c-kit), Spark

NIR™ 685 anti-human CD117 (c-kit) Antibody, APC/Fire™ 750 anti-human CD117 (c-kit), TotalSeq™-D0061 anti-human CD117 (c-kit), GMP APC anti-human CD117 (c-kit), GMP PE anti-human CD117 (c-kit)

Product Data

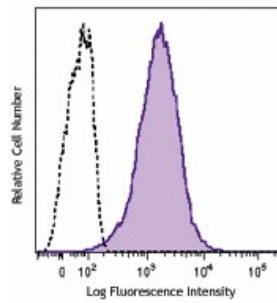

Human erythroleukemia cell line (HEL) was stained with CD117 (clone 104D2) Brilliant Violet 605™ (filled histogram) or mouse IgG1, κ Brilliant Violet 605™ isotype control (open histogram).

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587