

## CD4 Fluorophore Sampler Kit with Veri-Cells™ PBMC

<b>Catalog# / Size</b>	300562 / 1 kit
<b>Clone</b>	RPA-T4
<b>Regulatory Status</b>	RUO
<b>Other Names</b>	T4, staining index, spillover spreading
<b>Isotype</b>	Mouse IgG1, κ
<b>Description</b>	The CD4 Fluorophore Sampler Kit contains 18 different fluorophore formats of clone RPA-T4 with 10 tests of each antibody per vial. It can be useful to have the same clone in multiple fluorophore formats to enable instrument characterization and standardization. This kit also contains a vial of 25 tests of Veri-Cells™ PBMC.

CD4, also known as T4, is a 55 kD single-chain type I transmembrane glycoprotein expressed on most thymocytes, a subset of T cells, and monocytes/macrophages. CD4, a member of the Ig superfamily, recognizes antigens associated with MHC class II molecules, and participates in cell-cell interactions, thymic differentiation, and signal transduction. CD4 acts as a primary receptor for HIV, binding to HIV gp120. CD4 has also been shown to interact with IL-16.

### Product Details

---

<b>Verified Reactivity</b>	Human, Chimpanzee
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Mouse
<b>Formulation</b>	Brilliant Violet™ conjugates: Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA). All other conjugates: Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and 0.2% (w/v) BSA (origin USA).
<b>Storage &amp; Handling</b>	The kit should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a> . For flow cytometric staining, the suggested use of this reagent is 5 µl per million cells or 5 µl per 100 µl of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.  Alexa Fluor® and Pacific Blue™ are trademarks of Life Technologies Corporation.  <a href="#">View full statement regarding label licenses</a>  Brilliant Violet 421™, Brilliant Violet 510™, Brilliant Violet 570™, Brilliant Violet 605™, Brilliant Violet 650™, Brilliant Violet 711™, Brilliant Violet 785™ are trademarks of Sirigen Group Ltd.  <a href="#">Learn more about Brilliant Violet™.</a>  This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.
<b>Excitation Laser</b>	Violet Laser (405 nm) Blue Laser (488 nm) Green Laser (532 nm)/Yellow-Green Laser (561 nm) Red Laser (633 nm)
<b>Application Notes</b>	The RPA-T4 antibody binds to the D1 domain of CD4 (CDR1 and CDR3 epitopes) and can block

HIV gp120 binding and inhibit syncytia formation. Additional reported applications (for the relevant formats) include: immunohistochemistry of acetone-fixed frozen sections<sup>3,4,5</sup>, and blocking of T cell activation<sup>1,2</sup>. This clone was tested in-house and does not work on formalin fixed paraffin-embedded (FFPE) tissue. The LEAF™ purified antibody (Endotoxin <0.1 EU/μg, Azide-Free, 0.2 μm filtered) is recommended for functional assays (Cat. No. 300516).

## Application References

(PubMed link indicates BioLegend citation)

1. Knapp W, *et al.* 1989. *Leucocyte Typing IV*. Oxford University Press. New York. (Activ)
2. Moir S, *et al.* 1999. *J. Virol.* 73:7972. (Activ)
3. Deng MC, *et al.* 1995. *Circulation* 91:1647. (IHC)
4. Friedman T, *et al.* 1999. *J. Immunol.* 162:5256. (IHC)
5. Mack CL, *et al.* 2004. *Pediatr. Res.* 56:79. (IHC)
6. Lan RY, *et al.* 2006. *Hepatology* 43:729.
7. Zenaro E, *et al.* 2009. *J. Leukoc. Biol.* 86:1393. (FC) [PubMed](#)
8. Yoshino N, *et al.* 2000. *Exp. Anim. (Tokyo)* 49:97. (FC)

## Antigen Details

Cell Type	Dendritic cells, Tregs
Biology Area	Immunology
Molecular Family	CD Molecules
Gene ID	<a href="#">920</a>


## Related Protocols

[Protocol for using Veri-Cells™](#)


## Other Formats

APC anti-human CD4, Biotin anti-human CD4, FITC anti-human CD4, PE anti-human CD4, PE/Cyanine5 anti-human CD4, PE/Cyanine7 anti-human CD4, Purified anti-human CD4, APC/Cyanine7 anti-human CD4, Alexa Fluor® 488 anti-human CD4, Alexa Fluor® 647 anti-human CD4, Pacific Blue™ anti-human CD4, Brilliant Violet 421™ anti-human CD4, Alexa Fluor® 700 anti-human CD4, PerCP anti-human CD4, PerCP/Cyanine5.5 anti-human CD4, Brilliant Violet 570™ anti-human CD4, Brilliant Violet 650™ anti-human CD4, Purified anti-human CD4 (Maxpar® Ready), Alexa Fluor® 594 anti-human CD4, Brilliant Violet 510™ anti-human CD4, PE/Dazzle™ 594 anti-human CD4, Brilliant Violet 785™ anti-human CD4, Brilliant Violet 605™ anti-human CD4, Brilliant Violet 711™ anti-human CD4, APC/Fire™ 750 anti-human CD4, CD4 Fluorophore Sampler Kit, CD4 Fluorophore Sampler Kit with Veri-Cells™ PBMC, TotalSeq™-A0072 anti-human CD4, TotalSeq™-B0072 anti-human CD4, TotalSeq™-C0072 anti-human CD4, Ultra-LEAF™ Purified anti-human CD4, TotalSeq™-D0072 anti-human CD4

## Product Data

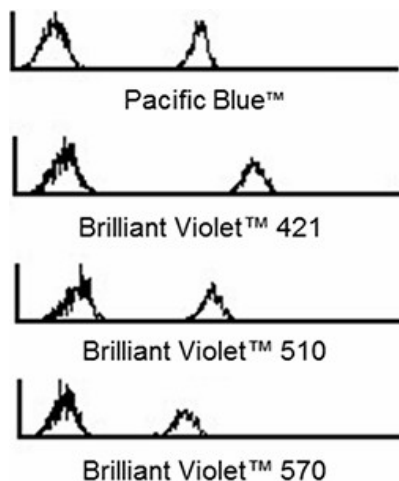

Veri-Cells™ PBMC were stained with anti-human CD4 (clone RPA-T4) in each fluorophore conjugate separately.

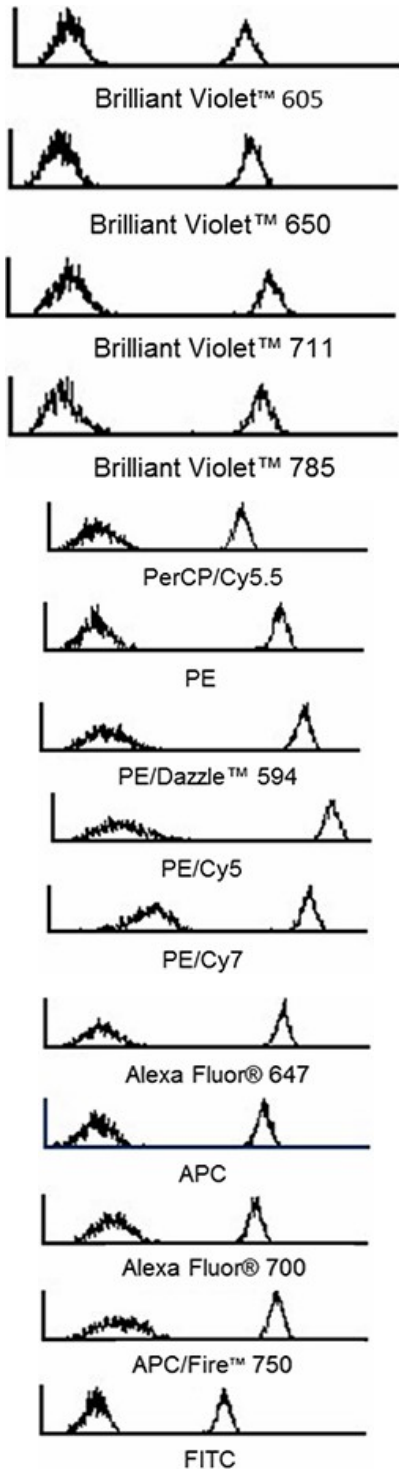

Veri-Cells™ PBMC were stained with anti-human CD4 (clone RPA-T4) in each fluorophore conjugate separately.

Veri-Cells™ PBMC were stained with anti-human CD4 (clone RPA-T4) in each fluorophore conjugate separately.

Veri-Cells™ PBMC were stained with anti-human CD4 (clone RPA-T4) in each fluorophore conjugate separately.

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

\*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, [www.biolegend.com/ordering#license](http://www.biolegend.com/ordering#license)). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
 Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587