

## PerCP anti-human CD4 Antibody

<b>Catalog# / Size</b>	300527 / 25 tests 300528 / 100 tests
<b>Clone</b>	RPA-T4
<b>Regulatory Status</b>	RUO
<b>Workshop</b>	IV T114
<b>Other Names</b>	T4
<b>Isotype</b>	Mouse IgG1, κ
<b>Description</b>	CD4, also known as T4, is a 55 kD single-chain type I transmembrane glycoprotein expressed on most thymocytes, a subset of T cells, and monocytes/macrophages. CD4, a member of the Ig superfamily, recognizes antigens associated with MHC class II molecules, and participates in cell-cell interactions, thymic differentiation, and signal transduction. CD4 acts as a primary receptor for HIV, binding to HIV gp120. CD4 has also been shown to interact with IL-16.

### Product Details

<b>Verified Reactivity</b>	Human
<b>Reported Reactivity</b>	Chimpanzee
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Mouse
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
<b>Preparation</b>	The antibody was purified by affinity chromatography, and conjugated with PerCP under optimal conditions.
<b>Concentration</b>	Lot-specific (to obtain lot-specific concentration, please enter the lot number in our <a href="#">Concentration and Expiration Lookup</a> or <a href="#">Certificate of Analysis</a> online tools.)
<b>Storage &amp; Handling</b>	The CD4 antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a> . For flow cytometric staining, the suggested use of this reagent is 5 µl per million cells in 100 µl staining volume or 5 µl per 100 µl of whole blood.  * PerCP has a maximum absorption of 482 nm and a maximum emission of 675 nm.
<b>Excitation Laser</b>	Blue Laser (488 nm)
<b>Application Notes</b>	The RPA-T4 antibody binds to the D1 domain of CD4 (CDR1 and CDR3 epitopes) and can block HIV gp120 binding and inhibit syncytia formation. Additional reported applications (for the relevant formats) include: immunohistochemistry of acetone-fixed frozen sections <sup>3,4,5</sup> , blocking of T cell activation <sup>1,2</sup> , and spatial biology (IBEX) <sup>10,11</sup> . This clone was tested in-house and does not work on formalin fixed paraffin-embedded (FFPE) tissue. The Ultra-LEAF™ purified antibody (Endotoxin < 0.01 EU/µg, Azide-Free, 0.2 µm filtered) is recommended for functional assays (Cat. No. 300569 - 300574).
<b>Application References</b>	<ol style="list-style-type: none"> <li>Knapp W, <i>et al.</i> 1989. Leucocyte Typing IV. Oxford University Press. New York. (Activ)</li> <li>Moir S, <i>et al.</i> 1999. <i>J. Virol.</i> 73:7972. (Activ)</li> <li>Deng MC, <i>et al.</i> 1995. <i>Circulation</i> 91:1647. (IHC)</li> <li>Friedman T, <i>et al.</i> 1999. <i>J. Immunol.</i> 162:5256. (IHC)</li> <li>Mack CL, <i>et al.</i> 2004. <i>Pediatr. Res.</i> 56:79. (IHC)</li> <li>Lan RY, <i>et al.</i> 2006. <i>Hepatology</i> 43:729.</li> <li>Zenaro E, <i>et al.</i> 2009. <i>J. Leukoc. Biol.</i> 86:1393. (FC) <a href="#">PubMed</a></li> </ol>
<b>(PubMed link indicates BioLegend citation)</b>	

8. Yoshino N, *et al.* 2000. *Exp. Anim. (Tokyo)* 49:97. (FC)
9. Stoeckius M, *et al.* 2017. *Nat. Methods.* 14:865. (PG)
10. Radtke AJ, *et al.* 2020. *Proc Natl Acad Sci USA.* 117:33455-33465. (SB) [PubMed](#)
11. Radtke AJ, *et al.* 2022. *Nat Protoc.* 17:378-401. (SB) [PubMed](#)

## Product Citations

1. Cao B, *et al.* 2022. *Nat Commun.* 13:6203. [PubMed](#)
2. De Hert E, *et al.* 2022. *Int J Mol Sci.* 23:. [PubMed](#)
3. Tran TM, *et al.* 2020. *Immunity.* 51(4):750-765. [PubMed](#)
4. Riner D, *et al.* 2016. *PLoS Negl Trop Dis.* 10:e0005180. [PubMed](#)
5. Petit N, *et al.* 2015. *PLoS One.* 10: 0138420. [PubMed](#)
6. Diao B, *et al.* 2020. *Front Immunol.* 1.032638889. [PubMed](#)
7. Zhao J, *et al.* 2021. *Front Immunol.* 12:658420. [PubMed](#)
8. Rajashekar JK, *et al.* 2021. *Cell Host Microbe.* 29(6):904-916.e6. [PubMed](#)
9. Boyle M, *et al.* 2015. *J Infect Dis.* 212: 416-425. [PubMed](#)
10. Bikker A, *et al.* 2012. *Ann Rheum Dis.* 71:1027. [PubMed](#)
11. Wouters K, *et al.* 2017. *Sci Rep.* 7:42665. [PubMed](#)
12. Vallejo-Gracia A, *et al.* 2020. *Nat Microbiol.* 1.002777778. [PubMed](#)
13. Cherneha M, *et al.* 2021. *Viruses.* 13:. [PubMed](#)
14. Smolders J, *et al.* 2010. *PLoS One.* 5:e15235. [PubMed](#)
15. Woodberry T, *et al.* 2017. *Infect Immun.* 85:e00986. [PubMed](#)
16. Ye H, *et al.* 2016. *Exp Mol Pathol.* 101:341-345. [PubMed](#)
17. Nakaya H, *et al.* 2016. *Proc Natl Acad Sci U S A.* 113: 1853 - 1858. [PubMed](#)
18. Verboven K, *et al.* 2018. *Sci Rep.* 8:4677. [PubMed](#)
19. Le J, *et al.* 2020. *Immunity.* 52(6):1105-1118.e9. [PubMed](#)
20. Hirschberger S, *et al.* 2021. *EMBO Mol Med.* 13:e14323. [PubMed](#)
21. FitzPatrick MEB, *et al.* 2021. *Cell Rep.* 34:108661. [PubMed](#)
22. Boehm D *et al.* 2017. *Cell host & microbe.* 21(5):569-579. [PubMed](#)
23. Conti L, *et al.* 2014. *J Immunol.* 188:1011. [PubMed](#)
24. Venugopal G, *et al.* 2017. *PLoS One.* . 10.1371/journal.pone.0188138. [PubMed](#)

## RRID

AB\_893327 (BioLegend Cat. No. 300527)  
 AB\_893321 (BioLegend Cat. No. 300528)

## Antigen Details

---

<b>Structure</b>	Ig superfamily, type I transmembrane glycoprotein, 55 kD
<b>Distribution</b>	T cell subset, majority of thymocytes, monocytes/macrophages
<b>Function</b>	MHC class II co-receptor, lymphocyte adhesion, thymic differentiation, HIV receptor
<b>Ligand/Receptor</b>	MHC class II molecules, HIV gp120, IL-16
<b>Cell Type</b>	Dendritic cells, Macrophages, Monocytes, T cells, Thymocytes, Tregs
<b>Biology Area</b>	Immunology
<b>Molecular Family</b>	CD Molecules
<b>Antigen References</b>	<ol style="list-style-type: none"> <li>1. Center D, <i>et al.</i> 1996. <i>Immunol. Today</i> 17:476.</li> <li>2. Gaubin M, <i>et al.</i> 1996. <i>Eur. J. Clin. Chem. Clin. Biochem.</i> 34:723.</li> </ol>
<b>Gene ID</b>	<a href="#">920</a>

## Related Protocols

---

[Cell Surface Flow Cytometry Staining Protocol](#)

## Other Formats


---

APC anti-human CD4, Biotin anti-human CD4, FITC anti-human CD4, PE anti-human CD4, PE/Cyanine5 anti-human CD4, PE/Cyanine7 anti-human CD4, Purified anti-human CD4, APC/Cyanine7 anti-human CD4, Alexa Fluor® 488 anti-human CD4, Alexa Fluor® 647 anti-human CD4, Pacific Blue™ anti-human CD4, Brilliant Violet 421™ anti-human CD4, Alexa Fluor® 700 anti-human CD4, PerCP anti-human CD4, PerCP/Cyanine5.5 anti-human CD4, Brilliant Violet 570™ anti-human CD4, Brilliant Violet 650™ anti-human CD4, Purified anti-human CD4 (Maxpar® Ready), Alexa Fluor® 594 anti-human CD4, Brilliant Violet 510™ anti-human CD4, PE/Dazzle™ 594 anti-human CD4, Brilliant Violet 785™ anti-human CD4, Brilliant Violet 605™ anti-human CD4, Brilliant Violet 711™ anti-human CD4, APC/Fire™ 750 anti-human CD4, CD4 Fluorophore Sampler Kit, CD4 Fluorophore Sampler Kit with Veri-Cells™

PBMC, TotalSeq™-A0072 anti-human CD4, TotalSeq™-B0072 anti-human CD4, TotalSeq™-C0072 anti-human CD4, Ultra-LEAF™ Purified anti-human CD4, TotalSeq™-D0072 anti-human CD4

## Product Data

---

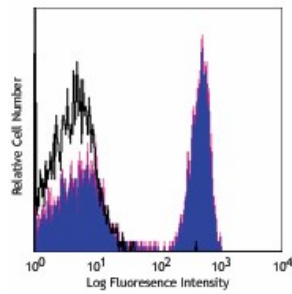

Human peripheral blood lymphocytes  
stained with RPA-T4 PerCP

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

\*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, [www.biolegend.com/ordering#license](http://www.biolegend.com/ordering#license)). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587