

Brilliant Violet 510™ anti-mouse CD45.2 Antibody

Catalog# / Size	109837 / 125 µL 109838 / 50 µg
Clone	104
Regulatory Status	RUO
Other Names	Ly-5.2, LCA
Isotype	Mouse (SJL) IgG2a, κ
Description	CD45.2 is an alloantigen of CD45, expressed by Ly5.2 bearing mouse strains (e.g., A, AKR, BALB/c, CBA/Ca, CBA/J, C3H/He, C57BL, C57BR, C57L, C58, DBA/1, DBA/2, NZB, SWR, 129). CD45, a member of the protein tyrosine phosphatase (PTP) family, is a 180-240 kD glycoprotein expressed on all hematopoietic cells except mature erythrocytes and platelets. There are multiple isoforms in the mouse that play key roles in TCR and BCR signal transduction. These isoforms are very specific to the activation and maturation states of the cell as well as specific cell type. The primary ligands for CD45 are galectin-1, CD2, CD3, CD4, TCR, CD22, and Thy-1.

Product Details

Verified Reactivity	Mouse
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	B10.S mouse thymocytes and splenocytes
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA).
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 510™ under optimal conditions.
Concentration	µg size: 0.2 mg/mL µL size: lot-specific (to obtain lot-specific concentration, please enter the lot number in our Concentration and Expiration Lookup or Certificate of Analysis online tools.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	<p>Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining using the µL size, the suggested use of this reagent is 5 µl per million cells in 100 µl staining volume or 5 µl per 100 µl of whole blood. For flow cytometric staining using the µg size, the suggested use of this reagent is ≤0.5 µg per million cells in 100 µl volume. It is recommended that the reagent be titrated for optimal performance for each application.</p> <p>Brilliant Violet 510™ excites at 405 nm and emits at 510 nm. The bandpass filter 510/50 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 510™ is a trademark of Sirigen Group Ltd.</p> <p>Learn more about Brilliant Violet™.</p> <p>This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.</p>
Excitation Laser	Violet Laser (405 nm)

Application Notes

The 104 antibody does not react with mouse cells expressing the CD45.1 alloantigen. Additional reported applications (for the relevant formats) include: immunoprecipitation⁴, *in vivo* and *in vitro* blocking of B cell responses^{1,2}, and immunohistochemical staining of acetone-fixed frozen sections³.

Application References

1. Yakura H, *et al.* 1983. *J. Exp. Med.* 157:1077. (Block)
2. Yakura H, *et al.* 1986. *J. Immunol.* 136:2729. (Block)
3. Suzuki K, *et al.* 2000. *Immunity* 13:691. (IHC)
4. Shen FW, *et al.* 1986. *Immunogenetics* 24:146. (IP)
5. Baldwin TA and Hogquist KA. 2007. *J. Immunol.* 179:837.
6. Pascal V, *et al.* 2007. *J. Immunol.* 179:1751.
7. Burman AC, *et al.* 2007. *Blood* 110:1064.
8. Kincaid EZ, *et al.* 2007. *J. Immunol.* 179:3187.
9. Phan TG, *et al.* 2007. *Nature Immunol.* 8:992.
10. Nakano-Yokomizo T, *et al.* 2011. *J. Exp Med.* 208:1661. [PubMed](#)
11. Wen T, *et al.* 2013. *PNAS.* 110:6067. [PubMed](#)
12. Kohlmeier JE, *et al.* 2008. *Immunity.* 29:101. (FC) [PubMed](#)

Product Citations

1. Wang L, *et al.* 2019. *Cell Rep.* 29:1848. [PubMed](#)
2. Han SJ *et al.* 2017. *Immunity.* 47(6):1154-1168. [PubMed](#)
3. Hirata Y *et al.* 2018. *Cell stem cell.* 22(3):445-453. [PubMed](#)
4. Delás MJ, *et al.* 2019. *Cell Rep.* 27:719. [PubMed](#)
5. Ng KK *et al.* 2018. *eLife.* 7 pii: e37851. [PubMed](#)
6. Zhang Z, *et al.* 2021. *Oncoimmunology.* 10:1912472. [PubMed](#)
7. Lam WY *et al.* 2018. *Cell reports.* 24(9):2479-2492. [PubMed](#)
8. Liu D *et al.* 2019. *Immunity.* 51(1):64-76. [PubMed](#)
9. Yang L, *et al.* 2021. *Cell Death Differ.* 28:2616. [PubMed](#)
10. Canton J, *et al.* 2021. *Nat Immunol.* 22:140. [PubMed](#)
11. Glaros V, *et al.* 2021. *Immunity.* 54:2005. [PubMed](#)
12. Milner JJ, *et al.* 2020. *Immunity.* 52(5):808-824.e7. [PubMed](#)
13. Gniadek TJ, *et al.* 2020. *J Immunother.* 43:217. [PubMed](#)
14. Jonas B, *et al.* 2016. *PLoS One.* 11: 0159189. [PubMed](#)
15. Izumi G, *et al.* 2021. *Nat Commun.* 12:5029. [PubMed](#)
16. van der Velden V, *et al.* 2017. *J Immunol Methods.* 10.1016/j.jim.2017.03.011. [PubMed](#)
17. Dutton EE, *et al.* 2018. *Wellcome Open Res.* 2:117. [PubMed](#)
18. Jing Y, *et al.* 2019. *J Allergy Clin Immunol.* 144:1377. [PubMed](#)
19. Yun JH, *et al.* 2021. *JCI Insight.* 6:. [PubMed](#)

RRID

AB_2561393 (BioLegend Cat. No. 109837)
AB_2650900 (BioLegend Cat. No. 109838)

Antigen Details

Structure	Protein tyrosine phosphatase (PTP) family, 180-240 kD
Distribution	All hematopoietic cells except mature erythrocytes and platelets of the CD45.2 strain of mice
Function	Phosphatase, T and B cell activation
Ligand/Receptor	Galectin-1, CD2, CD3, CD4
Biology Area	Cell Biology, Immunology, Inhibitory Molecules, Innate Immunity, Neuroscience, Neuroscience Cell Markers
Molecular Family	CD Molecules
Antigen References	1. Suzuki K, <i>et al.</i> 2000. <i>Immunity</i> 13:691.
Gene ID	19264


Related Protocols

[Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

Biotin anti-mouse CD45.2, FITC anti-mouse CD45.2, PE anti-mouse CD45.2, Purified anti-mouse CD45.2, APC anti-mouse CD45.2, Alexa Fluor® 488 anti-mouse CD45.2, Alexa Fluor® 647 anti-mouse CD45.2, Pacific Blue™ anti-mouse CD45.2, Alexa Fluor® 700 anti-mouse CD45.2, APC/Cyanine7 anti-mouse CD45.2, PerCP anti-mouse CD45.2, PerCP/Cyanine5.5 anti-mouse CD45.2, PE/Cyanine7 anti-mouse CD45.2, Brilliant Violet 421™ anti-mouse CD45.2, Brilliant Violet 570™ anti-mouse CD45.2, Brilliant Violet 650™ anti-mouse CD45.2, Brilliant Violet 510™ anti-mouse CD45.2, Brilliant Violet 785™ anti-mouse CD45.2, Brilliant Violet 605™ anti-mouse CD45.2, Purified anti-mouse CD45.2 (Maxpar® Ready), PE/Dazzle™ 594 anti-mouse CD45.2, Brilliant Violet 711™ anti-mouse CD45.2, Alexa Fluor® 594 anti-mouse CD45.2, APC/Fire™ 750 anti-mouse CD45.2, TotalSeq™-A0157 anti-mouse CD45.2, TotalSeq™-B0157 anti-mouse CD45.2, TotalSeq™-C0157 anti-mouse CD45.2, Brilliant Violet 750™ anti-mouse CD45.2, Spark Blue™ 550 anti-mouse CD45.2, Spark NIR™ 685 anti-mouse CD45.2

Product Data

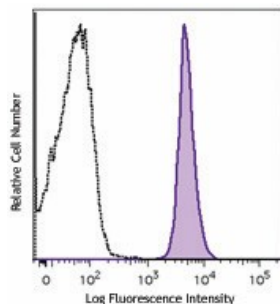

C57BL/6 mouse splenocytes were stained with CD45.2 (clone 104) Brilliant Violet 510™ (filled histogram) or mouse IgG2a, κ Brilliant Violet 510™ isotype control (open histogram).

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-BioLegend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587