

## PerCP/Cyanine5.5 anti-human CD11c Antibody

<b>Catalog# / Size</b>	371517 / 25 tests 371518 / 100 tests
<b>Clone</b>	S-HCL-3
<b>Regulatory Status</b>	RUO
<b>Other Names</b>	Integrin $\alpha$ X subunit, CR4, p150, ITGAX
<b>Isotype</b>	Mouse IgG2b, $\kappa$
<b>Description</b>	CD11c is a 145-150 kD type I transmembrane glycoprotein also known as integrin $\alpha$ X and CR4. CD11c non-covalently associates with integrin $\beta$ 2 (CD18) and is expressed on monocytes/macrophages, dendritic cells, granulocytes, NK cells, and subsets of T and B cells. CD11c has been reported to play a role in adhesion and CTL killing through its interactions with fibrinogen, CD54, and iC3b.

### Product Details

---

<b>Verified Reactivity</b>	Human
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Mouse
<b>Immunogen</b>	Spleen cells from patient diagnosed with hairy cell leukemia.
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
<b>Preparation</b>	The antibody was purified by affinity chromatography and conjugated with PerCP/Cyanine5.5 under optimal conditions.
<b>Concentration</b>	Lot-specific (to obtain lot-specific concentration, please enter the lot number in our <a href="#">Concentration and Expiration Lookup</a> or <a href="#">Certificate of Analysis</a> online tools.)
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a> . For flow cytometric staining, the suggested use of this reagent is 5 $\mu$ l per million cells in 100 $\mu$ l staining volume or 5 $\mu$ l per 100 $\mu$ l of whole blood.  * PerCP/Cyanine5.5 has a maximum absorption of 482 nm and a maximum emission of 690 nm.
<b>Excitation Laser</b>	Blue Laser (488 nm)
<b>Application Notes</b>	Additional reported applications (for the relevant formats) include: immunohistochemistry on frozen tissue sections <sup>1,2,3,4</sup> and immunoprecipitation <sup>1</sup> .
<b>Additional Product Notes</b>	BioLegend is in the process of converting the name PerCP/Cy5.5 to PerCP/Cyanine5.5. The dye molecule remains the same, so you should expect the same quality and performance from our PerCP/Cyanine5.5 products. Contact <a href="#">Technical Service</a> if you have any questions.
<b>Application References</b>	<ol style="list-style-type: none"><li>1. Schwarting R, <i>et al.</i> 1985. <i>Blood</i> 65:974.</li><li>2. Knowles DM, <i>et al.</i> 1990. <i>Am. J. Pathol.</i> 136:29.</li><li>3. Vandenabeele S, <i>et al.</i> 2001. <i>Blood</i> 97:1733.</li><li>4. Shaw JL, <i>et al.</i> 2011. <i>J. Reprod. Immunol.</i> 89:84.</li></ol>
<b>RRID</b>	AB_2650800 (BioLegend Cat. No. 371517) AB_2650801 (BioLegend Cat. No. 371518)

## Antigen Details

---

<b>Structure</b>	Integrin, type I transmembrane glycoprotein, associates with integrin $\beta 2$ (CD18), 145-150 kD.
<b>Distribution</b>	Myeloid, dendritic cells, NK cells, B cells and T cell subsets.
<b>Function</b>	Adhesion and CTL killing.
<b>Interaction</b>	Interacts with fibrinogen, CD54, and iC3b.
<b>Ligand/Receptor</b>	CD54, fibrinogen, iC3b, ICAM-1, ICAM-4
<b>Cell Type</b>	B cells, Dendritic cells, NK cells, T cells
<b>Biology Area</b>	Costimulatory Molecules, Immunology
<b>Molecular Family</b>	CD Molecules
<b>Antigen References</b>	1. Petty HR, Todd RF 3rd. 1996. <i>Immunol. Today</i> 17:209. 2. Springer T. 1994. <i>Cell</i> 76:301. 3. Ihanus E, <i>et al.</i> 2007. <i>Blood</i> 109:802-10.
<b>Gene ID</b>	<a href="#">3687</a>

## Related Protocols

---

[Cell Surface Flow Cytometry Staining Protocol](#)


## Other Formats

---

Purified anti-human CD11c, APC anti-human CD11c, PE anti-human CD11c, PE/Cyanine7 anti-human CD11c, Brilliant Violet 421™ anti-human CD11c, APC/Fire™ 750 anti-human CD11c, FITC anti-human CD11c, PerCP/Cyanine5.5 anti-human CD11c, Brilliant Violet 510™ anti-human CD11c, TotalSeq™-A0053 anti-human CD11c, TotalSeq™-C0053 anti-human CD11c, TotalSeq™-B0053 anti-human CD11c, Alexa Fluor® 647 anti-human CD11c Antibody, TotalSeq™-D0053 anti-human CD11c

## Product Data

---

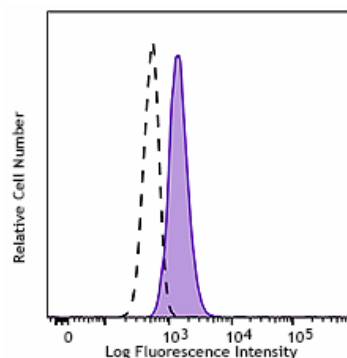

Human peripheral blood granulocytes were stained with CD11c (clone S-HCL-3) PerCP/Cyanine5.5 (filled histogram) or mouse IgG2b, ? PerCP/Cyanine5.5 isotype control (open histogram)

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

\*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, [www.biolegend.com/ordering#license](http://www.biolegend.com/ordering#license)). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587

