

## Brilliant Violet 510™ anti-human CD163 Antibody

<b>Catalog# / Size</b>	333627 / 25 tests 333628 / 100 tests
<b>Clone</b>	GHI/61
<b>Regulatory Status</b>	RUO
<b>Workshop</b>	VI M38
<b>Other Names</b>	GHI/61, M130, RM3/1, p155, Hemoglobin/Haptoglobin Complex Receptor, macrophage-associated antigen, ED2(rat)
<b>Isotype</b>	Mouse IgG1, κ
<b>Description</b>	CD163 is a member of the group B scavenger receptor cysteine-rich superfamily, also known as GHI/61, M130, RM3/1, p155, hemoglobin-haptoglobin complex receptor, or macrophage-associated antigen. It is a 134 kD (non-reduced)/155 kD (reduced) glycoprotein primarily expressed on macrophages, Kupffer cells, monocytes, a subset of dendritic cells, and a subset of hematopoietic stem/progenitor cells. CD163 binds to haptoglobin-hemoglobin complex and TWEAK, and plays a role in clearing hemoglobin and regulating cytokine production by macrophages. Membrane CD163 can be cleaved by metalloproteinases (MMP), resulting in a soluble form. Elevated serum level of sCD163 has been implicated in many kinds of inflammatory diseases.

### Product Details

---

<b>Verified Reactivity</b>	Human, Cynomolgus, Rhesus
<b>Reported Reactivity</b>	African Green
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Mouse
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA).
<b>Concentration</b>	Lot-specific (to obtain lot-specific concentration, please enter the lot number in our <a href="#">Concentration and Expiration Lookup</a> or <a href="#">Certificate of Analysis</a> online tools.)
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	<p>Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a>. For flow cytometric staining, the suggested use of this reagent is 5 µl per million cells in 100 µl staining volume or 5 µl per 100 µl of whole blood.</p> <p>Brilliant Violet 510™ excites at 405 nm and emits at 510 nm. The bandpass filter 510/50 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. <b>Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel.</b> Refer to your instrument manual or manufacturer for support. Brilliant Violet 510™ is a trademark of Sirigen Group Ltd.</p> <p><a href="#">Learn more about Brilliant Violet™.</a></p> <p>This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.</p>
<b>Excitation Laser</b>	Violet Laser (405 nm)
<b>Application Notes</b>	Clone GHI/61 binds to domain 7 of CD163. Additional reported applications (for the relevant

formats) include: immunocytochemical staining, immunoprecipitation, western blot<sup>1</sup>, and spatial biology (IBEX)<sup>6,7</sup>.

## Application References

(PubMed link indicates BioLegend citation)

1. Pulford K, *et al.* 1992. *Immunology* 75:588. (ICC, IP, WB)
2. Law SK, *et al.* 1993. *Eur. J. Immunol.* 23:2320.
3. Madsen M, *et al.* 2004. *J. Biol. Chem.* 279:51561.
4. Kim WK, *et al.* 2006. *Am. J. Pathol.* 168:822. (FC)
5. Buttari B, *et al.* 2011. *Atherosclerosis*. 215:316. [PubMed](#)
6. Radtke AJ, *et al.* 2020. *Proc Natl Acad Sci USA.* 117:33455-33465. (SB) [PubMed](#)
7. Radtke AJ, *et al.* 2022. *Nat Protoc.* 17:378-401. (SB) [PubMed](#)

## Product Citations

1. de Andrade LF, *et al.* 2019. *JCI Insight.* 4:e133103. [PubMed](#)
2. Hunegnaw R, *et al.* 2019. *Front Immunol.* 1.484027778. [PubMed](#)
3. Abd Hamid M *et al.* 2019. *Cancer Immunol Res.* 7(8):1293-1306. [PubMed](#)
4. Luoma AM, *et al.* 2020. *Cell.* 182(3):655-671.e22. [PubMed](#)

## RRID

AB\_2650631 (BioLegend Cat. No. 333627)  
AB\_2650632 (BioLegend Cat. No. 333628)

## Antigen Details

---

<b>Structure</b>	134 kD (non-reduced)/155 kD (reduced) glycoprotein, scavenger receptor superfamily
<b>Distribution</b>	Monocytes, macrophages, Kuffer cells, subset of dendritic cells, subset of hematopoietic stem/progenitor cells
<b>Function</b>	Clearance of haptoglobin-hemoglobin complex, regulation of cytokine production by macrophages
<b>Ligand/Receptor</b>	Haptoglobin-hemoglobin complex, TWEAK
<b>Cell Type</b>	Dendritic cells, Hematopoietic stem and progenitors, Macrophages, Monocytes
<b>Biology Area</b>	Cell Biology, Immunology, Innate Immunity, Neuroscience, Neuroscience Cell Markers
<b>Molecular Family</b>	CD Molecules
<b>Antigen References</b>	<ol style="list-style-type: none"><li>1. Roth J, <i>et al.</i> 1994 <i>Transplantation.</i> 57:127</li><li>2. Van den Heuvel MM, <i>et al.</i> 1999 <i>J. Leukoc. Biol.</i> 66:858</li><li>3. Sulahian TH, <i>et al.</i> 2000 <i>Cytokines</i> 12:1312</li><li>4. Fabrick BO, <i>et al.</i> 2007 <i>J. Neuroimmunol.</i> 187:179</li></ol>
<b>Gene ID</b>	<a href="#">9332</a>

## Related Protocols

---

[Cell Surface Flow Cytometry Staining Protocol](#)


## Other Formats

---

PerCP anti-human CD163, Purified anti-human CD163, Biotin anti-human CD163, PE anti-human CD163, PerCP/Cyanine5.5 anti-human CD163, APC anti-human CD163, Brilliant Violet 421™ anti-human CD163, PE/Cyanine7 anti-human CD163, Brilliant Violet 605™ anti-human CD163, FITC anti-human CD163, Alexa Fluor® 647 anti-human CD163, APC/Cyanine7 anti-human CD163, PE/Dazzle™ 594 anti-human CD163, Brilliant Violet 510™ anti-human CD163, Brilliant Violet 711™ anti-human CD163, Brilliant Violet 785™ anti-human CD163, APC/Fire™ 750 anti-human CD163, TotalSeq™-A0358 anti-human CD163, TotalSeq™-C0358 anti-human CD163, TotalSeq™-B0358 anti-human CD163, PE/Cyanine5 anti-human CD163, TotalSeq™-D0358 anti-human CD163

## Product Data

---

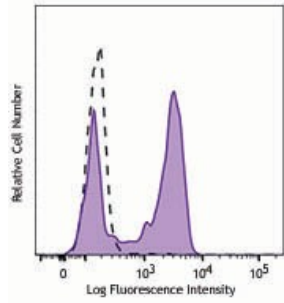

Human peripheral blood monocytes were stained with CD163 (clone GH1/61) Brilliant Violet 510™ (filled histogram) or mouse IgG1, κ Brilliant Violet 510™ isotype control (open histogram).

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

\*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, [www.biolegend.com/ordering#license](http://www.biolegend.com/ordering#license)). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587