

Brilliant Violet 605™ anti-human CD4 Antibody

Catalog# / Size	300555 / 25 tests 300556 / 100 tests
Clone	RPA-T4
Regulatory Status	RUO
Workshop	IV T114
Other Names	T4
Isotype	Mouse IgG1, κ
Description	CD4, also known as T4, is a 55 kD single-chain type I transmembrane glycoprotein expressed on most thymocytes, a subset of T cells, and monocytes/macrophages. CD4, a member of the Ig superfamily, recognizes antigens associated with MHC class II molecules, and participates in cell-cell interactions, thymic differentiation, and signal transduction. CD4 acts as a primary receptor for HIV, binding to HIV gp120. CD4 has also been shown to interact with IL-16.

Product Details

Verified Reactivity	Human
Reported Reactivity	Chimpanzee
Antibody Type	Monoclonal
Host Species	Mouse
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA).
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 605™ under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration, please enter the lot number in our Concentration and Expiration Lookup or Certificate of Analysis online tools.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	<p>Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 µl per million cells in 100 µl staining volume or 5 µl per 100 µl of whole blood.</p> <p>Brilliant Violet 605™ excites at 405 nm and emits at 603 nm. The bandpass filter 610/20 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 605™ is a trademark of Sirigen Group Ltd.</p> <p>Learn more about Brilliant Violet™.</p> <p>This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.</p>
Excitation Laser	Violet Laser (405 nm)
Application Notes	The RPA-T4 antibody binds to the D1 domain of CD4 (CDR1 and CDR3 epitopes) and can block HIV gp120 binding and inhibit syncytia formation. Additional reported applications (for the relevant formats) include: immunohistochemistry of acetone-fixed frozen sections ^{3,4,5} , blocking of T cell

activation^{1,2}, and spatial biology (IBEX)^{10,11}. This clone was tested in-house and does not work on formalin fixed paraffin-embedded (FFPE) tissue. The Ultra-LEAF™ purified antibody (Endotoxin < 0.01 EU/μg, Azide-Free, 0.2 μm filtered) is recommended for functional assays (Cat. No. 300569 - 300574).

Application References

(PubMed link indicates BioLegend citation)

1. Knapp W, *et al.* 1989. Leucocyte Typing IV. Oxford University Press. New York. (Activ)
2. Moir S, *et al.* 1999. *J. Virol.* 73:7972. (Activ)
3. Deng MC, *et al.* 1995. *Circulation* 91:1647. (IHC)
4. Friedman T, *et al.* 1999. *J. Immunol.* 162:5256. (IHC)
5. Mack CL, *et al.* 2004. *Pediatr. Res.* 56:79. (IHC)
6. Lan RY, *et al.* 2006. *Hepatology* 43:729.
7. Zenaro E, *et al.* 2009. *J. Leukoc. Biol.* 86:1393. (FC) [PubMed](#)
8. Yoshino N, *et al.* 2000. *Exp. Anim. (Tokyo)* 49:97. (FC)
9. Stoeckius M, *et al.* 2017. *Nat. Methods.* 14:865. (PG)
10. Radtke AJ, *et al.* 2020. *Proc Natl Acad Sci USA.* 117:33455-33465. (SB) [PubMed](#)
11. Radtke AJ, *et al.* 2022. *Nat Protoc.* 17:378-401. (SB) [PubMed](#)

Product Citations

1. Gerber AN, *et al.* 2021. *Cell Rep.* 37:109816. [PubMed](#)
2. Corleis B, *et al.* 2019. *Cell Rep.* 26:1409. [PubMed](#)
3. Braun J, *et al.* 2020. *Nature.* 587:270. [PubMed](#)
4. Álvaro de Mingo Pulido *et al.* 2018. *Cancer cell.* 33(1):60-74. [PubMed](#)
5. Kip P, *et al.* 2021. *Nutrients.* 13:. [PubMed](#)
6. Lindesmith LC, *et al.* 2020. *Cell Mol Gastroenterol Hepatol.* 0.586805556. [PubMed](#)
7. Juno JA, *et al.* 2020. *Nat Med.* 26:1428. [PubMed](#)
8. Bradley T *et al.* 2018. *Cell reports.* 25(1):107-117. [PubMed](#)
9. Csomós K, *et al.* 2022. *Nat Immunol.* 23:1256. [PubMed](#)
10. Olafsdottir TA, *et al.* 2022. *Commun Biol.* 5:914. [PubMed](#)

RRID

AB_2564390 (BioLegend Cat. No. 300555)
AB_2564391 (BioLegend Cat. No. 300556)

Antigen Details

Structure	Ig superfamily, type I transmembrane glycoprotein, 55 kD
Distribution	T cell subset, majority of thymocytes, monocytes/macrophages
Function	MHC class II co-receptor, lymphocyte adhesion, thymic differentiation, HIV receptor
Ligand/Receptor	MHC class II molecules, HIV gp120, IL-16
Cell Type	Dendritic cells, Macrophages, Monocytes, T cells, Thymocytes, Tregs
Biology Area	Immunology
Molecular Family	CD Molecules
Antigen References	<ol style="list-style-type: none">1. Center D, <i>et al.</i> 1996. <i>Immunol. Today</i> 17:476.2. Gaubin M, <i>et al.</i> 1996. <i>Eur. J. Clin. Chem. Clin. Biochem.</i> 34:723.
Gene ID	920


Related Protocols

[Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

APC anti-human CD4, Biotin anti-human CD4, FITC anti-human CD4, PE anti-human CD4, PE/Cyanine5 anti-human CD4, PE/Cyanine7 anti-human CD4, Purified anti-human CD4, APC/Cyanine7 anti-human CD4, Alexa Fluor® 488 anti-human CD4, Alexa Fluor® 647 anti-human CD4, Pacific Blue™ anti-human CD4, Brilliant Violet 421™ anti-human CD4, Alexa Fluor® 700 anti-human CD4, PerCP anti-human CD4, PerCP/Cyanine5.5 anti-human CD4, Brilliant Violet 570™ anti-human CD4, Brilliant Violet 650™ anti-human CD4, Purified anti-human CD4 (Maxpar® Ready), Alexa Fluor® 594 anti-human CD4, Brilliant Violet 510™ anti-human CD4, PE/Dazzle™ 594 anti-human CD4, Brilliant Violet 785™ anti-human CD4, Brilliant Violet 605™ anti-human CD4, Brilliant Violet 711™ anti-human CD4, APC/Fire™ 750 anti-human CD4, CD4 Fluorophore Sampler Kit, CD4 Fluorophore Sampler Kit with Veri-Cells™ PBMC, TotalSeq™-A0072 anti-human CD4, TotalSeq™-B0072 anti-human CD4, TotalSeq™-C0072 anti-human CD4, Ultra-LEAF™

Product Data

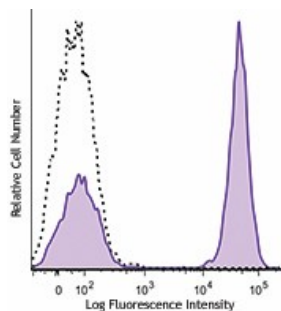

Human peripheral blood lymphocytes were stained with CD4 (clone RPA-T4) Brilliant Violet 605™ (filled histogram) or mouse IgG1, κ Brilliant Violet 605™ isotype control (open histogram).

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587