

Brilliant Violet 605™ anti-mouse CD69 Antibody

Catalog# / Size	104529 / 125 µL 104530 / 50 µg
Clone	H1.2F3
Regulatory Status	RUO
Other Names	Very Early Activation Antigen (VEA), AIM, EA1, MLR3, gp34/28
Isotype	Armenian Hamster IgG
Description	CD69 is a 60 kD type II membrane protein composed of a 27/33 kD disulfide-linked homodimer, also known as Very Early Activation Antigen (VEA), AIM, EA1, MLR3, and gp34/28. It is expressed on a subset of thymocytes and platelets. CD69 is rapidly induced on activated T and B cells, neutrophils, and NK cells. It is a C-type lectin, closely related to the NKR-P1 and Ly-49 NK cell activation molecules. CD69 is involved in the early events of cell activation and thymocyte positive selection.

Product Details

Verified Reactivity	Mouse
Antibody Type	Monoclonal
Host Species	Armenian Hamster
Immunogen	Mouse dendritic epidermal T cell line Y245
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA).
Preparation	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 605™ under optimal conditions.
Concentration	µg sizes: 0.2 mg/mL µL sizes: lot-specific (to obtain lot-specific concentration, please enter the lot number in our Concentration and Expiration Lookup or Certificate of Analysis online tools.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	<p>Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For immunofluorescent staining using the µg size, the suggested use of this reagent is ≤0.5 µg per million cells in 100 µl volume. For immunofluorescent staining using the µl size, the suggested use of this reagent is 5 µl per million cells in 100 µl staining volume or 5 µl per 100 µl of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.</p> <p>Brilliant Violet 605™ excites at 405 nm and emits at 603 nm. The bandpass filter 610/20 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 605™ is a trademark of Sirigen Group Ltd.</p> <p>Learn more about Brilliant Violet™.</p> <p>This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.</p>
Excitation Laser	Violet Laser (405 nm)

Application Notes The H1.2F3 antibody has been reported to augment T cell activation. Additional reported applications (for the relevant formats) include: *in vitro* T cell and NK cell activation¹⁻³, immunohistochemistry^{4,5}, and immunoprecipitation¹.

This antibody has been characterized in the literature as containing a lambda (?) light chain.

Application References

(PubMed link indicates BioLegend citation)

1. Yokoyama WM, *et al.* 1988. *J. Immunol.* 141:369. (IP)
2. Sobel ES, *et al.* 1993. *J. Immunol.* 150:673.
3. Karlhofer FM, *et al.* 1991. *J. Immunol.* 146:3662.
4. Zhou X, *et al.* 2005. *J. Biol. Chem.* 280:31240. (IHC)
5. Podd BS, *et al.* 2006. *J. Immunol.* 176:6532. (IHC)
6. Lawson BR, *et al.* 2007. *J. Immunol.* 178:5366.
7. Lee JW, *et al.* 2006. *Nature Immunol.* 8:181.
8. Epardaud M, *et al.* 2008. *Cancer Res.* 15:2972. [PubMed](#)
9. Jordan JM, *et al.* 2008. 76:3717. [PubMed](#)
10. Kenna TJ, *et al.* 2008. *Blood* 111:2091. [PubMed](#)
11. Ishikawa C, *et al.* 2013. *Biochim Biophys Acta.* 167:99. [PubMed](#)

Product Citations

1. Harsha Krovi S, *et al.* 2020. *Nat Commun.* 4.790277778. [PubMed](#)
2. Soon MSF, *et al.* 2020. *Nat Immunol.* 1.984027778. [PubMed](#)
3. Hodgins B, *et al.* 2019. *Immun Ageing.* 16:27. [PubMed](#)
4. Effer M, *et al.* 2022. *STAR Protoc.* 3:101038. [PubMed](#)
5. Aldon Y, *et al.* 2020. *J Immunol.* 204:903. [PubMed](#)
6. Toomer G, *et al.* 2022. *Viruses.* 14:. [PubMed](#)
7. Tan L, *et al.* 2019. *Cell Rep.* 27:3657. [PubMed](#)
8. Misumi I *et al.* 2019. *Cell Rep.* 27(2):514-524 . [PubMed](#)
9. Istaces N, *et al.* 2019. *Nat Commun.* 10:3306. [PubMed](#)
10. Boothby IC, *et al.* 2021. *Nature.* 599:667. [PubMed](#)
11. Hossain DMS, *et al.* 2018. *J Clin Invest.* 128:644. [PubMed](#)
12. Häselbarth L, *et al.* 2020. *EXCLI J.* 19:334. [PubMed](#)
13. Toubal A, *et al.* 2020. *Nat Commun.* 3755:11. [PubMed](#)
14. Kurniawan H, *et al.* 2020. *Cell Metabolism.* 31(5):920-936. [PubMed](#)
15. Schuldt N, *et al.* 2015. *PLoS One.* 10: 0145762. [PubMed](#)
16. Devi S, *et al.* 2021. *Immunity.* 54(6):1219-1230.e7. [PubMed](#)
17. Toshiro Hirai *et al.* 2019. *Immunity.* 50(5):1249-1261 . [PubMed](#)
18. Nelson CE *et al.* 2019. *Cell Rep.* 28(12):3092-3104 . [PubMed](#)
19. Highton AJ, *et al.* 2019. *Wellcome Open Res.* 4:78. [PubMed](#)
20. Hoover AR, *et al.* 2022. *Clin Transl Med.* 12:e937. [PubMed](#)
21. Quah B, *et al.* 2014. *J Vis Exp.* 88: 51627. [PubMed](#)
22. Laczkó D, *et al.* 2020. *Immunity.* 53:724. [PubMed](#)

RRID AB_11203710 (BioLegend Cat. No. 104529)
AB_2563062 (BioLegend Cat. No. 104530)

Antigen Details

Structure	C-type lectin, 27/33 kD
Distribution	Activated T cells and B cells, NK cells, granulocytes, thymocytes, platelets
Function	Lymphocyte activation
Cell Type	B cells, Granulocytes, NK cells, Platelets, T cells, Thymocytes, Tregs
Biology Area	Costimulatory Molecules, Immunology, Innate Immunity
Molecular Family	CD Molecules
Antigen References	<ol style="list-style-type: none">1. Barclay AN, <i>et al.</i> 1997. <i>The Leukocyte Antigen FactsBook</i> Academic Press.2. Testi R, <i>et al.</i> 1994. <i>Immunol. Today</i> 15:479.3. Moretta A, <i>et al.</i> 1991. <i>J. Exp. Med.</i> 174:1393.4. Yokoyama WM, <i>et al.</i> 1988. <i>J. Immunol.</i> 141:369.
Gene ID	12515


Related Protocols

[Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

Biotin anti-mouse CD69, FITC anti-mouse CD69, PE anti-mouse CD69, PE/Cyanine5 anti-mouse CD69, Purified anti-mouse CD69, PE/Cyanine7 anti-mouse CD69, APC anti-mouse CD69, Alexa Fluor® 488 anti-mouse CD69, Alexa Fluor® 647 anti-mouse CD69, PerCP anti-mouse CD69, PerCP/Cyanine5.5 anti-mouse CD69, Pacific Blue™ anti-mouse CD69, Brilliant Violet 421™ anti-mouse CD69, APC/Cyanine7 anti-mouse CD69, Brilliant Violet 605™ anti-mouse CD69, Brilliant Violet 510™ anti-mouse CD69, Purified anti-mouse CD69 (Maxpar® Ready), PE/Dazzle™ 594 anti-mouse CD69, Brilliant Violet 711™ anti-mouse CD69, Alexa Fluor® 700 anti-mouse CD69, Brilliant Violet 650™ anti-mouse CD69, Brilliant Violet 785™ anti-mouse CD69, TotalSeq™-A0197 anti-mouse CD69, APC/Fire™ 750 anti-mouse CD69, TotalSeq™-C0197 anti-mouse CD69, TotalSeq™-B0197 anti-mouse CD69, KIRAVIA Blue 520™ anti-mouse CD69, Spark NIR™ 685 anti-mouse CD69, Spark Red™ 718 anti-mouse CD69

Product Data

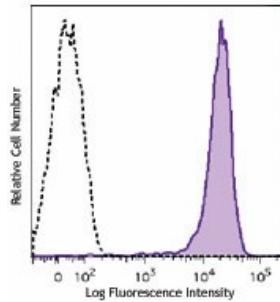

PMA+ionomycin-stimulated (6 hours)
C57BL/6 mouse splenocytes were
stained with CD69 (clone H1.2F3)
Brilliant Violet 605™ (filled histogram).

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587