

## PE anti-human Perforin Antibody

<b>Catalog# / Size</b>	308105 / 25 tests 308106 / 100 tests
<b>Clone</b>	dG9
<b>Regulatory Status</b>	RUO
<b>Other Names</b>	PRF1, P1, PFP, HPLH2
<b>Isotype</b>	Mouse IgG2b, κ
<b>Description</b>	Perforin is a 70 kD cytolytic protein that is expressed in the cytoplasmic granules of cytotoxic T lymphocytes (CTLs) and natural killer (NK) cells. Perforin is one of the major effector molecules used by cytotoxic T cells and NK cells to mediate targeted cell lysis.

### Product Details

---

<b>Verified Reactivity</b>	Human
<b>Reported Reactivity</b>	Cow
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Mouse
<b>Immunogen</b>	Purified granules from the human lymphoma cell line
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
<b>Preparation</b>	The antibody was purified by affinity chromatography, and conjugated with PE under optimal conditions.
<b>Concentration</b>	Lot-specific (to obtain lot-specific concentration, please enter the lot number in our <a href="#">Concentration and Expiration Lookup</a> or <a href="#">Certificate of Analysis</a> online tools.)
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">ICFC - Quality tested</a>
<b>Recommended Usage</b>	Each lot of this antibody is quality control tested by <a href="#">intracellular immunofluorescent staining with flow cytometric analysis</a> . For flow cytometric staining, the suggested use of this reagent is 5 µl per million cells in 100 µl staining volume or 5 µl per 100 µl of whole blood.
<b>Excitation Laser</b>	Blue Laser (488 nm) Green Laser (532 nm)/Yellow-Green Laser (561 nm)
<b>Application Notes</b>	Clone dG9 primarily recognizes perforin associated with cytotoxic granules <sup>9</sup> . Additional reported applications (for the relevant formats) include: immunoprecipitation, intracellular flow cytometric analysis and immunofluorescence microscopy <sup>5,7</sup> , and immunohistochemical staining of acetone-fixed frozen tissue sections and formalin-fixed paraffin-embedded tissue sections <sup>1,4</sup> .  Does not cross-react with mouse <sup>1</sup> .
<b>Application References</b>	<ol style="list-style-type: none"> <li>1. Hameed A, <i>et al.</i> 1992. <i>Am. J. Pathol.</i> 140:1025. (IHC)</li> <li>2. Schaerli P, <i>et al.</i> 2004. <i>J. Exp. Med.</i> 199:1265.</li> <li>3. Watanabe N, <i>et al.</i> 1997. <i>Blood</i> 90:3662.</li> <li>4. Mauad T, <i>et al.</i> 2004. <i>Pediatr. Pulmonol.</i> 38:233. (IHC)</li> <li>5. Barrat FJ, <i>et al.</i> 1999. <i>P. Natl. Acad. Sci. USA</i> 96:8645. (IF)</li> <li>6. Chen H, <i>et al.</i> 2005. <i>J. Immunol.</i> 175:591.</li> <li>7. Bryceson YT, <i>et al.</i> 2007. <i>Blood</i> doi:10.1182/blood-2007-02-074468. (IF)</li> <li>8. Wood SM, <i>et al.</i> 2009. <i>Blood</i> 114:4117. <a href="#">PubMed</a></li> <li>9. Makedonas G, <i>et al.</i> 2010. <i>PLoS Pathog.</i> 6:e1000798.</li> </ol>
<b>(PubMed link indicates BioLegend citation)</b>	

### Product Citations

1. Tang-Huau TL, *et al.* 2018. *Nat Commun.* 9:2570. [PubMed](#)
2. Seo IH, *et al.* 2021. *Immune Netw.* 21:e17. [PubMed](#)
3. Hejazi M, *et al.* 2022. *Front Immunol.* 12:798087. [PubMed](#)
4. Su S *et al.* 2018. *Cell.* 175(2):442-457. [PubMed](#)
5. Hegewisch-Solloa E, *et al.* 2021. *J Immunol.* 207:950. [PubMed](#)
6. Lettau M, *et al.* 2018. *Int Immunol.* 30:215. [PubMed](#)
7. Georg P, *et al.* 2022. *Cell.* 185:493. [PubMed](#)
8. Kuttruff S, *et al.* 2009. *Blood.* 113:358. [PubMed](#)
9. Van der Meer JM, *et al.* 2020. *Cancer Immunol Immunother.* . [PubMed](#)
10. Hu X, *et al.* 2020. *Neoplasia.* 1.290972222. [PubMed](#)
11. Van der Meer JMR, *et al.* 2021. *Oncoimmunology.* 10:1981049. [PubMed](#)
12. Liu CF, *et al.* 2017. *Front Immunol.* 1.445833333. [PubMed](#)
13. Bertino P, *et al.* 2014. *Vaccine.* 32:1670. [PubMed](#)

**RRID** AB\_314703 (BioLegend Cat. No. 308105)  
 AB\_314704 (BioLegend Cat. No. 308106)

## Antigen Details

<b>Structure</b>	70 kD
<b>Distribution</b>	CTL, NK (cytoplasmic granules)
<b>Function</b>	Mediates targeted cell lysis
<b>Cell Type</b>	NK cells, T cells
<b>Biology Area</b>	Cell Biology, Immunology, Innate Immunity, Neuroscience
<b>Molecular Family</b>	Cytokines/Chemokines
<b>Antigen References</b>	1. Lieberman J. 2003. <i>Nat. Rev. Immunol.</i> 3:361. 2. Trapani J, <i>et al.</i> 2002. <i>Nat. Rev. Immunol.</i> 2:735.
<b>Gene ID</b>	<a href="#">5551</a>

## Related Protocols


[Surface and Intracellular Cytokine Staining for Flow Cytometry - Video](#)

[Intracellular Flow Cytometry Staining Protocol](#)

## Other Formats

FITC anti-human Perforin, PE anti-human Perforin, Purified anti-human Perforin, Alexa Fluor® 488 anti-human Perforin, Alexa Fluor® 647 anti-human Perforin, APC anti-human Perforin, Pacific Blue™ anti-human Perforin, PerCP/Cyanine5.5 anti-human Perforin, Brilliant Violet 510™ anti-human Perforin, Brilliant Violet 421™ anti-human Perforin, Alexa Fluor® 594 anti-human Perforin, APC/Cyanine7 anti-human Perforin, PE/Cyanine7 anti-human Perforin, Brilliant Violet 711™ anti-human Perforin, PE/Dazzle™ 594 anti-human Perforin

## Product Data

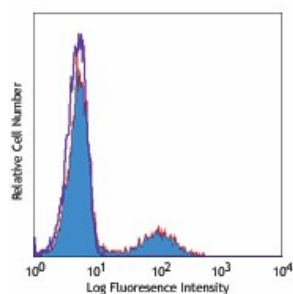

Human peripheral blood lymphocytes intracellularly stained with dG9 PE

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

\*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, [www.biolegend.com/ordering#license](http://www.biolegend.com/ordering#license)). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587