

Alexa Fluor® 647 anti-DYKDDDDK Tag Antibody

Catalog# / Size	637315 / 25 µg 637316 / 100 µg
Clone	L5
Regulatory Status	RUO
Other Names	FLAG tag
Isotype	Rat IgG2a, λ
Description	The DYKDDDDK tag, commonly referred to as Sigma®'s FLAG® Tag, is often used as a protein modification in order to simplify the labeling and detection of proteins. This unique amino acid sequence allows for specific antibody detection in western blotting, immunoprecipitation, and immunostaining techniques. Due to the short sequence, this modification is not likely to affect the structure or function of the modified proteins.

Product Details

Verified Reactivity	DYKDDDDK tag epitope
Antibody Type	Monoclonal
Host Species	Rat
Immunogen	DYKDDDDK-tagged mouse Langerin
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide.
Preparation	The antibody was purified by affinity chromatography and conjugated with Alexa Fluor® 647 under optimal conditions.
Concentration	0.5 mg/mL
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	ICC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunocytochemistry. For immunocytochemistry, a concentration range of 1:250-1:2500 (0.2 - 2 µg/mL) is recommended. It is recommended that the reagent be titrated for optimal performance for each application. * Alexa Fluor® 647 has a maximum emission of 668 nm when it is excited at 633 nm / 635 nm. Alexa Fluor® and Pacific Blue™ are trademarks of Life Technologies Corporation. View full statement regarding label licenses
Excitation Laser	Red Laser (633 nm)
Application Notes	The L5 clone has been demonstrated to have 2-8 fold better sensitivity in WB than another commonly used antibody clone, M2.
Application References	<ol style="list-style-type: none"> 1. Park SH, <i>et al.</i> 2008. <i>J Immunol Methods</i>. 331:27. 2. Moon SH, <i>et al.</i> 2010. <i>J. Biol Chem</i>. 285:12935. PubMed 3. Sasaki M, <i>et al.</i> 2011. <i>J. Biol Chem</i>. 286:39370. PubMed 4. Sonder SU, <i>et al.</i> 2012. <i>J Immunol</i>. 188:5906. PubMed 5. Jiang Y, <i>et al.</i> 2013. <i>Int Immunol</i>. 25:235. PubMed 6. Zuo X, <i>et al.</i> 2014. <i>PLoS One</i>. 9:84748. PubMed 7. Toyo-Oak K, <i>et al.</i> 2014. <i>J Neurosci</i>. 34:12168. PubMed
(PubMed link indicates BioLegend citation)	
Product Citations	<ol style="list-style-type: none"> 1. Weeks AM, <i>et al.</i> 2021. <i>Proc Natl Acad Sci U S A</i>. 118: . PubMed 2. Zost SJ, <i>et al.</i> 2020. <i>Nat Med</i>. 26:1422. PubMed

RRID AB_2716154 (BioLegend Cat. No. 637315)
 AB_2861049 (BioLegend Cat. No. 637316)

Antigen Details

Biology Area Cell Biology

Antigen References 1. Einhauer A. 2001. *J. Biochem. Biophys. Methods.* 49:455.
 2. Knappik A and Pluckthun A. 1994. *Biotechniques.* 17:754.

Gene ID NA


Related Protocols

[Immunocytochemistry Staining Protocol](#)


Other Formats

Purified anti-DYKDDDDK Tag, Anti-DYKDDDDK Tag (L5) Affinity Gel, APC anti-DYKDDDDK Tag, PE anti-DYKDDDDK Tag, Direct-Blot™ HRP anti-DYKDDDDK Tag, Alexa Fluor® 594 anti-DYKDDDDK Tag, Alexa Fluor® 647 anti-DYKDDDDK Tag, Alexa Fluor® 488 anti-DYKDDDDK Tag, PE/Cyanine7 anti-DYKDDDDK Tag, Brilliant Violet 421™ anti-DYKDDDDK Tag, PerCP/Cyanine5.5 anti-DYKDDDDK Tag, Ultra-LEAF™ Purified anti-DYKDDDDK Tag, PE/Dazzle™ 594 anti-DYKDDDDK Tag Antibody, TotalSeq™-B1129 anti-DYKDDDDK Tag, TotalSeq™-A1129 anti-DYKDDDDK Tag, TotalSeq™-C1129 anti-DYKDDDDK Tag

Product Data

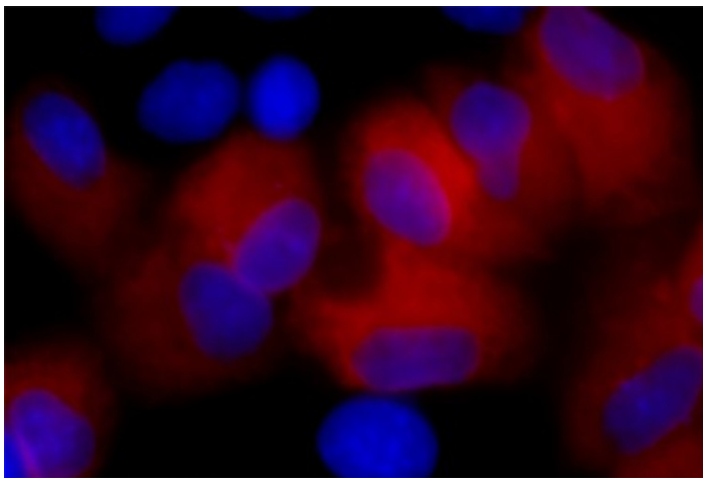

HeLa cells transiently transfected with DYKDDDDK tag fused protein were fixed with 4% paraformaldehyde (PFA) for 15 minutes, permeabilized with 0.5% Triton X-100 for 3 minutes, and blocked with 5% FBS for 60 minutes. Then the cells were intracellularly stained with 1:1000 diluted (0.5 µg/ml) Alexa Fluor® 647 anti-DYKDDDDK Tag antibody (red). Nuclei were counterstained with DAPI (blue, Cat. No. 422801). The image was captured with a 40X objective.

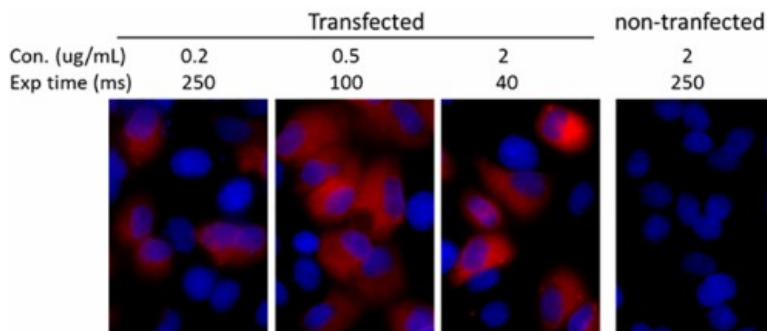

HeLa cells transiently transfected with DYKDDDDK tag fused protein or non-transfected were fixed with 4% paraformaldehyde (PFA) for 15 minutes, permeabilized with 0.5% Triton X-100 for 3 minutes, and blocked with 5% FBS for 60 minutes. Then the cells were intracellularly stained with 1:250 (2 µg/ml), 1:1000 (0.5 µg/ml) and 1:2500 diluted (0.2 µg/ml) Alexa Fluor® 647 anti-DYKDDDDK Tag antibody (red). Nuclei were counterstained with DAPI (blue, Cat. No. 422801). The image was captured with a 40X objective using KEYENCE BZ-X700 fluorescence microscope.

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587