

APC anti-mouse F4/80 Antibody

Catalog# / Size	123115 / 25 µg 123116 / 100 µg
Clone	BM8
Regulatory Status	RUO
Other Names	EMR1, Ly71
Isotype	Rat IgG2a, κ
Description	F4/80 is a 160 kD glycoprotein. It is characterized as a member of the epidermal growth factor (EGF)-transmembrane 7 (TM7) family. F4/80, also known as EMR1 or Ly71, has been widely used as a murine macrophage marker, which is expressed on the majority of tissue macrophages including peritoneal macrophages, macrophages in lung, gut, thymus and red pulp of spleen (but not on the macrophages located in T cell areas of the spleen, lymph node and Peyer's patch), Kuffer cells, Langerhans cells, and bone marrow stromal cells. F4/80 has also been shown on a subset of dendritic cells. The biological ligand of F4/80 has not been identified, but it has been reported that F4/80 is required for induction of CD8 ⁺ T cell-mediated peripheral tolerance.

Product Details

Verified Reactivity	Mouse
Antibody Type	Monoclonal
Host Species	Rat
Immunogen	Murine macrophages
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide.
Preparation	The antibody was purified by affinity chromatography, and conjugated with APC under optimal conditions.
Concentration	0.2 mg/ml
Storage & Handling	The F4/80 antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is ≤ 0.25 µg per 10 ⁶ cells in 100 µl. It is recommended that the reagent be titrated for optimal performance for each application.
Excitation Laser	Red Laser (633 nm)
Application Notes	Additional reported applications (for the relevant formats) include: immunohistochemical staining of acetone-fixed frozen sections ^{1,2} and formalin-fixed paraffin-embedded sections ^{6,7} , Western blotting, and spatial biology (IBEX) ^{12,13} .
Application References	<ol style="list-style-type: none"> Schaller E, <i>et al.</i> 2002. <i>Mol. Cell. Biol.</i> 22:8035. (IHC) Stevceva L, <i>et al.</i> 2001. <i>BMC Clin Pathol.</i> 1:3. (IHC) Kobayashi M, <i>et al.</i> 2008. <i>J. Leukoc. Biol.</i> 83:1354. PubMed Poeckel D, <i>et al.</i> 2009. <i>J. Biol Chem.</i> 284:21077. PubMed Glass AM, <i>et al.</i> 2013. <i>J. Immunol.</i> 190:4830. PubMed Koehm S, <i>et al.</i> 2007. <i>J. Allergy Clin. Immunol.</i> 120:570. (IHC) Rankin AL, <i>et al.</i> 2010. <i>J. Immunol.</i> 184:1526. (IHC) Sasi SP, <i>et al.</i> 2014. <i>J Biol Chem.</i> 289:14178. PubMed Thakus VS, <i>et al.</i> 2014. <i>Toxicol Lett.</i> 230:322. PubMed Watson NB, <i>et al.</i> 2015. <i>J Immunol.</i> 194:2796. PubMed Hirakawa H, <i>et al.</i> 2015. <i>PLoS One.</i> 10:119360. PubMed
(PubMed link indicates BioLegend citation)	

12. Radtke AJ, *et al.* 2020. *Proc Natl Acad Sci U S A.* 117:33455-65. (SB) [PubMed](#)
13. Radtke AJ, *et al.* 2022. *Nat Protoc.* 17:378-401. (SB) [PubMed](#)

Product Citations

1. Huang W, *et al.* 2015. *J Allergy Clin Immunol.* [PubMed](#)
2. Barrett N, *et al.* 2016. *Cell Rep.* 16: 1039-1054. [PubMed](#)
3. Liu Y, *et al.* 2021. *Cell Metabolism.* 33(6):1221-1233.e11. [PubMed](#)
4. Ying W, *et al.* 2019. *Cell Metab.* 29:457. [PubMed](#)
5. Saha D *et al.* 2017. *Cancer cell.* 32(2):253-267. [PubMed](#)
6. Zhao Z, *et al.* 2021. *Nat Commun.* 12:4355. [PubMed](#)
7. Marchelletta RR, *et al.* 2021. *J Clin Invest.* 131:. [PubMed](#)
8. Nguyen T, *et al.* 2021. *STAR Protoc.* 2:101004. [PubMed](#)
9. Grandjean CL, *et al.* 2021. *Sci Adv.* 7:. [PubMed](#)
10. Deng R, *et al.* 2022. *Adv Sci (Weinh).* 9:e2103343. [PubMed](#)
11. Yang D, *et al.* 2022. *Cell.* 185:1905. [PubMed](#)
12. Zheng J, *et al.* 2022. *Front Pharmacol.* 13:866228. [PubMed](#)
13. Furrer R, *et al.* 2021. *Sci Adv.* 7:eabi4852. [PubMed](#)
14. Li Z, *et al.* 2022. *Nat Commun.* 13:1845. [PubMed](#)
15. Platt DJ, *et al.* 2021. *Cell Reports.* 35(6):109113. [PubMed](#)
16. Du J, *et al.* 2020. *Developmental Cell.* 55(6):737-753.e7. [PubMed](#)
17. Gerwing M, *et al.* 2020. *Mol Imaging Biol.* 1.959027778. [PubMed](#)
18. Bousset L, *et al.* 2020. *PLoS Biol.* 18:e3000948. [PubMed](#)
19. Yu J, *et al.* 2020. *Sci Rep.* 10:21360. [PubMed](#)
20. Sharma-Chawla N, *et al.* 2016. *Infect Immun.* [PubMed](#)
21. Emal D, *et al.* 2016. *J Am Soc Nephrol.* 10.1681/ASN.2016030255. [PubMed](#)
22. Dale M, *et al.* 2016. *J Immunol.* 196: 4536 - 4543. [PubMed](#)
23. Kawano H, *et al.* 2016. *Int Immunol.* 28: 489 - 501. [PubMed](#)
24. Jaiswal M, *et al.* 2015. *Sci Rep.* 5: 15221. [PubMed](#)
25. Bieren J, *et al.* 2015. *J Immunol.* 194:1154. [PubMed](#)
26. Prantner D, *et al.* 2009. *Infect Immun.* 77:5334. [PubMed](#)
27. Cullis J, *et al.* 2017. *Cancer Immunol Res.* 5:182. [PubMed](#)
28. Kim SE, *et al.* 2019. *Nutrients.* 11. [PubMed](#)
29. Chang T, *et al.* 2017. *PLoS Pathog.* 10.1371/journal.ppat.1006485. [PubMed](#)
30. Yan X, *et al.* 2022. *Commun Biol.* 5:320. [PubMed](#)
31. Li Y, *et al.* 2022. *Int J Biol Sci.* 18:3697. [PubMed](#)
32. He R, *et al.* 2021. *J Immunol.* 206:2353. [PubMed](#)
33. Suresh R, *et al.* 2020. *J Immunother Cancer.* 8:. [PubMed](#)
34. Symeonidou V, *et al.* 2021. *Cell Rep.* 37:109900. [PubMed](#)
35. Zhang R, *et al.* 2021. *Cell Mol Immunol.* 18:1222. [PubMed](#)
36. Kurosawa M, *et al.* 2021. *Int J Mol Sci.* 22:. [PubMed](#)
37. He F, *et al.* 2021. *J Exp Med.* 218:. [PubMed](#)
38. Stubelius A, *et al.* 2012. *Immunobiology.* 217:751. [PubMed](#)
39. Sjöstrand M, *et al.* 2013. *J Immunol.* 191:3753. [PubMed](#)
40. Li X, *et al.* 2015. *Proc Natl Acad Sci U S A.* 112:6080. [PubMed](#)
41. Tan X, *et al.* 2016. *PLoS One.* 11: 0160985. [PubMed](#)
42. Yilmaz B, *et al.* 2021. *Cell Host Microbe.* 29(4):650-663.e9. [PubMed](#)
43. Dong L, *et al.* 2021. *Cancer Cell.* [PubMed](#)
44. Gonçalves NP, *et al.* 2017. *J Neuroinflammation.* 14:115. [PubMed](#)
45. Yates C, *et al.* 2017. *Stem Cell Res Ther.* 10.1186/s13287-017-0644-9. [PubMed](#)
46. Lancaster JN, *et al.* 2019. *Nat Commun.* 10:2220. [PubMed](#)
47. Chou YJ, *et al.* 2020. *Sci Rep.* 10:8422. [PubMed](#)
48. Lee H, *et al.* 2020. *Cell Metabolism.* 31(4):822-836. [PubMed](#)
49. Jiang W, *et al.* 2021. *Oncol Lett.* 22:625. [PubMed](#)
50. Liang S, *et al.* 2022. *Acta Pharm Sin B.* 12:2494. [PubMed](#)
51. Sepe JJ, *et al.* 2022. *JACC Basic Transl Sci.* 7:915. [PubMed](#)
52. Nam J, *et al.* 2021. *Adv Sci (Weinh).* 2002577:8. [PubMed](#)
53. de Mingo Pulido, *et al.* 2021. *Immunity.* 54(6):1154-1167.e7. [PubMed](#)
54. Zheng H, *et al.* 2021. *Frontiers in Cell and Developmental Biology.* 9:641527. [PubMed](#)
55. Xu Y, *et al.* 2020. *Mol Ther.* 28:202. [PubMed](#)
56. Mia MM, *et al.* 2020. *PLoS Biol.* 18:e3000941. [PubMed](#)
57. Xing J, *et al.* 2016. *Nat Immunol.* 17:1373-1380. [PubMed](#)
58. Zhu J, *et al.* 2016. *Sci Rep.* 6:27136. [PubMed](#)
59. Sun P, *et al.* 2015. *J Hepatol.* [PubMed](#)
60. Shen H, *et al.* 2015. *Infect Immun.* 83: 2694 - 2704. [PubMed](#)
61. Gritsina G, *et al.* 2015. *Mol Cancer Ther.* 14:1035. [PubMed](#)
62. Wong E, *et al.* 2019. *Cell Rep.* 29:3047. [PubMed](#)
63. Wu R, *et al.* 2019. *J Cell Mol Med.* 24:1684. [PubMed](#)
64. Álvaro de Mingo Pulido *et al.* 2018. *Cancer cell.* 33(1):60-74. [PubMed](#)
65. Neumann AM *et al.* 2017. *PloS one.* 12(11):e0187330. [PubMed](#)
66. Han J, *et al.* 2022. *FEBS J.* 289:417. [PubMed](#)
67. Zhao J, *et al.* 2022. *J Nanobiotechnology.* 20:62. [PubMed](#)
68. Liang F, *et al.* 2021. *Front Immunol.* 12:658715. [PubMed](#)
69. Rybalko V, *et al.* 2015. *PLoS One.* 10: 0145550. [PubMed](#)
70. Ohtsuki T, *et al.* 2016. *J Virol.* 90: 300 - 307. [PubMed](#)
71. Tajfirouz D, *et al.* 2017. *Virology.* 10.1016/j.virol.2017.02.022. [PubMed](#)
72. Krivanek J, *et al.* 2020. *Nat Commun.* 3.802777778. [PubMed](#)
73. Brestoff JR, *et al.* 2021. *Cell Metab.* 33:270. [PubMed](#)
74. Evren E, *et al.* 2020. *Immunity.* 54(2):259-275.e7. [PubMed](#)
75. Esmaili S, *et al.* 2021. *Cell Systems.* 12(5):432-445.e7. [PubMed](#)

76. Jeong D, *et al.* 2021. *eLife*. 0.416666666666667. [PubMed](#)
77. Molina-Sanchez P, *et al.* 2020. *Gastroenterology*. 2203:159. [PubMed](#)
78. Qi F, *et al.* 2017. *Molecular Immunology*. 10.1016/j.molimm.2017.05.008. [PubMed](#)
79. Stefan K, *et al.* 2017. *PLoS One*. . 10.1371/journal.pone.0188600. [PubMed](#)
80. Sinha M, *et al.* 2018. *Nat Commun*. 9:936. [PubMed](#)
81. Grohmann M *et al.* 2018. *Cell*. 175(5):1289-1306 . [PubMed](#)
82. Kosack L, *et al.* 2017. *Sci Rep*. 7:11289. [PubMed](#)
83. Christopher S Garris *et al.* 2018. *Immunity*. 49(6):1148-1161 . [PubMed](#)
84. Duarte N, *et al.* 2018. *Hepatol Commun*. 0.8333333333. [PubMed](#)
85. Zhong J, *et al.* 2018. *Sci Adv*. 4:eaas9864. [PubMed](#)
86. Sun Q, *et al.* 2022. *Cell Death Dis*. 13:724. [PubMed](#)
87. He C, *et al.* 2021. *Oncoimmunology*. 1897295:10. [PubMed](#)
88. Jiao H, *et al.* 2021. *Cell*. 184(11):2896-2910.e13. [PubMed](#)
89. Shang Y, *et al.* 2016. *Nat Immunol*. 10.1038/ni.3486. [PubMed](#)
90. Lasky C, *et al.* 2015. *Infect Immun*. 83: 2627 - 2635. [PubMed](#)
91. Peng L, *et al.* 2020. *Transl Lung Cancer Res*. 0.892361111. [PubMed](#)
92. Vannini N, *et al.* 2019. *Cell Stem Cell*. 24:405. [PubMed](#)
93. Nguyen CM, *et al.* 2019. *Diabetes*. 68:1499. [PubMed](#)
94. Dey A *et al.* 2019. *The EMBO journal*. 38(7) pii: e100293. [PubMed](#)
95. Cao W, *et al.* 2017. *Immunity*. 47:1182. [PubMed](#)
96. Wang Y, *et al.* 2022. *Front Immunol*. 13:767630. [PubMed](#)
97. Hou H, *et al.* 2021. *Front Cell Dev Biol*. 9:737003. [PubMed](#)
98. Kosaric N, *et al.* 2020. *Mol Ther*. 28:2007. [PubMed](#)
99. Ling X, *et al.* 2021. *ACS Nano*. 15:765. [PubMed](#)
100. Gordon RA, *et al.* 2020. *PLoS One*. 15:e0226396. [PubMed](#)
101. Helmi N, *et al.* 2015. *J Infect Dis*. 211:1677. [PubMed](#)
102. Tsoyi K, *et al.* 2015. *J Leukoc Biol*. 97:171. [PubMed](#)
103. Ma Y, *et al.* 2015. *Cardiovasc Res*. 106:421. [PubMed](#)
104. Ben-Shaan T, *et al.* 2016. *Nat Med*. 10.1038/nm.4133. [PubMed](#)
105. Karki R, *et al.* 2020. *JCI Insight*. 5:00. [PubMed](#)
106. Chi Z, *et al.* 2020. *Molecular Cell*. 80(1):43-58.e7. [PubMed](#)
107. Sun L, *et al.* 2020. *J Immunol*. 2177:204. [PubMed](#)
108. Müller A, *et al.* 2020. *J Clin Invest*. 5765:130. [PubMed](#)
109. Parada-Kusz M, *et al.* 2018. *Dis Model Mech*. 11:. [PubMed](#)
110. Luan HH *et al.* 2019. *Cell*. 178(5):1231-1244 . [PubMed](#)
111. Shimokawa C *et al.* 2017. *Immunity*. 46(5):863-874 . [PubMed](#)
112. Giri J, *et al.* 2020. *Cell Rep*. 30:1923. [PubMed](#)
113. Galvão I, *et al.* 2018. *Front Immunol*. 9:142. [PubMed](#)
114. Alexander RK, *et al.* 2020. *eLife*. 9:e54090.. [PubMed](#)
115. Prantner D, *et al.* 2011. *Infect Immun*. 79:3922. [PubMed](#)
116. Jeong SJ, *et al.* 2021. *Autophagy*. 17:3740. [PubMed](#)
117. Guo W, *et al.* 2022. *J Immunother Cancer*. 10:. [PubMed](#)
118. Shin JE, *et al.* 2021. *Heliyon*. 7:e08433. [PubMed](#)
119. Rao X, *et al.* 2022. *Cell Death Dis*. 13:891. [PubMed](#)
120. Carozza JA, *et al.* 2020. *Nat Cancer*. 184:1. [PubMed](#)
121. Agrawal V, *et al.* 2016. *J Leukoc Biol*. 100: 261 - 274. [PubMed](#)
122. Amash A, *et al.* 2016. *J Immunol*. 196: 3331-3340. [PubMed](#)
123. Gereke M, *et al.* 2012. *J Vis Exp*. 70: 4322. [PubMed](#)
124. Sullivan B, *et al.* 2015. *PLoS Pathog*. 11:1004588. [PubMed](#)
125. Oatley J, *et al.* 2009. *Development*. 136:1191. [PubMed](#)
126. Frazier T, *et al.* 2018. *Adipocyte*. 0.41875. [PubMed](#)
127. An W *et al.* 2019. *British journal of pharmacology*. 176(5):699-710 . [PubMed](#)
128. Yu Y, *et al.* 2019. *J Clin Invest*. 129:546. [PubMed](#)
129. Linde N, *et al.* 2018. *Nat Commun*. 9:21. [PubMed](#)
130. Anstine L, *et al.* 2017. *J Am Heart Assoc.*. 10.1161/JAHA.117.007097. [PubMed](#)
131. Zein HS, *et al.* 2022. *Sci Rep*. 12:4834. [PubMed](#)
132. Luo W, *et al.* 2022. *Front Immunol*. 13:816761. [PubMed](#)
133. Lv F, *et al.* 2022. *iScience*. 25:103639. [PubMed](#)
134. Wu L, *et al.* 2022. *Theranostics*. 12:842. [PubMed](#)
135. Stump CT, *et al.* 2021. *Open Biol*. 11:210245. [PubMed](#)
136. Jiao Y, *et al.* 2021. *Crit Care*. 25:356. [PubMed](#)
137. Gao H, *et al.* 2021. *Cells*. 10: [PubMed](#)
138. Glauser D, *et al.* 2013. *J Virol*. 87:10828. [PubMed](#)
139. Bretscher P, *et al.* 2015. *EMBO Mol Med*. 7:593. [PubMed](#)
140. Hara-Chikuma H 2016. *Mol Cell Biol*. 36: 1206-1218. [PubMed](#)
141. Hu J, *et al.* 2020. *Theranostics*. 7.154166667. [PubMed](#)
142. DeVito NC, *et al.* 2021. *Cell Reports*. 35(5):109071. [PubMed](#)
143. Gozgit JM, *et al.* 2021. *Cancer Cell*. . [PubMed](#)
144. Chen F *et al.* 2018. *Cell reports*. 25(10):2775-2783 . [PubMed](#)
145. Liao T, *et al.* 2017. *Front Immunol*. 8:1334. [PubMed](#)
146. Deng M, *et al.* 2020. *Nat Commun*. 11:2193. [PubMed](#)
147. Mao FY, *et al.* 2021. *Cell Mol Gastroenterol Hepatol*. 12:395. [PubMed](#)
148. Shi R, *et al.* 2022. *Theranostics*. 12:875. [PubMed](#)
149. Du Y, *et al.* 2022. *Nat Commun*. 13:231. [PubMed](#)
150. Tartey S, *et al.* 2021. *Cancer Research*. 81(9):2358-2372. [PubMed](#)
151. Cheng Y, *et al.* 2020. *PLoS Pathog*. 16:e1008569. [PubMed](#)
152. Wang W, *et al.* 2018. *Oncogenesis*. 7:97. [PubMed](#)
153. Xie Z, *et al.* 2018. *J Am Heart Assoc*. 7:e007442. [PubMed](#)
154. Pinho S, *et al.* 2022. *Nat Cell Biol*. 24:290. [PubMed](#)
155. Lin S, *et al.* 2022. *Nat Commun*. 13:3301. [PubMed](#)

156. Ye P, *et al.* 2022. *Front Cardiovasc Med.* 8:810477. [PubMed](#)
157. Zhang Y, *et al.* 2021. *J Immunol Res.* 2021:5599439. [PubMed](#)
158. Fairfax K, *et al.* 2013. *J Immunol.* 191:4202. [PubMed](#)
159. Ericson J, *et al.* 2014. *PLoS One.* 9:108553. [PubMed](#)
160. Toda G, *et al.* 2021. *STAR Protoc.* 2:100246. [PubMed](#)
161. Guldner IH, *et al.* 2020. *Cell.* 183(5):1234-1248.e25. [PubMed](#)
162. Engler AE, *et al.* 2020. *Cell Reports.* 33(13):108553. [PubMed](#)
163. Fraccarollo D, *et al.* 2021. *eLife.* 0.4166666666666667. [PubMed](#)
164. Wang T, *et al.* 2018. *Front Physiol.* 9:209. [PubMed](#)
165. Xiang W, *et al.* 2018. *Nat Commun.* 9:2574. [PubMed](#)
166. Moon HG *et al.* 2018. *Immunity.* 49(2):275-287. [PubMed](#)
167. LaFleur MW, *et al.* 2019. *Nat Commun.* 10:1668. [PubMed](#)
168. Hu J, *et al.* 2019. *Mol Ther Nucleic Acids.* 16:650. [PubMed](#)
169. Matundan H, *et al.* 2019. *J Virol.* 93. [PubMed](#)
170. Yu L, *et al.* 2019. *Elife.* 8:e46677. [PubMed](#)
171. Ruiz de Azua I, *et al.* 2017. *J Clin Invest.* 127:4148. [PubMed](#)
172. Li Y, *et al.* 2020. *Cell Reports.* 27(10):2809-2816.e3. [PubMed](#)
173. Pfefferlé M, *et al.* 2021. *Front Immunol.* 12:680855. [PubMed](#)
174. Wara AK, *et al.* 2020. *Cell Rep.* 33:108550. [PubMed](#)
175. Fu Z, *et al.* 2021. *Cell Res.* 31:631. [PubMed](#)
176. Yang Y, *et al.* 2022. *Elife.* 11. [PubMed](#)
177. Frenis K, *et al.* 2021. *Basic Res Cardiol.* 31:116. [PubMed](#)
178. Lei J, *et al.* 2021. *Theranostics.* 4251:11. [PubMed](#)
179. Dewi IMW, *et al.* 2021. *Cell Reports Medicine.* 2(5):100289. [PubMed](#)
180. Oguri Y, *et al.* 2020. *Cell.* 182(3):563-577.e20. [PubMed](#)
181. Ni J, *et al.* 2020. *Immunity.* 52(6):1075-1087.e8. [PubMed](#)
182. Madenspacher JH, *et al.* 2020. *JCI Insight.* 5:00. [PubMed](#)
183. Kostadinova E, *et al.* 2016. *Sci Rep.* 6:30943. [PubMed](#)
184. Tan X, *et al.* 2015. *Sci Rep.* 5: 15702. [PubMed](#)
185. Kuwahara G, *et al.* 2014. *PLoS One.* 9:106987. [PubMed](#)
186. Patankar Y, *et al.* 2013. *Infect Immun.* 81:2043. [PubMed](#)
187. Kyburz A, *et al.* 2019. *J Allergy Clin Immunol.* 143:1496. [PubMed](#)
188. Liu C *et al.* 2019. *Immunity.* 51(2):381-397. [PubMed](#)
189. Simula L *et al.* 2018. *Cell reports.* 25(11):3059-3073. [PubMed](#)
190. Tomay F, *et al.* 2019. *J Transl Med.* 17:237. [PubMed](#)
191. Zhou J, *et al.* 2019. *Nat Commun.* 10:2427. [PubMed](#)
192. Brown SM *et al.* 2017. *Endocrinology.* 158(10):3592-3604. [PubMed](#)
193. Nam GH, *et al.* 2018. *Nat Commun.* 9:2165. [PubMed](#)
194. Zheng X, *et al.* 2018. *J Clin Invest.* 128:2104. [PubMed](#)
195. Chung EJ, *et al.* 2022. *Aging (Albany NY).* 14:7692. [PubMed](#)
196. Puigdelloses M, *et al.* 2021. *J Immunother Cancer.* 9. [PubMed](#)
197. Lang V, *et al.* 2021. *Elife.* 10. [PubMed](#)
198. Jaggi U, *et al.* 2021. *PLoS Pathog.* 17:e1009999. [PubMed](#)
199. Rashid MH, *et al.* 2021. *Oncol Rep.* 45:1171. [PubMed](#)
200. Tapia-Abellán A, *et al.* 2021. *Sci Adv.* 7:eabf4468. [PubMed](#)
201. Henrich IC, *et al.* 2021. *Cancer Res.* 81:2171. [PubMed](#)
202. Kim K, *et al.* 2014. *Proc Natl Acad Sci U S A.* 111:11774. [PubMed](#)
203. Hu J, *et al.* 2015. *J Hepatol.* Available online 31 August 2015. [PubMed](#)
204. Terrazas M 2014. *Int J Parasit.* 44 613. [PubMed](#)
205. Yang C, *et al.* 2016. *Sci Rep.* 6:34643. [PubMed](#)
206. Jiang Q, *et al.* 2020. *JCI Insight.* 5:00. [PubMed](#)
207. Kesavardhana S, *et al.* 2020. *J Biol Chem.* 295:8325. [PubMed](#)
208. Cao Q, *et al.* 2019. *Theranostics.* 5.576388889. [PubMed](#)
209. Ye ZH, *et al.* 2021. *Transl Oncol.* 14:101162. [PubMed](#)
210. Ito T, *et al.* 2021. *Nat Commun.* 12:2232. [PubMed](#)
211. Wuggenig P, *et al.* 2020. *Commun Biol.* 3:130. [PubMed](#)
212. Xiao Y, *et al.* 2022. *Nat Commun.* 13:758. [PubMed](#)
213. Ali S, *et al.* 2021. *PLoS One.* 16:e0246646. [PubMed](#)
214. Yu W, *et al.* 2022. *Nat Commun.* 13:3544. [PubMed](#)
215. Li Q, *et al.* 2022. *Front Cardiovasc Med.* 9:872102. [PubMed](#)
216. Pandit M, *et al.* 2021. *Eur J Immunol.* 51:1461. [PubMed](#)
217. Li D, *et al.* 2022. *Emerg Microbes Infect.* 11:2248. [PubMed](#)
218. Ishikura S, *et al.* 2021. *J Cell Biochem.* [PubMed](#)
219. Wang G, *et al.* 2021. *Cell Host Microbe.* 29(5):777-791.e6. [PubMed](#)
220. Desai P, *et al.* 2021. *Cell.* 184(5):1214-1231.e16. [PubMed](#)
221. Chen Z, *et al.* 2016. *Nat Commun.* 7:11596. [PubMed](#)
222. Deng W, *et al.* 2015. *Clin Immunol.* [PubMed](#)
223. Honda H, *et al.* 2014. *J Leukoc Biol.* 96:1087. [PubMed](#)
224. Shaw O, *et al.* 2014. *Rheumatology.* 53:1901. [PubMed](#)
225. Rosenbaum SR, *et al.* 2020. *Cell Rep.* 30:510. [PubMed](#)
226. Yuan X, *et al.* 2017. *Elife.* 6:e29540. [PubMed](#)
227. Petkevicius K, *et al.* 2019. *Elife.* 8:e47990. [PubMed](#)
228. Cianciaruso C, *et al.* 2020. *Cell Reports.* 27(10):3062-3080.e11. [PubMed](#)
229. Engblom C, *et al.* 2017. *Science.* 358:6367. [PubMed](#)
230. Zhao F *et al.* 2018. *Immunity.* 48(1):147-160. [PubMed](#)
231. Lima-Junior D, *et al.* 2017. *J Immunol.* 10.4049/jimmunol.1700258. [PubMed](#)
232. Pascual-Pasto G, *et al.* 2022. *J Immunother Cancer.* 10. [PubMed](#)
233. Matsumura T, *et al.* 2022. *Nat Commun.* 13:7064. [PubMed](#)
234. Chen F, *et al.* 2022. *Cell Rep.* 38:110215. [PubMed](#)
235. Wang S, *et al.* 2021. *Commun Biol.* 22. [PubMed](#)

236. Bent EH, *et al.* 2021. Nat Commun. 12:6218. [PubMed](#)
237. Calle P, *et al.* 2021. Cells. 10:. [PubMed](#)
238. Horsle B 2013. Development. 140:1517. [PubMed](#)
239. Calderon B, *et al.* 2015. J Exp Med. 212: 1497-1512. [PubMed](#)
240. Ding H, *et al.* 2016. Nat Commun. 7:11533. [PubMed](#)
241. Gamradt P, *et al.* 2016. PLoS Pathog. 12:e1006032. [PubMed](#)
242. Wang X, *et al.* 2016. J Immunol. 197(12):4750-4761. [PubMed](#)
243. Zhang H, *et al.* 2017. Immunol Cell Biol. 10.1038/icb.2017.11. [PubMed](#)
244. Tian Y, *et al.* 2020. JCI Insight. 5:00. [PubMed](#)
245. Kerkhofs D, *et al.* 2020. Theranostics. 7.022222222. [PubMed](#)
246. Toda G, *et al.* 2020. Molecular Cell. 79(1):43-53.e4. [PubMed](#)
247. Karki R, *et al.* 2020. bioRxiv. . [PubMed](#)
248. Di Genua C, *et al.* 2019. Haematologica. 104:2215. [PubMed](#)
249. Hua X, *et al.* 2021. Front Immunol. 12:706027. [PubMed](#)
250. Liu J, *et al.* 2021. Front Immunol. 12:733808. [PubMed](#)
251. Feng L, *et al.* 2022. Signal Transduct Target Ther. 7:21. [PubMed](#)
252. Zhang MS, *et al.* 2022. Nat Commun. 13:954. [PubMed](#)
253. Pan Y, *et al.* 2021. Br J Pharmacol. 178:4907. [PubMed](#)
254. Chen C, *et al.* 2020. Cell Rep. 2136:30. [PubMed](#)
255. Li K, *et al.* 2020. J Cell Mol Med. 24:13115. [PubMed](#)
256. Karki R, *et al.* 2021. Cell. 184:149. [PubMed](#)
257. Cao X, *et al.* 2017. Neurosci Lett. 10.1016/j.neulet.2017.02.045. [PubMed](#)
258. Jeurissen M, *et al.* 2016. PLoS One. 11:e0167199. [PubMed](#)
259. Montero-Barrera D, *et al.* 2015. BioMed Res Int. 2015 615865. [PubMed](#)
260. Neimert-Andersson T, *et al.* 2011. Vaccine. 29:8965. [PubMed](#)
261. Sun L, *et al.* 2009. FASEB J. 3.602777778. [PubMed](#)
262. Esposito E, *et al.* 2019. Nat Commun. 4.101388889. [PubMed](#)
263. Ji L *et al.* 2019. Immunity. 51(2):272-284 . [PubMed](#)
264. Zhu K, *et al.* 2022. EMBO Rep. 23:e54499. [PubMed](#)
265. Ding R, *et al.* 2022. J Neuroinflammation. 19:137. [PubMed](#)
266. Liu R, *et al.* 2020. Eur J Vasc Endovasc Surg. 60:254. [PubMed](#)
267. Zhang W, *et al.* 2021. Front Immunol. 12:768435. [PubMed](#)
268. Panagi M, *et al.* 2020. Theranostics. 10:1910. [PubMed](#)
269. Zhou H, *et al.* 2021. Cell Death Discov. 7:332. [PubMed](#)
270. Gu H, *et al.* 2021. Elife. 10:. [PubMed](#)
271. Nishio H, *et al.* 2020. Regen Ther. 14:262. [PubMed](#)
272. Mirsoian A, *et al.* 2014. J Exp Med. 211:2373. [PubMed](#)
273. Yin X, *et al.* 2015. J Immunol. 194:4029. [PubMed](#)
274. Kuo P, *et al.* 2016. J Am Heart Assoc. 5: 002610. [PubMed](#)
275. Santopaolo M, *et al.* 2021. Front Immunol. 609406:12. [PubMed](#)
276. Wang Z *et al.* 2019. Cell reports. 26(5):1273-1285 . [PubMed](#)
277. Nishino T, *et al.* 2018. Arterioscler Thromb Vasc Biol. 38:2460. [PubMed](#)
278. Meng Michelle Xu *et al.* 2017. Immunity. 47(2):363-373 . [PubMed](#)
279. Singh KS, *et al.* 2020. Nat Commun. 0.786805556. [PubMed](#)
280. Chao Y, *et al.* 2020. Sci Adv. 6:eaaz4204. [PubMed](#)
281. Chan JY, *et al.* 2019. Diabetologia. 62:993. [PubMed](#)
282. Chen L, *et al.* 2021. Theranostics. 11:6668. [PubMed](#)
283. Zhang A, *et al.* 2021. Theranostics. 11:3839. [PubMed](#)
284. Adzika GK, *et al.* 2021. Front Cell Dev Biol. 9:719351. [PubMed](#)
285. Liu SS, *et al.* 2021. Immunity. .: [PubMed](#)
286. Freed-Pastor WA, *et al.* 2021. Cancer Cell. .: [PubMed](#)
287. Cao H, *et al.* 2020. Front Pharmacol. 11:578943. [PubMed](#)
288. Ngkelo A, *et al.* 2016. J Exp Med. 213: 1353 - 1374. [PubMed](#)
289. Sockolosky J, *et al.* 2016. Proc Natl Acad Sci U S A. 113: 2646 - 2654. [PubMed](#)
290. Wan X, Thomas J, Unanue E 2016. J Exp Med. 213: 967 - 978. [PubMed](#)
291. Ishizuka M, *et al.* 2015. PLoS One. 10: 0137978. [PubMed](#)
292. Pechous R, *et al.* 2015. MBio. 6:2302. [PubMed](#)
293. Yin T, *et al.* 2015. J Immunol. 194:429. [PubMed](#)
294. Klezovich-Bénard M, *et al.* 2012. PLoS One. 8:e1002481. [PubMed](#)
295. Karpurapu M, *et al.* 2011. Blood. 118:5255. [PubMed](#)
296. Yan K, *et al.* 2020. Cell Death Dis. 1.013888889. [PubMed](#)
297. Eghbalzadeh K, *et al.* 2019. Front Immunol. 2.022916667. [PubMed](#)
298. Ouadah Y, *et al.* 2019. Cell. 179:403. [PubMed](#)
299. Ling X, *et al.* 2021. Biomaterials. 270:120690. [PubMed](#)
300. Lau P, *et al.* 2022. Cell Mol Immunol. .: [PubMed](#)
301. Zhou Y, *et al.* 2022. Theranostics. 12:5488. [PubMed](#)
302. Sheng D, *et al.* 2022. J Immunother Cancer. 10:. [PubMed](#)
303. Teng F, *et al.* 2021. Cell Rep. 37:110051. [PubMed](#)
304. Adzika GK, *et al.* 2021. Front Cardiovasc Med. 8:719805. [PubMed](#)
305. Wang X, *et al.* 2021. Theranostics. 11:209. [PubMed](#)
306. Collins C, *et al.* 2014. Proc Natl Acad Sci U S A. 111:9899. [PubMed](#)
307. Bramwell K, *et al.* 2015. J Immunol. 195: 1647 - 1656. [PubMed](#)
308. Legrain S, *et al.* 2015. J Immunol. 195: 4171 - 4175. [PubMed](#)
309. Schuijs M, *et al.* 2016. PLoS One. 11: 016188. [PubMed](#)
310. Lerbs T, *et al.* 2020. JCI Insight. 5:00. [PubMed](#)
311. Wang J, *et al.* 2020. Cell. 183(7):1867-1883.e26. [PubMed](#)
312. Wu L, *et al.* 2020. Cancer Immunol Res. 710:8. [PubMed](#)
313. Fu W, *et al.* 2020. J Immunol Res. 8820355:2020. [PubMed](#)
314. Saul S, *et al.* 2019. J Leukoc Biol. N/A. [PubMed](#)
315. Koupenova M, *et al.* 2019. iScience. 19:916. [PubMed](#)

316. Yuan Y, *et al.* 2019. Clin Sci (Lond). 133:1759. [PubMed](#)
 317. Hong Y, *et al.* 2019. J Extracell Vesicles. 8:1670893. [PubMed](#)
 318. Boguslawski KM, *et al.* 2020. Sci Adv. 6:eaax7515. [PubMed](#)
 319. Kim M, *et al.* 2012. J Immunol. 188:4158. [PubMed](#)
 320. Gardner A, *et al.* 2022. J Immunother Cancer. 10:. [PubMed](#)
 321. Tang L, *et al.* 2022. Front Immunol. 12:770402. [PubMed](#)
 322. Caporarello N, *et al.* 2022. Nat Commun. 13:4170. [PubMed](#)
 323. Zhang B, *et al.* 2021. Front Neurol. 607370:12. [PubMed](#)
 324. Kapralov AA, *et al.* 2020. Nat Chem Biol. 278:16. [PubMed](#)
 325. Xiang W, *et al.* 2020. Signal Transduct Target Ther. 0.374305556. [PubMed](#)
 326. Badgular DC, *et al.* 2020. PLoS Pathog. 16:e1009016. [PubMed](#)
 327. Steinmetz O, *et al.* 2016. J Infect Dis. 10.1093/infdis/jiv763. [PubMed](#)
 328. Malik M, *et al.* 2012. Proc Natl Acad Sci U S A. 109:6130. [PubMed](#)
 329. Kim SH, *et al.* 2020. Neoplasia. 1.3375. [PubMed](#)
 330. Maas SLN, *et al.* 2020. J Neuroinflammation. 17:120. [PubMed](#)
 331. Thomson CA, *et al.* 2018. J Immunol. 201:215. [PubMed](#)
 332. Sparber F *et al.* 2019. British journal of pharmacology. 176(11):1728-1744 . [PubMed](#)
 333. Lalor R, *et al.* 2019. Nutrients. 11. [PubMed](#)
 334. Fearon AE, *et al.* 2021. iScience. 24:103143. [PubMed](#)
 335. Chung EJ, *et al.* 2021. Int J Radiat Oncol Biol Phys. 110:526. [PubMed](#)
 336. Murray CW, *et al.* 2022. Nat Commun. 13:1090. [PubMed](#)
 337. Qu J, *et al.* 2021. Int J Biol Sci. 17:2756. [PubMed](#)
 338. Liu ZZ, *et al.* 2021. Front Immunol. 12:636861. [PubMed](#)

RRID AB_893493 (BioLegend Cat. No. 123115)
 AB_893481 (BioLegend Cat. No. 123116)

Antigen Details

Structure	EGF-TM7 family member, 160 kD glycoprotein
Distribution	Majority of tissue macrophages including peritoneal macrophages, macrophages in lung, gut, thymus and red pulp of spleen, Kuffer cells, Langerhans cells, bone marrow stromal cells, and a subset of dendritic cells
Function	Induction of immunological tolerance
Cell Type	Dendritic cells, Langerhans cells, Macrophages, Tregs
Biology Area	Cell Biology, Immunology, Innate Immunity, Neuroinflammation, Neuroscience
Antigen References	<ol style="list-style-type: none"> 1. Austy JM and Gordon S. 1981. <i>Eur. J. Immunol.</i> 11:805. 2. Hume DA, <i>et al.</i> 1983. <i>J. Exp. Med.</i> 158:1522. 3. Ruedl C, <i>et al.</i> 1996. <i>Eur. J. Immunol.</i> 26:1801. 4. McKnight AJ, <i>et al.</i> 1996. <i>J. Biol. Chem.</i> 271:486. 5. Lin HH, <i>et al.</i> 2005. <i>J. Exp. Med.</i> 201:1615.
Gene ID	13733

Related Protocols


[Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

Brilliant Violet 605™ anti-mouse F4/80, Purified anti-mouse F4/80, Biotin anti-mouse F4/80, FITC anti-mouse F4/80, PE anti-mouse F4/80, PE/Cyanine5 anti-mouse F4/80, PE/Cyanine7 anti-mouse F4/80, APC anti-mouse F4/80, APC/Cyanine7 anti-mouse F4/80, Alexa Fluor® 488 anti-mouse F4/80, Alexa Fluor® 647 anti-mouse F4/80, Pacific Blue™ anti-mouse F4/80, PerCP anti-mouse F4/80, PerCP/Cyanine5.5 anti-mouse F4/80, Alexa Fluor® 700 anti-mouse F4/80, Brilliant Violet 421™ anti-mouse F4/80, Brilliant Violet 510™ anti-mouse F4/80, Alexa Fluor® 594 anti-mouse F4/80, Brilliant Violet 785™ anti-mouse F4/80, Purified anti-mouse F4/80 (Maxpar® Ready), PE/Dazzle™ 594 anti-mouse F4/80, Brilliant Violet 650™ anti-mouse F4/80, Brilliant Violet 711™ anti-mouse F4/80, APC/Fire™ 750 anti-mouse F4/80, TotalSeq™-A0114 anti-mouse F4/80, TotalSeq™-B0114 anti-mouse F4/80, TotalSeq™-C0114 anti-mouse F4/80, Spark YG™ 570 anti-mouse F4/80, KIRAVIA Blue 520™ anti-mouse F4/80, Ultra-LEAF™ Purified anti-mouse F4/80, APC/Fire™ 810 anti-mouse F4/80, Spark NIR™ 685 anti-mouse F4/80

Product Data

Thioglycolate-elicited BALB/c mouse
peritoneal macrophages stained with
BM8 APC

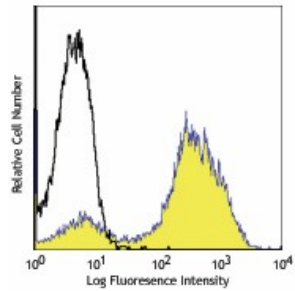

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587