

Brilliant Violet 605™ anti-human CD83 Antibody

| | |
|--------------------------|---|
| Catalog# / Size | 305335 / 25 tests 305336 / 100 tests |
| Clone | HB15e |
| Regulatory Status | RUO |
| Other Names | HB15 |
| Isotype | Mouse IgG1, κ |
| Description | CD83 is a 43 kD single chain type I glycoprotein also known as HB15. A member of the immunoglobulin superfamily, CD83 is expressed on a subset of dendritic cells, Langerhans cells, and weakly on activated lymphocytes. Although CD83 is thought to play a role in antigen presentation and/or lymphocyte activation, the precise function of this protein is unknown. CD83 is considered to be a useful marker for mature dendritic cells. |

Product Details

| | |
|-------------------------------|---|
| Verified Reactivity | Human |
| Reported Reactivity | African Green, Baboon, Pigtailed Macaque, Rhesus |
| Antibody Type | Monoclonal |
| Host Species | Mouse |
| Formulation | Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA). |
| Preparation | The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 605™ under optimal conditions. |
| Concentration | Lot-specific (to obtain lot-specific concentration, please enter the lot number in our Concentration and Expiration Lookup or Certificate of Analysis online tools.) |
| Storage & Handling | The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze. |
| Application | FC - Quality tested |
| Recommended Usage | <p>Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 µl per million cells in 100 µl staining volume or 5 µl per 100 µl of whole blood.</p> <p>Brilliant Violet 605™ excites at 405 nm and emits at 603 nm. The bandpass filter 610/20 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel. Refer to your instrument manual or manufacturer for support. Brilliant Violet 605™ is a trademark of Sirigen Group Ltd.</p> <p>Learn more about Brilliant Violet™.</p> <p>This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.</p> |
| Excitation Laser | Violet Laser (405 nm) |
| Application Notes | Additional reported applications (for the relevant formats) include: immunohistochemical staining of acetone-fixed frozen tissue sections ⁴ . |
| Application References | 1. Knapp W, <i>et al.</i> 1989. Leucocyte Typing IV. Oxford University Press New York. |

(PubMed link indicates BioLegend citation)

2. Zhou L, *et al.* 1995. *J. Immunol.* 154:3821.
3. Cao W, *et al.* 2005. *Biochem. J.* 385:85.
4. Lore K, *et al.* 2002. *AIDS* 16:683. (IHC)
5. Cho H, *et al.* 2007. *Physiol Genomics* doi:10.1152/physiolgenomics.00051.2006

RRID AB_2687382 (BioLegend Cat. No. 305335)
AB_2687383 (BioLegend Cat. No. 305336)

Antigen Details

| | |
|---------------------------|---|
| Structure | Ig superfamily, single chain transmembrane glycoprotein, 43 kD |
| Distribution | Dendritic cells, Langerhan cells, activated B and T cells |
| Cell Type | B cells, Dendritic cells, Langerhans cells, T cells |
| Biology Area | Costimulatory Molecules, Immunology |
| Molecular Family | CD Molecules |
| Antigen References | 1. Kozlow E, <i>et al.</i> 1993. <i>Blood</i> 81:454. 2. Zhou L, <i>et al.</i> 1992. <i>J. Immunol.</i> 149:735. 3. Zhou L, <i>et al.</i> 1995. <i>Blood</i> 86:3295. |
| Gene ID | 9308 |


Related Protocols

[Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

APC anti-human CD83, Biotin anti-human CD83, FITC anti-human CD83, PE anti-human CD83, PE/Cyanine5 anti-human CD83, Purified anti-human CD83, Alexa Fluor® 488 anti-human CD83, Alexa Fluor® 647 anti-human CD83, PerCP/Cyanine5.5 anti-human CD83, Brilliant Violet 421™ anti-human CD83, PE/Cyanine7 anti-human CD83, PE/Dazzle™ 594 anti-human CD83, APC/Cyanine7 anti-human CD83, Brilliant Violet 711™ anti-human CD83, APC/Fire™ 750 anti-human CD83, Brilliant Violet 605™ anti-human CD83, Brilliant Violet 785™ anti-human CD83, TotalSeq™-A0359 anti-human CD83, TotalSeq™-C0359 anti-human CD83, TotalSeq™-B0359 anti-human CD83, TotalSeq™-D0359 anti-human CD83

Product Data

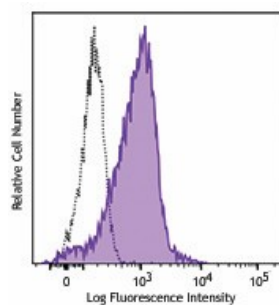

Human monocyte-derived dendritic cells (induced with GM-CSF and IL-4 and TNF- α) stained with CD83 Brilliant Violet 605™ (Clone HB15e, filled histogram) or Mouse IgG1, κ Brilliant Violet 605™ isotype control (open histogram).

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587