

## Brilliant Violet 650™ anti-human CD137 (4-1BB) Antibody

<b>Catalog# / Size</b>	309827 / 25 tests 309828 / 100 tests
<b>Clone</b>	4B4-1
<b>Regulatory Status</b>	RUO
<b>Workshop</b>	VI C-7
<b>Other Names</b>	4-1BB, ILA, CD137, TNFRSF9
<b>Isotype</b>	Mouse IgG1, κ
<b>Description</b>	CD137 is a 39 kD transmembrane protein also known as 4-1BB. It is expressed on activated T cells. CD137 is a type I membrane protein and a member of the tumor necrosis factor receptor superfamily. CD137 appears to be important for T cell proliferation and survival, and induces monocyte activation through its interaction with 4-1BB ligand.

### Product Details

---

<b>Verified Reactivity</b>	Human
<b>Reported Reactivity</b>	Chimpanzee, Baboon, Cynomolgus, Rhesus
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Mouse
<b>Immunogen</b>	Ectodomain of recombinant human 4-1BB fusion protein
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA).
<b>Preparation</b>	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 650™ under optimal conditions.
<b>Concentration</b>	Lot-specific (to obtain lot-specific concentration, please enter the lot number in our <a href="#">Concentration and Expiration Lookup</a> or <a href="#">Certificate of Analysis</a> online tools.)
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	<p>Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a>. For flow cytometric staining, the suggested use of this reagent is 5 µl per million cells in 100 µl staining volume or 5 µl per 100 µl of whole blood.</p> <p>Brilliant Violet 650™ excites at 405 nm and emits at 645 nm. The bandpass filter 660/20 nm is recommended for detection, although filter optimization may be required depending on other fluorophores used. <b>Be sure to verify that your cytometer configuration and software setup are appropriate for detecting this channel.</b> Refer to your instrument manual or manufacturer for support. Brilliant Violet 650™ is a trademark of Sirigen Group Ltd.</p> <p><a href="#">Learn more about Brilliant Violet™.</a></p> <p>This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.</p>
<b>Excitation Laser</b>	Violet Laser (405 nm)
<b>Application Notes</b>	Additional reported applications (for the relevant formats) include: immunoprecipitation <sup>1,4</sup> , inhibition of cytokine production <sup>2,3</sup> , and ELISA. For most successful immunofluorescent staining results, it

may be important to maximize signal over background by using a relatively bright fluorochrome-antibody conjugate (Cat. No. 309804) or by using a high sensitivity, three-layer staining technique (e.g., including a biotinylated anti-mouse IgG second step (Cat. No. 405303), followed by Streptavidin-PE (Cat. No. 405204)).

#### Application References

(PubMed link indicates BioLegend citation)

1. Garni-Wagner B, *et al.* 1996. *Cell. Immunol.* 169:91. (IP)
2. Salih HR, *et al.* 2000. *J. Immunol.* 165:2903. (FA)
3. Kienzle G, *et al.* 2000. *Int. Immunol.* 12:73. (FA)
4. Langstein J, *et al.* 1998. *J. Immunol.* 160:2488. (IP)

#### Product Citations

1. Nielsen CM, *et al.* 2021. *Cell Reports Medicine.* 2(3):100207. [PubMed](#)
2. Minassian AM, *et al.* 2021. *Med (N Y).* 2:701. [PubMed](#)

#### RRID

AB\_2572192 (BioLegend Cat. No. 309827)  
AB\_2572193 (BioLegend Cat. No. 309828)

## Antigen Details

---

**Structure** TNFR superfamily, type I transmembrane protein, 30 kD

**Distribution** Activated T cells

**Function** T cell costimulation

**Ligand/Receptor** 4-1BB ligand

**Cell Type** T cells

**Biology Area** Costimulatory Molecules, Immunology

**Molecular Family** CD Molecules

#### Antigen References

1. Gruss H, *et al.* 1995. *Blood* 85:3378.
2. Sica G, *et al.* 2000. *Adv. Exp. Med. Biol.* 465:355.
3. Alderson M, *et al.* 1994. *Eur. J. Immunol.* 24:2219.
4. Schwarz H, *et al.* 1996. *Blood* 87:2839.

#### Gene ID

[3604](#)

## Related Protocols

---

[Cell Surface Flow Cytometry Staining Protocol](#)


## Other Formats

---

Purified anti-human CD137 (4-1BB), PE anti-human CD137 (4-1BB), Biotin anti-human CD137 (4-1BB), PE/Cyanine5 anti-human CD137 (4-1BB), APC anti-human CD137 (4-1BB), PerCP/Cyanine5.5 anti-human CD137 (4-1BB), Alexa Fluor® 700 anti-human CD137 (4-1BB), PE/Cyanine7 anti-human CD137 (4-1BB), Brilliant Violet 421™ anti-human CD137 (4-1BB), APC/Cyanine7 anti-human CD137 (4-1BB), Brilliant Violet 605™ anti-human CD137 (4-1BB), Alexa Fluor® 647 anti-human CD137 (4-1BB), PE/Dazzle™ 594 anti-human CD137 (4-1BB), Brilliant Violet 650™ anti-human CD137 (4-1BB), Brilliant Violet 711™ anti-human CD137 (4-1BB), APC/Fire™ 750 anti-human CD137 (4-1BB), TotalSeq™-A0355 anti-human CD137 (4-1BB), TotalSeq™-B0355 anti-human CD137 (4-1BB), TotalSeq™-C0355 anti-human CD137 (4-1BB), Ultra-LEAF™ Purified anti-human CD137 (4-1BB), Brilliant Violet 750™ anti-human CD137 (4-1BB), TotalSeq™-D0355 anti-human CD137 (4-1BB)

## Product Data

---

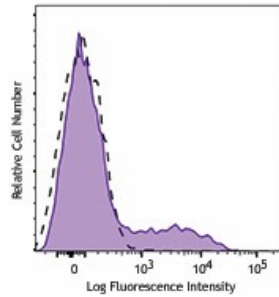

PHA-stimulated (3 days) human peripheral blood lymphocytes were stained with CD137 (clone 4B4-1) Brilliant Violet 650™ (filled histogram) or mouse IgG1, κ Brilliant Violet 650™ isotype control (open histogram).

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

\*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, [www.biolegend.com/ordering#license](http://www.biolegend.com/ordering#license)). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587