

DAPI (4',6-Diamidino-2-Phenylindole, Dilactate)

Catalog# / Size 422801 / 10 mg

Regulatory Status RUO

Other Names 4',6-Diamidino-2-Phenylindole, Dilactate

Description DAPI (4',6-Diamidino-2-Phenylindole, Dilactate), a blue fluorescent nucleic acid stain, preferentially binds dsDNA and associates with the minor groove AT clusters. The molecular weight for DAPI is 457.5. Association of DAPI with dsDNA increases fluorescence approximately 20-fold, with an emission maximum of 460 nm. The increase in fluorescence is believed to be due to displacement of water molecules from the minor groove and DAPI. DAPI can also bind RNA and emits at a longer wavelength of 500 nm, with only a 20% increase in fluorescence.

DAPI is a popular nuclear stain for identification of cell cycle and specifically stains nuclei but not cytoplasm. The dilactate form is more water soluble compared to the dihydrochloride form.

DAPI can be excited by a UV laser, xenon, or a mercury-arc lamp.

Product Details

Formulation Lyophilized

Preparation To make 5 mg/ml (10.9 mM) solution, dissolve entire contents of vial in 2 ml of deionized water.

Storage & Handling Store between 2°C and 8°C, protected from light. For long term storage, aliquot and store at -20°C.

Application [ICFC](#), [ICC](#)

Recommended Usage This product can be used at 3 µM for flow cytometry and 300 nM for immunofluorescence. It is recommended to titrate the dye for optimal performance.

Excitation Laser Ultraviolet Laser (355 nm)

Application Notes **Protocols for Counterstaining with DAPI:**

Immunofluorescence Staining

Note: Cells may be fixed by method of choice, such as with 4% paraformaldehyde. DAPI staining is done after staining for other markers. Fixation and permeabilization of cells is not necessary to counterstain with DAPI.

1. Fix cells using method of choice.
2. Incubate the cells with phosphate buffered saline (PBS) for 15 minutes.
3. Dilute the DAPI solution (we recommend to 300 nM) with PBS. Remove excess PBS from the slide and cover with DAPI solution, making sure the cells are completely covered.
4. Incubate for up to 5 minutes. Rinse the slide several times to remove all free DAPI. It is recommended to use a mounting medium with antifade reagent to reduce fluorescence quenching. The slide is ready to view under the microscope with appropriate filters.

Nuclear Staining for Flow Cytometry

Note: Cells may be fixed with 4% paraformaldehyde or absolute ethanol.

1. Fix cells with absolute ethanol or 4% paraformaldehyde.
2. Centrifuge, discard the supernatant, and add 5 ml of PBS, allowing cells to rehydrate for 15 minutes.
3. Dilute the DAPI to 3 µM (recommended) in staining buffer (such as BioLegend's FOXP3 Perm Buffer, Cat. No. 421402).
4. Centrifuge cell suspension from step 3 above, discard supernatant, and add 1 ml of DAPI solution; then incubate at room temperature for 15 minutes.
5. Analyze by flow cytometry using an instrument with appropriate lasers.

Additional Product Notes View more applications data for this product in our [Scientific Poster Library](#).

Application References

1. Solomonov I, *et al.* 2014. *Sci Rep.* 4:5987. [PubMed](#)
2. Sodhi CP, *et al.* 2015. *J Immunol.* 194:4931. [PubMed](#)

(PubMed link indicates BioLegend citation)

Product Citations

1. Sierra F, *et al.* 2017. *Immunity.* 47:374. [PubMed](#)
2. Davila ML, *et al.* 2022. *iScience.* 25:104426. [PubMed](#)
3. Grubwieser P, *et al.* 2022. *Front Cell Infect Microbiol.* 12:875543. [PubMed](#)
4. Sountoulidis A, *et al.* 2020. *PLoS Biol.* 18:e3000675. [PubMed](#)
5. Mester B, *et al.* 2015. *PLoS One.* 10: e0140432. [PubMed](#)
6. Wu N, *et al.* 2020. *Cell Rep.* 30:1129. [PubMed](#)
7. Velilla J, *et al.* 2019. *Neurol Genet.* 5:e312. [PubMed](#)
8. Jiao Q, *et al.* 2017. *Clin Exp Immunol.* 188:420. [PubMed](#)
9. Rodrigues M, *et al.* 2022. *Biochem J.* 479:1891. [PubMed](#)
10. Gabriely G, *et al.* 2021. *iScience.* 24:103347. [PubMed](#)
11. Effern M, *et al.* 2022. *STAR Protoc.* 3:101038. [PubMed](#)
12. Wilson K, *et al.* 2016. *Cell Death Dis.* 7: e2197. [PubMed](#)
13. Kobayashi T, *et al.* 2019. *Cell.* 176:982. [PubMed](#)
14. Anguiano T, *et al.* 2020. *Toxicol Sci.* 177:494. [PubMed](#)
15. Afshar Y, *et al.* 2021. *Nat Commun.* 12:4408. [PubMed](#)
16. Goodnough LH, *et al.* 2022. *Front Endocrinol (Lausanne).* 13:924927. [PubMed](#)
17. Nguyen DT, *et al.* 2022. *Cells.* 11: . [PubMed](#)
18. Baldominos P, *et al.* 2022. *STAR Protoc.* 3:101795. [PubMed](#)
19. Gamboa CM, *et al.* 2021. *STAR Protoc.* 100345:2. [PubMed](#)
20. Sugawara E, *et al.* 2020. *Autophagy.* 1.323611111. [PubMed](#)
21. Elias A, *et al.* 2016. *Sci Rep.* 6: 24971. [PubMed](#)
22. Miedema A, *et al.* 2022. *Acta Neuropathol Commun.* 10:8. [PubMed](#)
23. Emmons TR, *et al.* 2021. *Cancer Immunol Res.* 9:790. [PubMed](#)
24. Valkov N, *et al.* 2021. *Life Sci Alliance.* 4: . [PubMed](#)
25. Bannard O, *et al.* 2016. *J Exp Med.* 13: 993 - 1009. [PubMed](#)
26. Sumitomo Y, *et al.* 2016. *Blood.* 128: 1614 - 1624. [PubMed](#)
27. Paris AJ, *et al.* 2020. *Nat Cell Biol.* 1.747916667. [PubMed](#)
28. Brunner JS, *et al.* 2020. *Nat Commun.* 0.757638889. [PubMed](#)
29. Fukushima T, *et al.* 2019. *Cell Rep.* 29:4144. [PubMed](#)
30. Hopp AK, *et al.* 2021. *Mol Cell.* 81:340. [PubMed](#)
31. Zhang M, *et al.* 2017. *Oncol Lett.* 3.069444444. [PubMed](#)
32. Greenblatt SM, *et al.* 2018. *Cancer Cell.* 33:1111. [PubMed](#)
33. Alexander MS, *et al.* 2021. *Muscle Nerve.* 63:928. [PubMed](#)
34. Hendricks A, *et al.* 2021. *Front Oncol.* 11:646885. [PubMed](#)
35. Hiramatsu-Asano S, *et al.* 2021. *Front Immunol.* 11:616141. [PubMed](#)
36. Ambrosi TH, *et al.* 2021. *Elife.* 10: . [PubMed](#)
37. DeDreu J, *et al.* 2020. *FASEB J.* 34:9316. [PubMed](#)
38. Allam O, *et al.* 2018. *PLoS One.* 13:e0194185. [PubMed](#)
39. Wilson AS, *et al.* 2022. *Nat Commun.* 13:528. [PubMed](#)
40. Deng J, *et al.* 2021. *Breast Cancer Res Treat.* 185:85. [PubMed](#)
41. Koh P, *et al.* 2016. *Sci Data.* 3:160109. [PubMed](#)
42. Sodhi C, *et al.* 2015. *J Immunol.* 194:4931. [PubMed](#)
43. Kanda K, *et al.* 2018. *JCI Insight.* 3:e91316. [PubMed](#)
44. Fujiwara H, *et al.* 2021. *Nat Immunol.* 22:1440. [PubMed](#)
45. Ravi A, *et al.* 2021. *Dev Cell.* 56:1661. [PubMed](#)
46. Wu CL, *et al.* 2021. *Nat Commun.* 12:362. [PubMed](#)
47. Poulos M, *et al.* 2016. *Nat Commun.* 7:13829. [PubMed](#)
48. Noah AC, *et al.* 2020. *J Appl Physiol (1985).* 473:128. [PubMed](#)
49. Wang ECE, *et al.* 2019. *Cell Stem Cell.* 24:654. [PubMed](#)
50. So EY, *et al.* 2021. *Am J Physiol Cell Physiol.* 321:C569. [PubMed](#)
51. Hendricks A, *et al.* 2020. *Cancers (Basel).* 12:00. [PubMed](#)
52. Aronovich A, *et al.* 2021. *J Invest Dermatol.* 141:619. [PubMed](#)
53. Teo YV, *et al.* 2019. *Cell Rep.* 27:997. [PubMed](#)
54. Hayama Y, *et al.* 2018. *J Immunol.* 200:3790. [PubMed](#)
55. Fredriksson-Lidman K, *et al.* 2017. *PLoS One.* 10.1371/journal.pone.0185509. [PubMed](#)
56. Wennerberg E, *et al.* 2022. *Sci Transl Med.* 14:eabe8195. [PubMed](#)
57. Kumar A, *et al.* 2021. *STAR Protoc.* 2:100874. [PubMed](#)
58. Faisal S, *et al.* 2016. *Sci Rep.* 6:39530. [PubMed](#)
59. Caeser R, *et al.* 2021. *Nat Protoc.* 16:2499. [PubMed](#)
60. Thompson M, *et al.* 2015. *J Immunol.* 195: 5637 - 5647. [PubMed](#)
61. Kundu K, *et al.* 2019. *Cancer Immunol Res.* 1.069444444. [PubMed](#)
62. Pagan JD, *et al.* 2018. *Cell.* 172:564. [PubMed](#)
63. Zhang X, *et al.* 2022. *JCI Insight.* 7: . [PubMed](#)
64. Zhu M, *et al.* 2021. *Eneuro.* 8: . [PubMed](#)
65. Flores-Pliego A, *et al.* 2021. *Cells.* 10: . [PubMed](#)
66. Shani O, *et al.* 2021. *Elife.* 10: . [PubMed](#)
67. Moser J, *et al.* 2016. *Sci Rep.* 6:35387. [PubMed](#)
68. Mould KJ, *et al.* 2019. *JCI Insight.* 4:00. [PubMed](#)
69. Heath O, *et al.* 2021. *Cancer Immunol Res.* 9:665. [PubMed](#)
70. Dasinger JH, *et al.* 2021. *Pregnancy Hypertens.* 24:126. [PubMed](#)
71. Ren X, *et al.* 2022. *STAR Protoc.* 3:101818. [PubMed](#)
72. Correa A, *et al.* 2015. *PLoS One.* 10: 0143391. [PubMed](#)
73. Yanting C, *et al.* 2018. *J Dairy Sci.* 101:1601. [PubMed](#)
74. Kurena B, *et al.* 2017. *Front Immunol.* 1.490972222. [PubMed](#)

75. Pinho-Ribeiro FA, *et al.* 2018. *Cell*. 173:1083. [PubMed](#)
76. Aziz M, *et al.* 2017. *J Immunol*. 10.4049/jimmunol.1602056. [PubMed](#)
77. Pascual-Pasto G, *et al.* 2022. *J Immunother Cancer*. 10:. [PubMed](#)
78. Lehrke MJ, *et al.* 2021. *Elife*. 10:. [PubMed](#)
79. Wang Y, *et al.* 2021. *EMBO Mol Med*. 13:e14089. [PubMed](#)
80. Chiang M, *et al.* 2019. *Cell Rep*. 28:3212. [PubMed](#)
81. Valdespino-Vázquez MY, *et al.* 2021. *J Med Virol*. 93:4480. [PubMed](#)
82. Kissiov DU, *et al.* 2022. *Elife*. 11:. [PubMed](#)
83. Isobe T, *et al.* 2022. *Nat Commun*. 13:4501. [PubMed](#)
84. Caporarello N, *et al.* 2020. *Aging Cell*. 19:e13196. [PubMed](#)
85. Vergoni B, *et al.* 2016. *Diabetes*. 65: 3062 - 3074. [PubMed](#)
86. Dent JR, *et al.* 2019. *J Appl Physiol (1985)*. 127:1117. [PubMed](#)
87. Schmidt A, *et al.* 2017. *Toxicol Lett*. 10.1016/j.toxlet.2017.08.008. [PubMed](#)
88. Canals Hernaez D, *et al.* 2022. *Front Oncol*. 12:856424. [PubMed](#)
89. Mascharak S, *et al.* 2022. *Cell Stem Cell*. 29:315. [PubMed](#)
90. Rodriguez A, *et al.* 2021. *Cell Stem Cell*. 28:33. [PubMed](#)
91. Pang J, *et al.* 2021. *STAR Protoc*. 2:100879. [PubMed](#)
92. Haschak M, *et al.* 2021. *Aging (Albany NY)*. 13:16938. [PubMed](#)
93. Phan TS, *et al.* 2021. *Sci Adv*. 7:. [PubMed](#)
94. Stancil IT, *et al.* 2021. *Nat Commun*. 12:4566. [PubMed](#)
95. Margaroli C, *et al.* 2021. *Cell Rep Med*. :100242. [PubMed](#)
96. Burns F, *et al.* 2015. *Aquatic Toxicology*. 164: 52-60. [PubMed](#)
97. Stevens M, *et al.* 2016. *Nat Biotechnol*. 34:1161-1167. [PubMed](#)
98. Shiels J, *et al.* 2021. *Sci Rep*. 0.745138889. [PubMed](#)
99. Kohn M, *et al.* 2021. *Front Immunol*. 626627:12. [PubMed](#)
100. Jain A, *et al.* 2020. *Nat Immunol*. 0.920138889. [PubMed](#)
101. Zheng C, *et al.* 2022. *Acta Pharm Sin B*. 12:3726. [PubMed](#)
102. Anisul M, *et al.* 2021. *Elife*. 10:. [PubMed](#)
103. Pelham-Webb B, *et al.* 2021. *Mol Cell*. 81:1732. [PubMed](#)
104. Chmielewski M and Abken H 2017. *Cell Rep*. 10.1016/j.celrep.2017.11.063. [PubMed](#)
105. Hassanein W, *et al.* 2017. *Organogenesis*. 13:16. [PubMed](#)
106. Aimi F, *et al.* 2021. *Mod Pathol*. 34:1468. [PubMed](#)
107. Barcia Durán JG, *et al.* 2021. *Commun Biol*. 4:406. [PubMed](#)
108. Karimzadeh A, *et al.* 2021. *Front Cell Dev Biol*. 9:734176. [PubMed](#)
109. Reeves KM, *et al.* 2021. *Clin Cancer Res*. Online ahead of print. [PubMed](#)
110. Chen HY, *et al.* 2020. *J Immunol*. 1158:204. [PubMed](#)
111. Gosztyla C, *et al.* 2020. *Tissue Eng Part C Methods*. 26:67. [PubMed](#)
112. Köhler VK, *et al.* 2020. *Int J Mol Sci*. 21:00. [PubMed](#)
113. Solomonov I, *et al.* 2014. *Sci Rep*. 4:5987. [PubMed](#)
114. Srijiwangsa P, *et al.* 2018. *BMC Pharmacol Toxicol*. 19:32. [PubMed](#)
115. Ladd MR, *et al.* 2018. *Tissue Eng Part C Methods*. 24:331. [PubMed](#)
116. Bennett AH, *et al.* 2018. *PLoS Genet*. 14:e1007226. [PubMed](#)
117. Doorley L, *et al.* 2017. *Immunobiology*. 10.1016/j.imbio.2017.08.002. [PubMed](#)
118. Guo P, *et al.* 2022. *Nat Cell Biol*. 24:99. [PubMed](#)
119. Tominaga K, *et al.* 2022. *iScience*. 25:103667. [PubMed](#)


Antigen Details

Biology Area Cell Biology, Cell Cycle/DNA Replication, Neuroscience


Antigen References 1. Langlois RG, *et al.* 1980. *Chromosoma*. 77:229.
2. Kubista M, *et al.* 1987. *Biochemistry* 26:4545.

Gene ID NA


Product Data

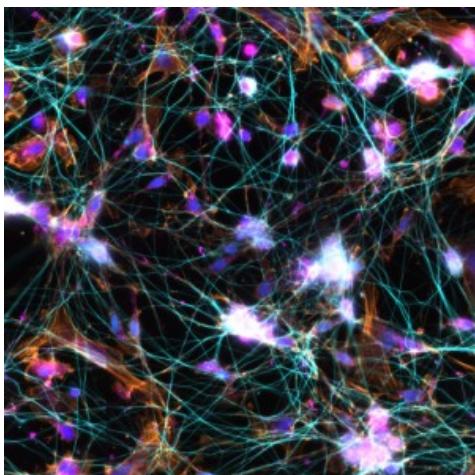


Methanol fixed superior cervical ganglion neurons were stained with anti-tubulin $\beta 3$ (pink) and anti-TRIM46 (aqua), then stained with secondary antibodies, phalloidin (orange) and DAPI (blue). Image generously submitted to the 2017 Cell Life Imaging Competition by Ankita Patil from Drexel University.

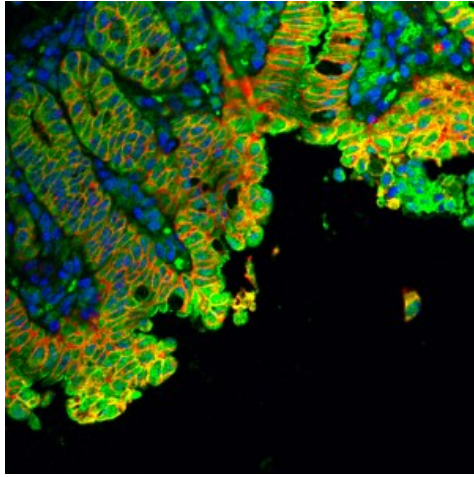


Frozen mouse intestine stained with purified anti-Ankyrin-B (green), E-Cadherin (red) and DAPI (blue) followed by anti-goat Alexa Fluor® 488 and anti-rabbit Alexa Fluor® 568. Image generously submitted to the 2017 Cell Life Imaging Competition by Kathleen Ignatoski from University of Michigan.

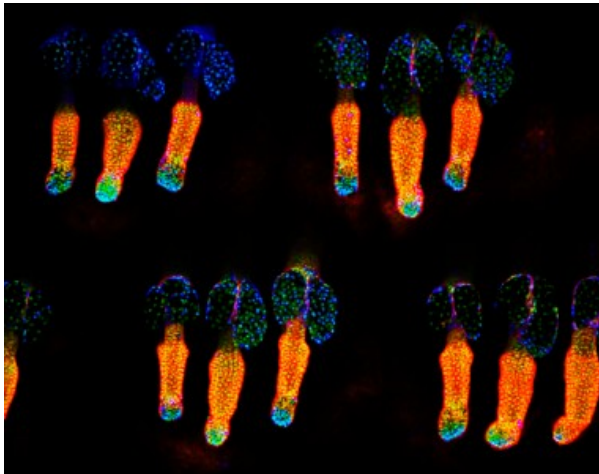

Whole mount mouse tail epidermis was stained with purified anti-Keratin 14 (green, clone Poly19053), purified anti-Keratin 15 (orange, clone Poly18339) and DAPI (blue). Image generously submitted to the 2017 Cell Life Imaging Competition by Kif Liakath.

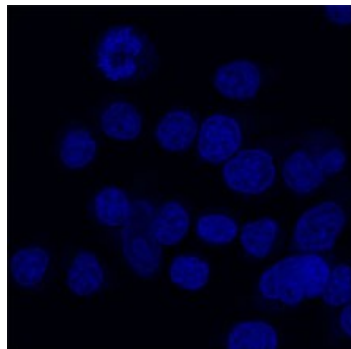

Human cervical cancer cell line, HeLa, was fixed with 4% paraformaldehyde and stained with 300 nM DAPI.

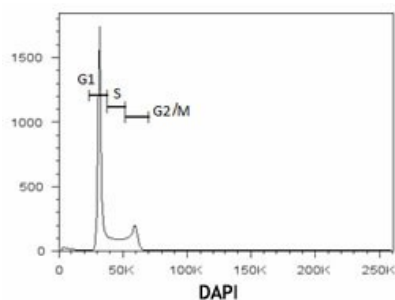

Human cervical cancer cell line, HeLa, was fixed with 4% paraformaldehyde and stained with 300 nM DAPI.

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587