

FITC anti-human CD16 Antibody

Catalog# / Size	302005 / 25 tests 302006 / 100 tests
Clone	3G8
Regulatory Status	RUO
Workshop	V NK80
Other Names	FcγRIII, Fc gamma receptor, Fc gamma receptor 3
Isotype	Mouse IgG1, κ
Description	CD16 is known as low affinity IgG receptor III (FcγRIII). It is expressed as two distinct forms (CD16a and CD16b). CD16a (FcγRIIIA) is a 50-65 kD polypeptide-anchored transmembrane protein. It is expressed on the surface of NK cells, activated monocytes, macrophages, and placental trophoblasts in humans. CD16b (FcγRIIIB) is a 48 kD glycosylphosphatidylinositol (GPI)-anchored protein. Its extracellular domain is over 95% homologous to that of CD16a, and it is expressed specifically on neutrophils. CD16 binds aggregated IgG or IgG-antigen complex which functions in NK cell activation, phagocytosis, and antibody-dependent cell-mediated cytotoxicity (ADCC).

Product Details

Verified Reactivity	Human, Cynomolgus, Rhesus
Reported Reactivity	African Green, Baboon, Capuchin Monkey, Chimpanzee, Common Marmoset, Pigtailed Macaque, Sooty Mangabey, Squirrel Monkey
Antibody Type	Monoclonal
Host Species	Mouse
Immunogen	Human PMN cells
Formulation	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA)
Preparation	The antibody was purified by affinity chromatography, and conjugated with FITC under optimal conditions.
Concentration	Lot-specific (to obtain lot-specific concentration, please enter the lot number in our Concentration and Expiration Lookup or Certificate of Analysis online tools.)
Storage & Handling	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze.
Application	FC - Quality tested
Recommended Usage	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis . For flow cytometric staining, the suggested use of this reagent is 5 μl per million cells in 100 μl staining volume or 5 μl per 100 μl of whole blood.
Excitation Laser	Blue Laser (488 nm)
Application Notes	<p>The 3G8 antibody clone blocks neutrophil phagocytosis and stimulates NK cell proliferation. It has been reported that this clone interacts with the FcγRIIIa and FcγRIIIb receptors causing neutrophil activation and aggregation¹⁸. Due to this phenomenon staining in whole blood may cause a reduction in the number of granulocytes or alter their scatter profile.</p> <p>Additional reported applications (for the relevant formats) include: immunohistochemical staining of acetone-fixed frozen tissue sections⁶, immunoprecipitation³, stimulation of NK cell proliferation⁴, blocking of phagocytosis⁵, and blocking of immunoglobulin binding to FcγRIII^{7,8}. The Ultra-LEAF™ purified antibody (Endotoxin < 0.01 EU/μg, Azide-Free, 0.2 μm filtered) is recommended for functional assays (Cat. No. 302049, 302050, 302057, 302058).</p>

Application References

**(PubMed link indicates
BioLegend citation)**

1. Knapp W, *et al.* Eds. 1989. Leucocyte Typing IV. Oxford University Press. New York.
2. Schlossman S, *et al.* Eds. 1995. Leucocyte Typing V. Oxford University Press. New York.
3. Edberg J, *et al.* 1997. *J. Immunol.* 159:3849. (IP)
4. Hoshino S, *et al.* 1991. *Blood* 78:3232. (Stim)
5. Tamm A, *et al.* 1996. *Immunol.* 157:1576. (Block)
6. Da Silva DM, *et al.* 2001. *Int. Immunol.* 13:633. (IHC)
7. Holl V, *et al.* 2004. *J. Immunol.* 173:6274. (Block)
8. Hober D, *et al.* 2002. *J. Gen. Virol.* 83:2169. (Block)
9. Brainard DM, *et al.* 2009. *J. Virol.* 83:7305. [PubMed](#)
10. Smed-Sørensen A, *et al.* 2008. *Blood* 111:5037. (Block) [PubMed](#)
11. Timmerman KL, *et al.* 2008. *J. Leukoc. Biol.* 84:1271. (FC) [PubMed](#)
12. Yoshino N, *et al.* 2000. *Exp. Anim. (Tokyo)* 49:97. (FC)
13. Rout N, *et al.* 2010. *PLoS One* 5:e9787. (FC)
14. Kim WK, *et al.* 2006. *Am. J. Pathol.* 168:822. (FC)
15. Boltz A, *et al.* 2011. *J. Biol Chem.* 286:21896. [PubMed](#)
16. Wu Z, *et al.* 2013. *J. Virol.* 87:7717. [PubMed](#)
17. Peterson VM, *et al.* 2017. *Nat. Biotechnol.* 35:936. (PG)
18. Vossebeld PJ, *et al.* 1997. *Biochem J.* 323:87-94 (Stim)

Product Citations

1. Babania O, *et al.* 2021. *J Clin Lab Anal.* 35:e23984. [PubMed](#)
2. Sagnella SM, *et al.* 2020. *Cancer Cell.* 37(3):354-370. [PubMed](#)
3. Okada T, *et al.* 2010. *Am J Pathol.* 176:2309. [PubMed](#)
4. Rosendahl Huber A, *et al.* 2022. *STAR Protoc.* 3:101361. [PubMed](#)
5. Brandsma AM, *et al.* 2021. *Blood Cancer Discov.* 2:484. [PubMed](#)
6. Tripathi H, *et al.* 2020. *J Mol Cell Cardiol.* 149:95. [PubMed](#)
7. Tong M, *et al.* 2019. *J Immunol.* 203:500. [PubMed](#)
8. Poonia B, Pauza C 2014. *PLoS One.* 9:88884. [PubMed](#)
9. Kerr SC, *et al.* 2020. *Clin Exp Allergy.* 50:904. [PubMed](#)
10. Laustsen A, *et al.* 2021. *Elife.* 10:. [PubMed](#)
11. Hasaart KAL, *et al.* 2020. *Sci Rep.* 10:12991. [PubMed](#)
12. AC Belkina, JE Snyder-Cappione 2017. *Cytometry A.* 91:175-179. [PubMed](#)
13. Søndergaard J, *et al.* 2014. *Int Immunol.* 26) : 705-16. [PubMed](#)
14. Seamon K, *et al.* 2020. *Front Immunol.* 11:396. [PubMed](#)
15. Evren E, *et al.* 2020. *Immunity.* 54(2):259-275.e7. [PubMed](#)
16. Xu Y, *et al.* 2021. *Stem Cell Res Ther.* 12:112. [PubMed](#)
17. Wu HL, *et al.* 2018. *J Immunol.* 200:49. [PubMed](#)
18. Wang Y, *et al.* 2019. *J Transl Med.* 17:93. [PubMed](#)
19. Gärtner K, *et al.* 2018. *Cancer Med.* 11:39. [PubMed](#)
20. Kobayashi Y, *et al.* 2020. *Int J Oncol.* 999:56. [PubMed](#)
21. Hunn M, *et al.* 2012. *Clin Cancer Res.* 18:6446. [PubMed](#)
22. Tong M, *et al.* 2021. *J Immunol.* 206:1039. [PubMed](#)
23. Brandsma AM, *et al.* 2021. *Cancer Discov.* . [PubMed](#)
24. Wang H, *et al.* 2022. *J Neuroinflammation.* 19:171. [PubMed](#)
25. Herter JM, *et al.* 2022. *Strahlenther Onkol.* Online ahead of print. [PubMed](#)
26. Zhang D, *et al.* 2020. *Oncoimmunology.* 9:1744921. [PubMed](#)
27. Yokoyama K, *et al.* 2022. *Front Immunol.* 13:1016914. [PubMed](#)
28. Evans RDR, *et al.* 2020. *Nat Commun.* 3.491666667. [PubMed](#)
29. Mulfaul K, *et al.* 2020. *Cell Rep.* 30:2209. [PubMed](#)
30. Ugawa M, *et al.* 2021. *Elife.* 10:. [PubMed](#)
31. Osorio FG *et al.* 2018. *Cell reports.* 25(9):2308-2316 . [PubMed](#)
32. Souriant S, *et al.* 2019. *Cell Rep.* 26:3586. [PubMed](#)
33. Walter F, *et al.* 2020. *PLoS One.* 15:e0239369. [PubMed](#)
34. Fisch D, *et al.* 2022. *Pathog Dis.* 79:. [PubMed](#)
35. Guo M, *et al.* 2021. *Front Oncol.* 11:632540. [PubMed](#)
36. Witkowski MT, *et al.* 2020. *Cancer Cell.* 37:867. [PubMed](#)
37. Zee BM, *et al.* 2021. *iScience.* 24(6):102651. [PubMed](#)
38. Sefik E, *et al.* 2021. *Nat Biotechnol.* . [PubMed](#)
39. Kacherovsky N, *et al.* 2019. *Nat Biomed Eng.* 0.66875. [PubMed](#)
40. Liu J, *et al.* 2022. *iScience.* 25:105509. [PubMed](#)
41. Hazenberg MD, *et al.* 2019. *Blood Adv.* 2.659722222. [PubMed](#)
42. Cheng Y, *et al.* 2021. *Immunity.* 54(8):1825-1840.e7. [PubMed](#)
43. Tong M, *et al.* 2021. *J Reprod Immunol.* 145:103306. [PubMed](#)
44. Yoshida R, *et al.* 2022. *Cancer Res.* .: [PubMed](#)
45. Ren C, *et al.* 2020. *J Immunol.* 204:1012. [PubMed](#)
46. Nowak W, *et al.* 2020. *EBioMedicine.* 50:290-305.. [PubMed](#)
47. Zhao NQ, *et al.* 2020. *PLoS One.* 15:e0238347. [PubMed](#)
48. Pahl JHW, *et al.* 2018. *Cancer Immunol Res.* 0.609027778. [PubMed](#)
49. Xue D, *et al.* 2013. *Transplant Proc.* 45:3279. [PubMed](#)
50. Hasaart KAL, *et al.* 2022. *iScience.* 25:103736. [PubMed](#)
51. Theurich S *et al.* 2017. *Cell metabolism.* 26(1):171-184 . [PubMed](#)
52. Bal S, *et al.* 2016. *Nat Immunol.* 10.1038/ni.3444. [PubMed](#)
53. Johnson RK, *et al.* 2022. *Sci Rep.* 12:19920. [PubMed](#)
54. Reitingner C, *et al.* 2022. *Front Immunol.* 13:970290. [PubMed](#)
55. Wells DK, *et al.* 2020. *Cell.* 183:818. [PubMed](#)
56. Vanoni G, *et al.* 2021. *eLife.* 10:00. [PubMed](#)
57. Wang J, *et al.* 2020. *Cell.* 183(7):1867-1883.e26. [PubMed](#)
58. Geng S, *et al.* 2016. *Nat Commun.* 7:13436. [PubMed](#)

RRID AB_314205 (BioLegend Cat. No. 302005)
AB_314206 (BioLegend Cat. No. 302006)

Antigen Details

Structure	Ig superfamily, transmembrane form (50-65 kD) or GPI-linked form (48 kD)
Distribution	NK cells, activated monocytes, macrophages, neutrophils
Function	Low affinity IgG Fc receptor, phagocytosis, ADCC
Ligand/Receptor	Aggregated IgG, IgG-antigen complex
Cell Type	Dendritic cells, Macrophages, Monocytes, Neutrophils, NK cells
Biology Area	Immunology, Innate Immunity
Molecular Family	CD Molecules, Fc Receptors
Antigen References	1. Fleit H, <i>et al.</i> 1982. <i>P. Natl. Acad. Sci. USA</i> 79:3275. 2. Stroncek D, <i>et al.</i> 1991. <i>Blood</i> 77:1572. 3. Wirthmueller U, <i>et al.</i> 1992. <i>J. Exp. Med.</i> 175:1381.
Gene ID	2214


Related Protocols

[Cell Surface Flow Cytometry Staining Protocol](#)

Other Formats

APC anti-human CD16, Biotin anti-human CD16, FITC anti-human CD16, Brilliant Violet 711™ anti-human CD16, PE anti-human CD16, PE/Cyanine5 anti-human CD16, Purified anti-human CD16, APC/Cyanine7 anti-human CD16, PE/Cyanine7 anti-human CD16, Alexa Fluor® 488 anti-human CD16, Alexa Fluor® 647 anti-human CD16, Pacific Blue™ anti-human CD16, Alexa Fluor® 700 anti-human CD16, PerCP/Cyanine5.5 anti-human CD16, PerCP anti-human CD16, Brilliant Violet 421™ anti-human CD16, Brilliant Violet 570™ anti-human CD16, Brilliant Violet 605™ anti-human CD16, Brilliant Violet 650™ anti-human CD16, Brilliant Violet 785™ anti-human CD16, Brilliant Violet 510™ anti-human CD16, Ultra-LEAF™ Purified anti-human CD16, Purified anti-human CD16 (Maxpar® Ready), PE/Dazzle™ 594 anti-human CD16, APC/Fire™ 750 anti-human CD16, TotalSeq™-A0083 anti-human CD16, TotalSeq™-B0083 anti-human CD16, TotalSeq™-C0083 anti-human CD16, PE/Fire™ 640 anti-human CD16, Spark YG™ 581 anti-human CD16, TotalSeq™-D0083 anti-human CD16, APC/Fire™ 810 anti-human CD16, GMP APC anti-human CD16, GMP PE/Dazzle™ 594 anti-human CD16, GMP PE anti-human CD16, Spark Red™ 718 anti-human CD16, GMP Pacific Blue™ anti-human CD16, GMP FITC anti-human CD16

Product Data

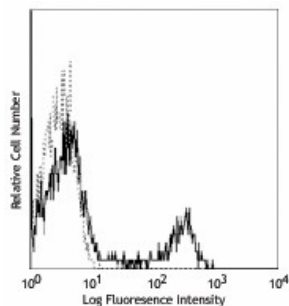

Human peripheral blood lymphocytes
stained with 3G8 FITC

For research use only. Not for diagnostic use. Not for resale. BioLegend will not be held responsible for patent infringement or other violations that may occur with the use of our products.

*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, www.biolegend.com/ordering#license). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally

similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. Unless otherwise indicated, these products are for research use only and are not intended for human or animal diagnostic, therapeutic or commercial use.

BioLegend Inc., 8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587